

END USER LICENSING AGREEMENT

Microsoft, Inc.

You have acquired a PATHFINDER™ which includes software licensed by Prentke Romich Company from one or more software licensors (Microsoft Licensing Inc. or its affiliates and SMART Modular Technologies (Systems) Inc.). Such software products, as well as associated media, printed materials and "online" or electronic documentation (software) are protected by copyright laws and international copyright treaties, as well as other intellectual property laws and treaties.

IF YOU DO NOT AGREE TO THIS END USER LICENSE AGREEMENT (EULA), DO NOT USE THE DEVICE OR COPY THE SOFTWARE. INSTEAD, PROMPTLY CONTACT PRENTKE ROMICH COMPANY FOR INSTRUCTIONS ON RETURN OF THE UNUSED DEVICE(S) FOR A REFUND. **ANY USE OF THE SOFTWARE, INCLUDING BUT NOT LIMITED TO USE ON THE DEVICE, WILL CONSTITUTE YOUR AGREEMENT TO THIS EULA (OR RATIFICATION OF ANY PREVIOUS CONSENT).**

Microsoft GRANT OF LICENSE.

The SOFTWARE is licensed, not sold. This EULA grants you the following rights to the SOFTWARE:

- You may use the SOFTWARE only on the PATHFINDER.
- **NOT FAULT TOLERANT.** THE SOFTWARE IS NOT FAULT TOLERANT. PRC HAS INDEPENDENTLY DETERMINED HOW TO USE THE SOFTWARE AND MS HAS RELIED UPON PRC TO CONDUCT SUFFICIENT TESTING TO DETERMINE THAT THE SOFTWARE IS SUITABLE FOR SUCH USE.
- **Note on Java Support.** The SOFTWARE may contain support for programs written in Java. Java technology is not fault tolerant and is not designed, manufactured, or intended for use or resale as online control equipment in hazardous environments requiring fail-safe performance, such as in the operation of nuclear facilities, aircraft navigation or communication systems, air traffic control, direct life support machines, or weapons systems, in which the failure of Java technology could lead directly to death, personal injury or severe physical or environmental damage. Sun Microsystems, Inc. has contractually obligated MS to make this disclaimer.
- No Liability for Certain damages. **EXCEPT AS PROHIBITED BY LAW, MS SHALL HAVE NO LIABILITY FOR ANY INDIRECT, SPECIAL CONSEQUENTIAL OR INCIDENTAL DAMAGES ARISING FROM OR IN CONNECTION WITH THE USE OR PERFORMANCE OF THE SOFTWARE. THIS LIMITATION SHALL APPLY EVEN IF ANY REMEDY FAILS OF ITS ESSENTIAL PURPOSE. IN NO EVENT SHALL MS BE LIABLE FOR ANY AMOUNT IN EXCESS OF U.S. TWO HUNDRED FIFTY DOLLARS (U.S. \$250.00).**
- Limitations on Reverse Engineering, Decompilation and Disassembly. You may not reverse engineer, decompile, or disassemble the SOFTWARE, except and only

to the extent that such activity is expressly permitted by applicable law notwithstanding this limitation.

- **SOFTWARE TRANSFER ALLOWED BUT WITH RESTRICTIONS.** You may permanently transfer rights under this EULA only as part of a permanent sale or transfer of the Device, and only if the recipient agrees to this EULA. If the SOFTWARE is an upgrade, any transfer must also include all prior versions of SOFTWARE.

EXPORT RESTRICTIONS. Export of SOFTWARE from the United States is regulated by the Export Administration Regulations (EAR, 15 CFR 730-744) of the U.S. Commerce Department, Bureau of Export Administration (BXA). You agree to comply with the EAR in the export or re-export of the SOFTWARE: (1) to any country to which the U.S. has embargoed or restricted the export of goods and services, which as of May 1999 include, but are not necessarily limited to Cuba, Iran, Iraq, Libya, North Korea, Sudan, Syria, and the Federal Republic of Yugoslavia (including Serbia, but not Montenegro), or to any national of any such country, wherever located, who intends to transmit or transport the SOFTWARE back to such country, (ii) to any person or entity who you know or have reason to know will utilize the SOFTWARE or portion thereof in the design, development or production of nuclear, chemical or biological weapons; or (iii) to any person or entity who has been prohibited from participating in the U.S. export transactions by any federal agency of the U.S. government. You warrant and represent that neither the BXA nor any other U.S. federal agency has suspended, revoked or denied your export privileges. For additional information see <http://www.microsoft.com/exporting/>.

END USER LICENSING AGREEMENT

SMART Modular Technologies (Systems) Inc.

- You may use DECTalk RUNTIME SOFTWARE TECHNOLOGY **only** on the PATHFINDER™
- You may not COPY DECTalk RUNTIME SOFTWARE except for archival use.
- You may not reverse engineer, decompile, or disassemble the DECTalk RUNTIME SOFTWARE.
- The Title to DECTalk RUNTIME SOFTWARE, which is a component of the PATHFINDER, does not pass to the END USER.
- SMART disclaims any liability for any damages, whether direct, indirect, incidental, special or consequential arising from the use of the DECTalk RUNTIME SOFTWARE.
- At the termination or expiration of the END USER AGREEMENT, you must discontinue to use, and destroy or return to Prentke Romich Company, all PATHFINDERS and all archival or other copies of DECTalk RUNTIME SOFTWARE that contain any part of DECTalk TECHNOLOGY.
- The END USER has no right to transfer or distribute the DECTalk RUNTIME SOFTWARE to others.

Pathfinder

Operator's Manual

Setting Up and Using Pathfinder's Basic Tools

PRENTKE ROMICH COMPANY

Corporate Headquarters:

1022 Heyl Rd. • Wooster, Ohio 44691

Telephone: (330) 262-1984

Fax: (330) 263-4829

E-Mail Address: info@prentrom.com

Web Site Address: www.prentrom.com

Worldwide Headquarters

Prentke Romich International, Ltd.

Whitegates • Swinstead • Grantham • Lincs. • NG33 4PA • UK.

Fax: 44 1476 550872

E-Mail Address: info@prentromint.com



The cover photograph of multiple generations of stars in the Tarantula Nebula
was taken by the Hubble Space Telescope.

The picture was downloaded from the Space Telescope Science Institute/Association of Universities
for Research in Astronomy, Inc. web site on the internet.

☞ Data can be uploaded from this unit to a personal computer for storage on a hard disk, floppy diskette, or other storage medium. Be sure to make multiple copies of all important data to protect against accidental loss.

Pathfinder, Liberator, DeltaTalker, Language Activity Monitor (LAM), Memory Transfer Interface (MTI), PRC Infrared Receiver (IRR) are trademarks of Prentke Romich Company. Use of the name "Liberator" with permission from Liberator, Ltd.

Unity is a registered trademark of Semantic Compaction Systems. Unity uses the methodology derived from the concepts covered by U.S Patent Numbers 4,661,916; 5,097,425; 5,210,689 and subsequent patents and applications of Bruce Baker.

Unity is a registered trademark, in the USA, UK and Canada, of Semantic Compaction Systems.

MinTalk is a trademark of Semantic Compaction Systems.

DECtalk software is a trademark of FORCE Computers, Inc.

Windows, Microsoft Word, Access Pack and Microsoft Works are registered trademarks of Microsoft, Inc.

Corel WordPerfect is a registered trademark of Corel, Inc.

Access DOS is a registered trademark of IBM.

KE:NX is a registered trademark of Don Johnston Company

Apple, Macintosh, PowerBook, AppleTalk and SerialKeys are registered trademarks of Apple Computer, Inc.

IBM is a registered trademark of International Business Machines.

T-TAM (Trace Transparent Access Module) was developed by the Trace Research and Development Center, University of Wisconsin-Madison.

HeadMaster is a trademark of Origin Instruments, Inc.

Tracker 2000 is a trademark of Madenta, Inc.

DECtalk song programming instructions, charts and speech tips are used with permission of Digital Equipment Corporation.

©Prentke Romich Company 2000, 2001

All rights reserved. Under copyright laws this manual may not be copied, in whole or in part, without the written consent of the Prentke Romich Company.

ELECTRONIC INTERFERENCE INFORMATION

This device complies with Part 15 of the FCC Rules. Operation is subject to two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning! Changes or modifications to this product which are not authorized by Prentke Romich Company could exceed FCC limits and negate your authority to use this product.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. You can determine whether this product is causing interference in your radio or television by turning this product off. If the interference stops, it was probably caused by this product or one of its accessories. You can attempt to correct the interference by using one or more of the following measures:

- Turn the television or radio antenna until the interference stops.
- Move this product to one side or the other of the television or radio.
- Move this product farther away from the television or radio.
- Plug this product into an outlet that is on a different circuit from the television or radio; that is, this product should be controlled by different circuit breakers/fuses from the television or radio.
- If necessary, contact a Prentke Romich Company service technician for assistance.

This Class B digital apparatus meets all requirements of the Canadian Interference-Causing Equipment Regulations. Cet appareil numérique de la classe B respect toutes les exigences du Règlement sur le matériel brouilleur du Canada.

Prentke Romich Company • 1022 Heyl Rd. • Wooster, OH 44691

ACKNOWLEDGEMENTS
Pathfinder Manual: Basics

13113L3.0

CE ENVIRONMENTAL INFORMATION: Use: 0° - 45° C; Storage: (-20°) - 50° C; 95% maximum relative humidity @ ambient temperatures less than 40° C.



Attention! Consult accompanying documents.



Type BF Applied Part: F type applied part complying with the specified requirements of 60601-1 to provide protection against electric shock, particularly regarding allowable leakage current.

Not intended to be an emergency call device or sole communication device.

Turning the Pathfinder On for the First Time

To conserve batteries, the entire unit was turned OFF at the factory before it was shipped. To use it, you must turn the unit ON.

To Turn the Pathfinder ON

On the left front of the case:

Press the  (I/O) and the  (TOOLS) buttons down at the same time. When the display lights up and you see the Pathfinder name, release the buttons.

Recommendation

To Avoid Static Discharge

You can avoid static discharge by turning the Pathfinder OFF before you plug anything into or unplug anything from the connectors on the case—this includes switches, headpointers, and **especially** the battery charger.

When you follow this recommendation you are helping to preserve the device's memory by safeguarding its electronic circuits.

Table of Contents

GETTING ACQUAINTED WITH THE PATHFINDER

What's on the Front of the Pathfinder?.....	11
What's on the Top of the Pathfinder?	17
What's on the Back of the Pathfinder?	19
What's on the Right and Left Sides of the Pathfinder?	21
About the Batteries	23
About the Backlight Shut-Off	25
About Minspeak.....	27
Icon Prediction Shows You Where Messages are Stored.....	29
Examples of Icon Associations.....	31
About Selecting Your Own Icons to Represent Messages	33
About Activities	35
Locating Keys on the Keyboard	37
About the Toolbox.....	39
Using Menus in the Toolbox	41
Change the Speaker and Earphone Volume	43

SETTING UP A USER AREA AND A SELECTION TECHNIQUE

The USER AREA Menu	47
Select a User Area.....	49
About Overlays	51
Rename a User Area.....	52
Copy a User Area	53
Erase a User Area.....	54
About Computer Access Mode.....	55
The Selection Techniques.....	57
Using the SELECTION TECHNIQUE Menu	59
Keyboard.....	61
Infrared Headpointer	63
More About Your IR Headpointer	65
1-Switch.....	67
2-Switch.....	69
Joystick	71
Quarter Row/Column Scan	73

Some Options in the SELECTION TECHNIQUE MENU	75
<i>IR Headpointing Options:</i>	
Averaging.....	77
Horizontal and Vertical Adjustment	79
Pointer Color.....	81
Battery Information.....	82
Connecting the HeadMouse® to the Pathfinder	83
Connecting the Madentec Tracker 2000® to the Pathfinder	84
<i>Some Scanning Options</i>	
1-Switch—Scan Mode 1	85
1-Switch—Scan Mode 2	86
2-Switch—Scan Mode 1	87
2-Switch—Scan Mode 2	88
Scan Direction	89
Scanning Speed	90
Acceptance Time	91
Release Time.....	92
Activation Delay	93
Auto Repeat Rate.....	94
Predictive Selection	95
Modify Predictive Keys.....	97
Scan Start Position: Center or Top Left	99
About Linear Scan Cut-Off	100
WORD PREDICTION	
About Minspeak and Word Prediction.....	103
The Benefits of Minspeak and Word Prediction Together	105
Turn Word Prediction On and Off.....	107
Spell and Select Words.....	109
Set Up Word Prediction in the WORD PREDICTION Menu	111
Predict in Minspeak	112
Recency/Frequency.....	113
About Predictor Keys.....	115
The Word Prediction Dictionary.....	117
Add a Word to the Word Prediction Dictionary	118
Delete a Word from the Word Prediction Dictionary.....	118
Review Words in the Word Prediction Dictionary	119
About the Computer "Backspace" Command.....	121
The Word Prediction Display Area.....	122
STORING MESSAGES	
Store Digitized (Recorded) or Synthesized (Text-to-Speech) Messages	125
About Core Vocabulary and Activity Row Vocabulary	127

About SPELL/MINSPEAK, GO TO SPELL, GO TO MINSPEAK and GO TO COMPUTER Tools .	128
Plan Before You Begin to Store	129
Store a Recorded Message	131
Erase a Recorded Message	132
Store a Synthesized Message	133
Edit a Synthesized Message	134
Erase a Synthesized Message	135
Some Helpful Storing Hints	136
Storing and Retrieving Using the 64 and 32 Location Overlays	137
About Themes	139
About Using INSERT TOOL to Store Theme Tools on your Overlay	140
Assign PICK 1 (or 2) ICON THEME to a Key or Store it under an Icon Sequence	141
Assign SET 1 (or 2) ICON THEME to a Key or Store it under an Icon Sequence	143
About EXIT THEME	145
About Single Meaning Icons	147
About Creating an Activity	151
Create a New Activity—An Example	153
Dynamic Activities	157
Erase a Dynamic Activity	159
Review All Dynamic Activities	159
Modify an Activity	160
Linking Activities/Creating Sub-Categories	161
Copy Activities from one User Area to Another	163
Using the REVIEW VOCABULARY Menu	165

PERSONALIZE YOUR OVERLAY

About Personalizing Your Overlay	171
The ASSIGN KEY Menu	173
About Copying Tools from the Toolbox to your overlay	175
Assign a Tool to a Key in Minspeak Mode	176
Assign a Minspeak Icon to a Key	177
Assign a Tool to Spell Mode	178
Assign a Tool to Computer Mode	179
Erase a Tool from a Key	181
Copying a Tool and Storing it Under and Icon Sequence	182
About Macros	183
Assign a Macro to a Key	184
Store a Macro Under an Icon Sequence	185
Erase a Macro	185
Swap Keys	186
Hide/Show Keys	187
Auditory Prompts and the AUDITORY PROMPTS Menu	189
Prompt Type: Key or Content	190

Key Prompts and Category Prompts	191
Mid-Sequence Prompts	193
Some Options in the FEEDBACK Menu	195
Spell One Word.....	197
Auto Spell	198
Icon Tutor.....	199
System Colors.....	200
The MAINTENANCE MENU.....	203
Reset Individual Menu Options to their Default Settings.....	205

APPENDICES


Appendix A: All Tools List and Explanations of Tools.....	209
Appendix B: Data Logging: Using the LAM	224
Appendix C: How to Use the <i>Review to Output Option</i> in Various Pathfinder Menus.....	232
Appendix D: Pre-Stored Songs in the Pathfinder.....	236
To Copy a Pre-Stored Song onto your Overlay.....	237
Appendix E: How to Program and Store your own Songs on your Overlay	239
English Vowel Phonemes	253
English Consonant Phonemes.....	254
Pitch Notes and Pitch Values.....	255
Pitch Value Code	257
Time Values of Notes.....	259
Appendix F: Helpful DECTalk Speech Tips	261
Words with Alternate Pronunciations.....	263
Abbreviations Recognized by DECTalk Speech.....	265
Produce a Long "A" in a Sentence.....	268
Creating Pauses	269
How DECTalk Speech Handles Numbers	269
Appendix G: Dependent Icons.....	271
Appendix H: External Speaker and Microphone Specifications.....	275
Appendix I: Creating Macros that Contain Menu Options.....	276
Appendix J: Computer Commands for IBM and Macintosh Keyboards and Keypads.....	279
Appendix K: Save Your Setup Settings.....	283
Appendix L: Troubleshooting	285
INDEX	293

Getting Acquainted with the Pathfinder



KEYBOARD AND TOOLBOX

Front of the Pathfinder: Turned Off

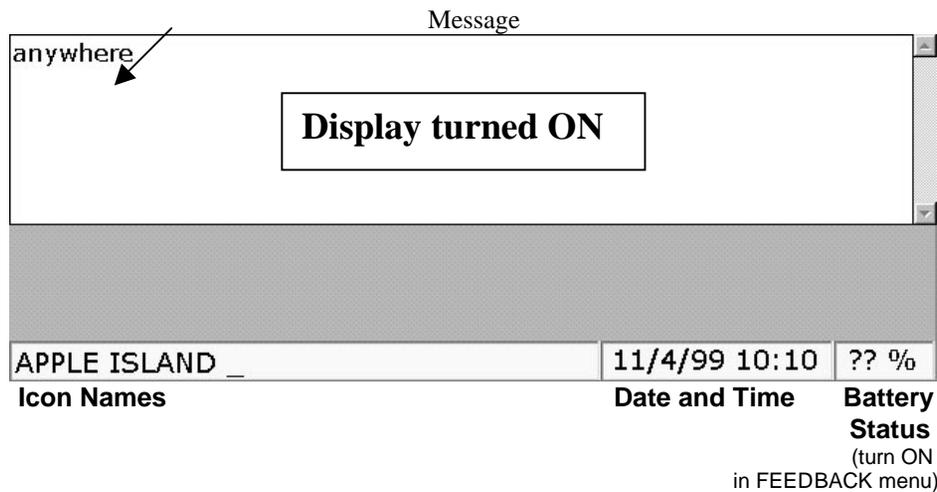

What's on the Front of the Pathfinder?

Display Screen

The display screen is where you will see text. You can choose large, medium, or small sized text. The size of the text you select dictates how many lines of text will fit on the screen.

You can press on the display to make menu selections or you can use the keyboard.

Directly below the display screen is the status bar. This is where you can see icon names and check for directions or instructions about what you are doing. If you are using a Notebook, the status bar will display the icon of a small notebook.

The blank area below the display area is where you will see your Activity Row or Word Prediction choices when those options are On.

Buttons

On the Front Left Side:



I/O: Press this button to turn the Pathfinder On or Off.




Tools: Press this to open and close the Toolbox.



System Status: Press this to see the Pathfinder's current settings.

There is also a TOOLBOX option on this screen that allows you to access all the Toolbox menus from the display screen without removing the keyguard and overlay and going to the Toolbox.

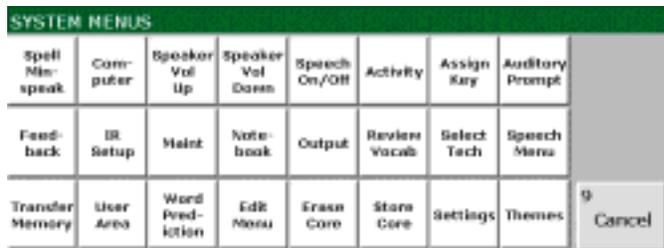

User Defined Key: This key defaults to CHOOSE ACTIVITY. If you do not want CHOOSE ACTIVITY, you can program it with whatever command you want. The command will be set in all three modes—Minspeak, Spell and Computer.



Battery Indicator: When the LED in the battery box comes on, charge the battery. When you plug in the charger, the LED in the lightning volt will flash. A rapid flash indicates fast-charging; a slow flash indicates that the battery is charged but it will continue to trickle-charge as long as the charger is plugged in.

On the Front Right Side



User Areas: The areas are numbered 1-6 and each number has an LED (Light Emitting Diode) next to it. The area you are currently using will be lit.

ABC: When this LED is lit, you are in *Spell* mode and can spell messages.



● **Minspeak Symbol:** When this LED is lit, you are in *Minspeak* mode and can retrieve stored messages.



● **Computer Icon:** When this LED is lit, you are in *Computer* mode and can use the pre-set computer keys for computer emulation.




● **Shift/Caps Lock:** Press Caps Lock on your overlay and this LED will light. Press SHIFT and this LED will flash for one activation.



● **Microphone** The recording microphone is in the upper right corner of the Pathfinder case. Speak into the microphone when you record a message.



● **Headpointing IR Window:** Infrared signals are received from your *headpointer* through this window.

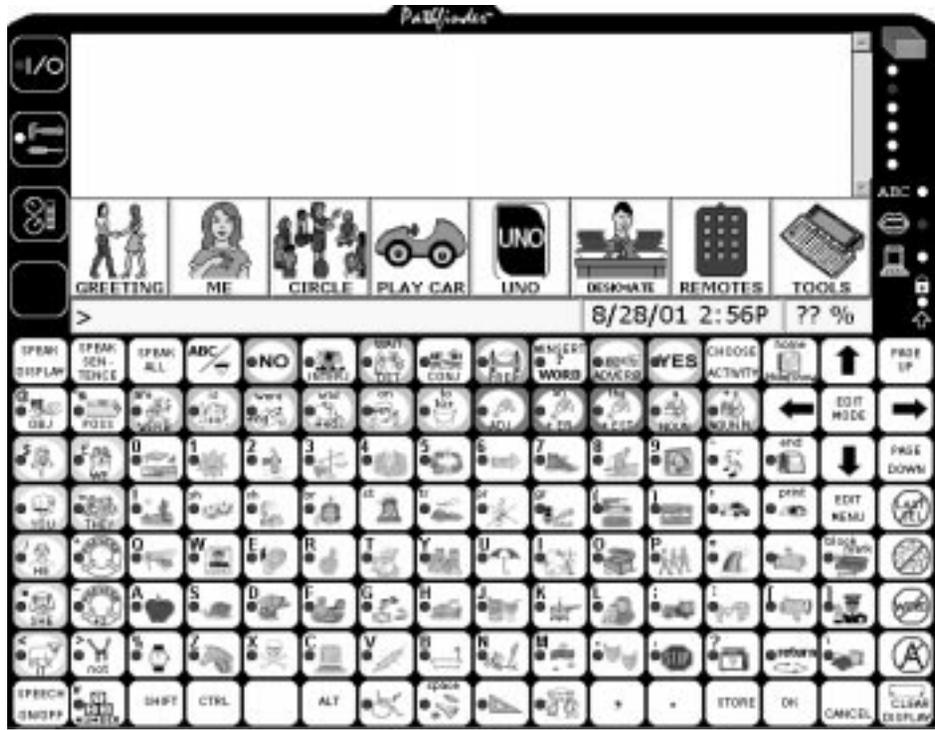

Unity Default Overlay

GETTING ACQUAINTED

Pathfinder Manual: Basics

The Keyboard:

Keyboard and LEDs: The Pathfinder comes with 128 keys, each with an LED (Light Emitting Diode) that will light when the key is activated.

Keyguard: The keyguard fits over the keyboard and outlines the keys on the keyboard. It is made of durable plastic. Remove the keyguard by sliding it up from the bottom of the case and lifting up.


Overlay Protector: This is a clear sheet of vinyl that fits under the keyguard and on top of the overlay. It helps protect the overlay from spills, crumbs, drooling, etc.

Overlay: The Pathfinder is shipped with a Unity overlay in place. This overlay has icons but no printed words on it. It will work with the Unity 1-Hit or the Unity Sequenced-Hit vocabularies. The overlay package also includes a Unity I-Hit overlay, Unity 32-1 and 2 hit overlays and a computer mouse access overlay that you place on top of your Unity overlay.

Toolbox: Remove the keyguard, overlay protector and overlay to see the Toolbox. The Toolbox contains "tools" that you can use to increase the Pathfinder's usefulness. The names of the tools are printed on the keys. The Toolbox has 128 keys and can be accessed the same way as the keyboard.

➡ There may be times when you would like to go to a menu without having to remove the keyguard and overlay and opening the Toolbox. You can do this by using the SYSTEM STATUS key on the top left front of the Pathfinder.

Activate the **SYSTEM STATUS** key, . On the display screen, activate the **TOOLBOX** key. Select the menu you want. You can also activate the **STATUS** key twice to go directly to the Toolbox.

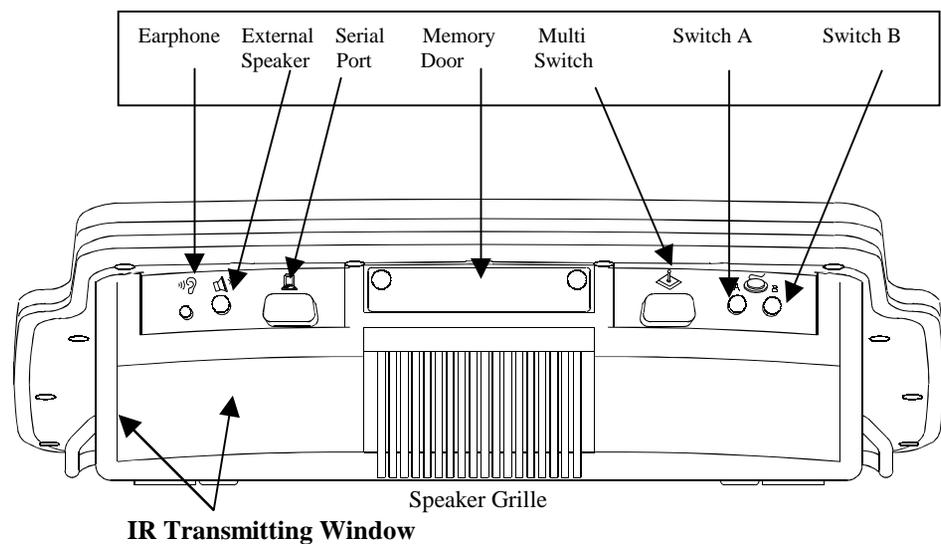

What's On the Top of the Pathfinder?

The back of the Pathfinder contains the connectors for things you can plug into the Pathfinder.

Earphone jack: Plug an earphone in here. See **Appendix H** for earphone specifications.

External Speaker: Plug an external speaker in here. See **Appendix H** for speaker specifications.


Memory Card Door: Underneath this door is the PCMCIA card that contains the system software for your Pathfinder. Unless you have received a new system software card from PRC and need to install it, this door should remain securely in place at all times.

Serial Port: If you want to make a memory transfer, plug in a computer or plug in a serial printer, use the computer serial port. ➡ If you are going to use a multi-switch (e.g., a joystick) plug it into the **Multi-Switch** connector.

Switch 1 and Switch 2: Plug a switch into one of these ports.

Speaker Grille: The internal speaker for the Pathfinder is under this grille. This area should not be covered up if you want optimum speaker volume and speech.

IR Transmitting Window: Infrared signals are sent out to other devices and appliances through this area on the Pathfinder.

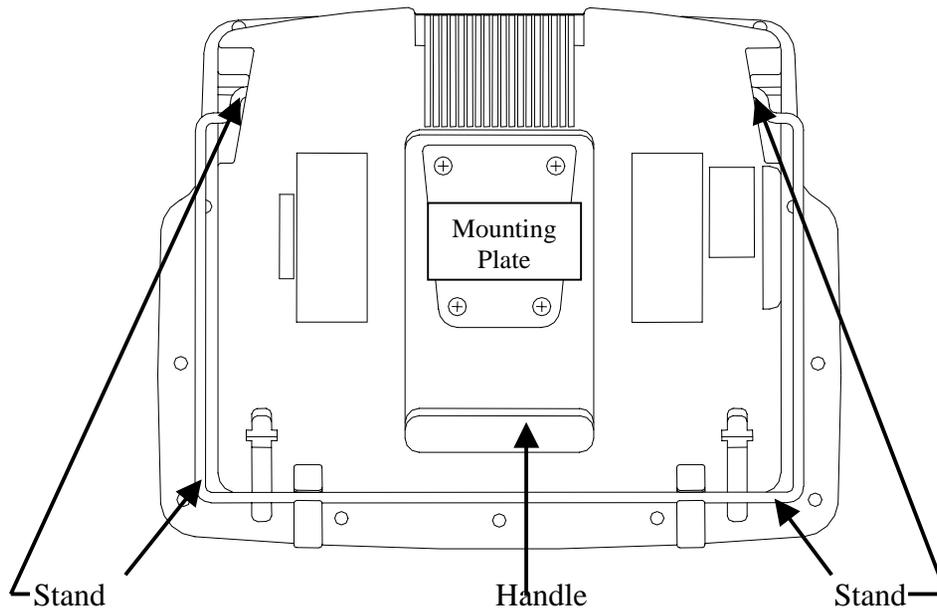

When you are using the stand, you can change its position/angle.

Carefully remove the sides of the stand from the positional holes they are currently in. Replace them in another set of holes (there are 3 holes on each side).

 **Warning!**

It is possible to pinch your fingers while re-positioning the stand. Remove the sides of the stand from their positional holes with care.

What's on the Back of the Pathfinder?

Turn the Pathfinder face down. You will see:


Wheelchair/Table Stand Mounting Plate: This plate slips into the PRC wheelchair mount or the PRC table stand to hold the Pathfinder upright.

Handle: The Pathfinder has a molded handle on the back.


Built-In Table Stand: The Pathfinder has its own table stand. Pull it out from the case. It can be set at one of three different angles by repositioning the stand in the different positioning holes.

Carry Strap Holders: Connect the hooks on your carry strap through these holders. Unscrew the hook, slip it around the holder and screw the hook back together again. See the label on the strap bag.

Left side of Pathfinder

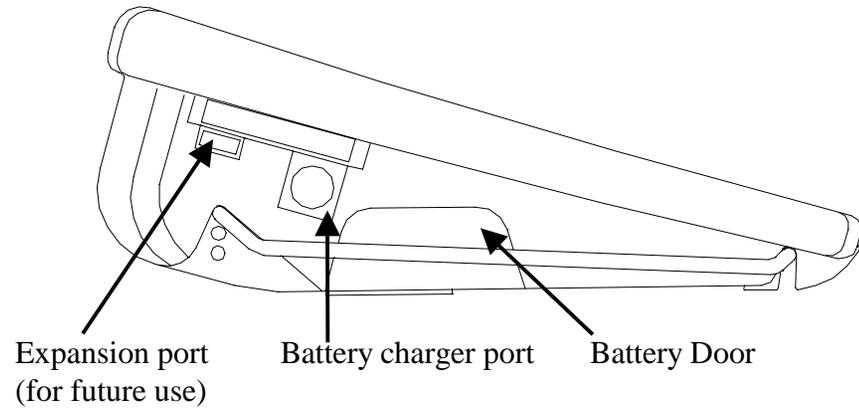

Right Side of Pathfinder

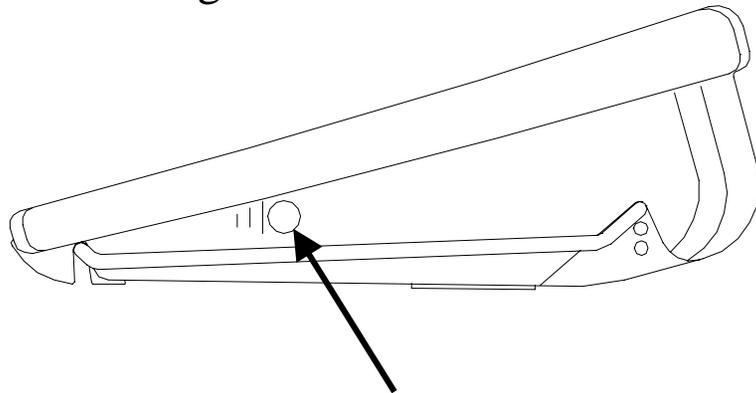

Learning IR Window

What's on the Left and Right Sides of the Pathfinder?

On the left side of the Pathfinder:

Expansion Port: This port is for future development and does not work at this time.


Battery Door: The battery pack is under this door. If you need to replace your battery pack, you will remove this door, pull the old battery out, slide the new battery in and replace the battery door.

Battery Charger Port: Plug the battery charger that came with your Pathfinder into this port.  See warning below.

 **WARNING!** Use *only* the battery charger that came with your Pathfinder. Other chargers may damage the batteries or the Pathfinder.

On the Right Side of the Pathfinder

IR Learning Window: When you teach remote control signals to the Pathfinder you will point the infrared area of your remote control at this area on the Pathfinder. The easiest way is simply to lay the remote control on a table next to the Pathfinder with the remote control IR area pointing to the IR learning window on the Pathfinder. Just press the remote keys to send the correct signals to the Pathfinder. See the *Power Tools* manual for complete learning infrared instructions.

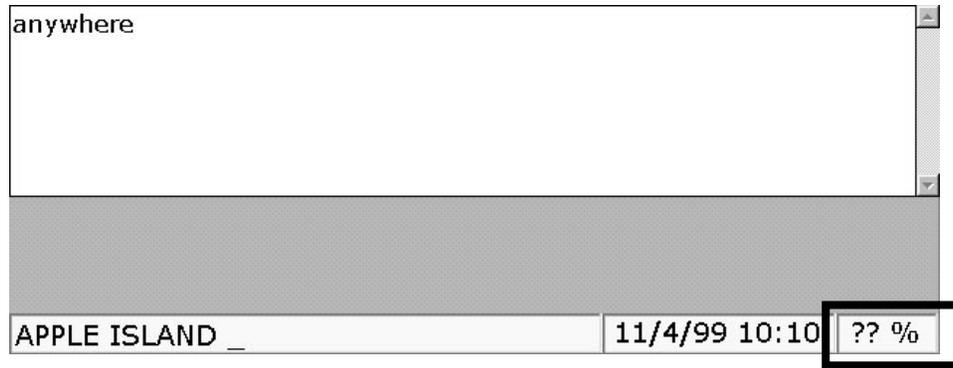

**Battery
Status
Area**

Battery Status Area on the Display Screen

About the Pathfinder Batteries



Always use the battery charger that came with the Pathfinder.

Any other charger may damage your batteries.

If you turn **Battery Status On** in the FEEDBACK menu, you will see a percentage of the full battery charge displayed in the battery status area (which is the far right area) in the status bar at the bottom of your display screen. For example, you might see, "78%" which means your batteries are 78% full.

The battery status area turns blue and displays the word, **Charger**, when the charger is plugged in and the batteries are charging. When the batteries are fully charged, the battery status area turns green and displays the word, **Full**.

When the batteries get low, the Pathfinder sounds a **Low Battery Warning**. This is a two-toned running beep. The battery status area on the display turns red and displays the amount of charge remaining in the batteries, for example, 14%. The **Low Battery** LED on the front of the case also turns on.  When the Low Battery warning sounds, you can no longer do any storing.

The Low Battery warning starts when battery capacity drops below 15%. The warning beeps are repeated at every 2% drop in battery capacity. When the battery capacity drops to 4%, the warning is repeated every 20 seconds. When the batteries are dead, the Pathfinder powers down completely and you cannot use it.

 When you hear the Low Battery warning, plug in the battery charger. **When the Pathfinder is in use during charging**, it will take 8-10 hours to fully charge a dead battery. **When the Pathfinder is not in use during charging**, it will take 4-5 hours to fully charge a dead battery.

A fully charged battery should last for 8-10 hours of normal use.

MAINTENANCE MENU			1		
1	System Lock <input type="text" value="OFF"/>	4	Set Time/Date	7	More Items
2	Backlight Shut-off <input type="text" value="10 Minutes"/>	5	Time/Date Format	8	OK
3	Erase All Vocabulary	6	Delete Imported Icon	9	Cancel



About the Backlight Shut-Off

Summary:

The Backlight Shut-Off can be set for either 5 or 10 minutes. If 5 (or 10) minutes pass with no key activations, the Pathfinder will turn the display backlight off to conserve battery energy. Activate a key or switch to wake it back up. The Backlight Shut-Off option is in the MAINTENANCE MENU in the Toolbox.

The Pathfinder automatically shuts down the backlight on your display if you don't make any key activations for 5 minutes. This helps to conserve the batteries. Just activate a key or switch to turn the display on again. After 30 minutes of no key activations the pathfinder will go to sleep. To wake it up, press any key or switch. When the Pathfinder powers down, it automatically saves any storing or programming you have done.

If you would prefer a longer wait before the backlight powers down, you can set the Backlight Shut-Off for 10 minutes. Open the MAINTENANCE MENU in the Toolbox and change the *Backlight Shut-Off* option to 10 minutes. (See pages 39-41 for information about using the Toolbox.) You can also choose to turn *Backlight Shut-Off* to *Off*. This means the Pathfinder will never power down by itself.

➤ **If You Use an Infrared Headpointer:** an infrared headpointer will always be able to wake up the system as long as the Pathfinder has not been turned off using the I/O button on the front of the case.



➤ **To turn the Pathfinder completely Off,** press the I/O button on the left front of the case. If you turn the Pathfinder off with this button, you must turn it back on by pressing this button again.



I



WANT

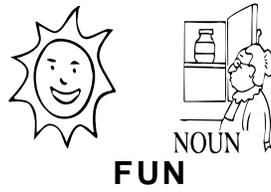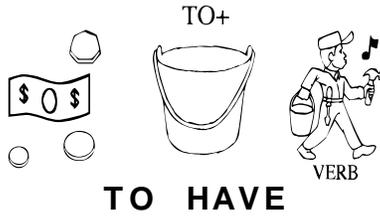

Minspeak Icon Sequences

About Minspeak


The Pathfinder uses a patented vocabulary-coding method called Minspeak to help you communicate.

Minspeak is a means of coding vocabulary to a small set of icons that are rich in meaning. Minspeak's use of icons taps into a person's natural tendency to associate multiple meanings to pictures. By coding a large vocabulary of words, phrases and whole messages using this small set of icons, Minspeak offers many advantages:

- ◆ it is logical in how it represents a person's language;
- ◆ it is easy to learn because it fits into a person's natural ability to associate many meanings to one picture;
- ◆ it promotes language development by providing a rich vocabulary of words that you can combine into sentences;
- ◆ because the small set of icons remains the same on a single overlay, you can, with practice, learn to recall your vocabulary automatically;
- ◆ it promotes further language development by offering many opportunities for language-learning experiences.


These advantages lead to the ultimate benefit of using Minspeak: it allows for and promotes independent, effective communication.

ICON PREDICTION

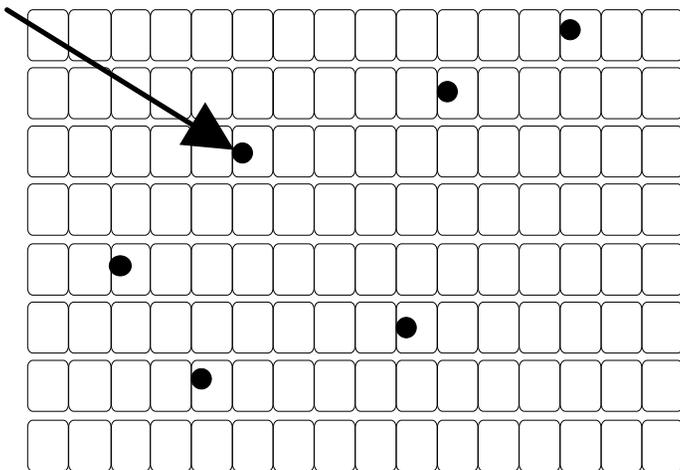

Icons that begin messages are lit.

Select an icon.

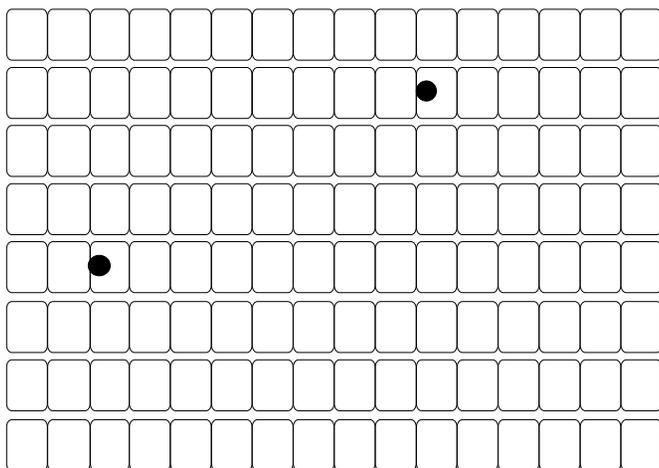

Activate an icon that is lit.

Pathfinder speaks the message.

**Icon Prediction Lights only those Icons where
Messages are Stored**

Icon Prediction Shows Where You Have Messages Stored

Summary:

When Icon Prediction is ON, the icons where messages are stored are lit, other icons are not lit. Choose a **Low** or **High** intensity for the lights.

If you have activated any keys on the keyboard, you have probably noticed that some LEDs on the keyboard remain lit while other LEDs go out. This is called "Icon Prediction". The icons with lighted LEDs tell you that they are possible selections for messages that are part of the vocabulary that was pre-stored in the Pathfinder software. The default for Icon Prediction is On.


When you turn **Icon Prediction On** (in the FEEDBACK menu) you can select a **Low**, **Medium** or **High** intensity for the LEDs. If you select **Low**, the LEDs will not be overly bright. If you select **High**, the LEDs will be very bright.

If you press a key that is lit, the name of that icon will appear in the status bar area of the display. Other keys on the keyboard will either go out or remain lit. If you select another lit key, its icon will appear in the status bar area and the device will speak the message that is stored under the sequence you just activated.

Once a message is spoken, the keyboard returns to lighting all the icons. When you press a key again, some keys will remain lit and some will go out.

If you are using scanning you will want to adjust the Scan LED Intensity so that it is easy for you to see the key you are scanning or pointing to.

Icon Prediction is a good tool to use to help you learn your vocabulary.

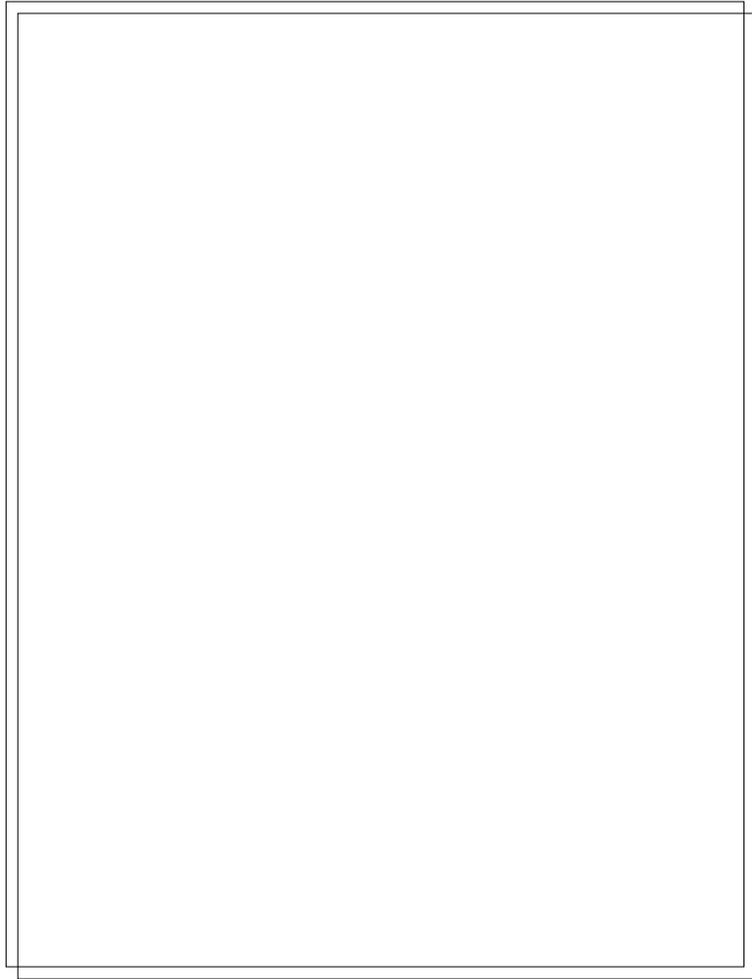

Icon Associations

An Example of Icon Associations

Minspeak uses **icon combinations** to represent a vocabulary. For example, suppose you want to store the word, "Where." Minspeak allows you to choose several icons to represent the word. If you choose a two-icon sequence, you will find that you can speak the word with only two key activations. This results in a 60% savings in key activations over spelling out the word.

Minspeak also makes it easy for the user to remember which locations he/she must activate to speak a message. Here is where the power of the icons becomes apparent. When storing a message you pick icons that will remind you of its contents. With the word "anywhere", you might



choose the sequence,  plus  . In this case the icons were chosen

because  represents the letter "A" and  represents an island (or 'someplace').

You will notice that in the future when you want to speak the word "anywhere" you can simply activate the icons that you selected while storing the message. Not only will this require fewer activations, but with practice it will become automatic to recall. If the association of "where" with the map of an island is not understood to start with, there is a good teaching activity built into this kind of Minspeak sequence to help a person understand and to create a rewarding language experience.

It is important to realize that the Minspeak icons are rich in meaning. The meanings associated with the icons can be the meanings you learn from your vocabulary program and they can also have very individual meanings. For instance, the icon,  might be used in different sequences to mean 'umbrella', 'weather', 'under', 'rain', 'whether', or any number of other things that the icon makes you think of.

ICONS SHOULD BE
easy to learn
easy to remember
easy to teach

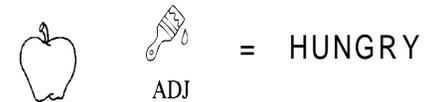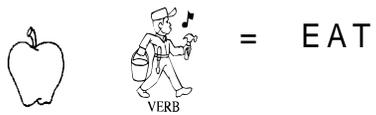

Choosing Your Own Icons

About Selecting Your Own Minspeak Icons to Represent a Message

Summary:


When choosing your own icons for messages consider the following things: the "rules" for your vocabulary program; all the associations an icon may have; the category your message may belong in; how easy the association between icon and word will be to learn, to remember and to teach; what icons might be combined to represent a wide range of messages.

There is a Pathfinder Unity overlay in place on your device. We suggest you use this overlay for practice.

A big consideration in choosing icon sequences is that they are logical, easy to learn and easy to remember.

If you are using a vocabulary program you will want your own sequences to work the same way as the sequences work in your program. To achieve this you will need to understand the "rules" of your vocabulary program. Spend a little time thinking about multiple associations for icons. This will help you to understand how various icons can be combined to represent a word. Also spend some time thinking about what kind of a 'category' your message might fit into. Organizing a vocabulary into categories is something human beings do naturally so having a vocabulary program that is stored in the same way can help you learn and remember your sequences more easily. ➡ For more information about this, see your *Unity* manual. Also see the **STORING MESSAGES** chapter in this manual beginning on page 125.

If you are not using a vocabulary program, you will want to achieve the same goals, i.e., logical sequences that are easy to learn and remember. To achieve this you will want to devise a good system for storing vocabulary. Picking icons is part of this system.

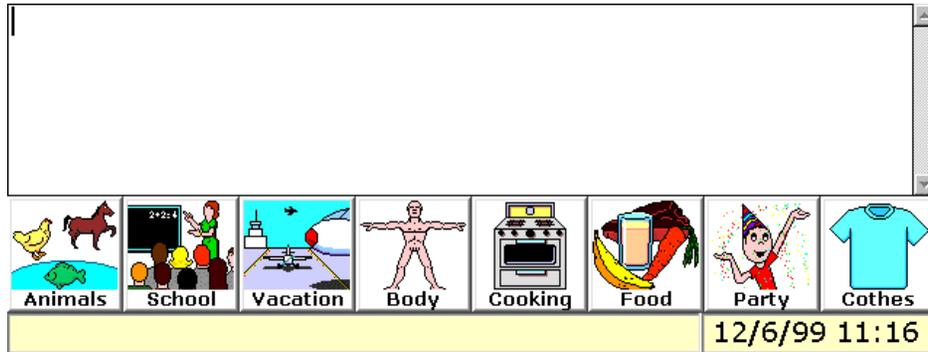

The Activity Row Showing a List of Activities
(the Activities on your Pathfinder may be different)

About Activities

When Activities are turned On (in the ACTIVITIES menu), you will see a row of keys/icons at the bottom of your display screen. This is the Activity Row. Activities can be used in a number of ways. You can create activities that allow you to access words that are centered around a particular theme, such as "School," "Work," "Party," etc. Or, as some of the Pathfinder Unity overlays do, you can have an activity row that contains the picture producing nouns from your overlay, or all the action words, or whatever you wish.

Activities can be accessed two ways:

Dynamic: these activities are brought up by activating an icon or an icon sequence on your overlay. There are pre-stored dynamic activities in the Unity vocabulary, but their icons are not shown in the Activity Row.

List of Activities: you can create activities that appear in a row (or rows) on your display screen when Activities are turned On. You can then choose an Activity from the row.

➡ The User Defined Key (on the left front of the Pathfinder case) defaults to CHOOSE ACTIVITY. If you do not assign something else to this key, you may use it to select activities from the Activity row.

➡ You can also assign the CHOOSE ACTIVITY tool to a key or store it under an icon sequence, if you would like to. See pages 173-182 to learn about this.

➡ The Unity overlay also has CHOOSE ACTIVITY assigned to key M1.

Key Location Codes

C O L U M N S

	A	B	C	D	E	F	G	H	I	J	K	L	M	N	O	P
R																
	1															
	2															
O	3															
4							■									
W	5															
	6															
S	7															
	8															

Location **G-4**
Column **G** Row **4**

	A	B	C	D	E	F	G	H	I	J	K	L	M	N	O	P
1																
2															■	
3																
4																
5																
6																
7																
8																

Location **O-2**
Column **O** Row **2**

Key Coordinates


Use Key Coordinates to Locate Keys on an Overlay

There are 16 columns of keys going across a 128-location overlay. There are 8 rows going down the overlay. These numbers are also true in the Toolbox.

The columns are labeled from "A" (far left column) to "P" (far right column). The rows are numbered from "1" (the top row) to "8" (bottom row). Using these guides as coordinates, you can find any key on the overlay.

For example, key A1 is the top key in the far left column. Key P8 is the bottom key in the far right column. Key H4 is near the middle of the keyboard.

Key coordinates are used frequently in this manual to help you locate keys easily.

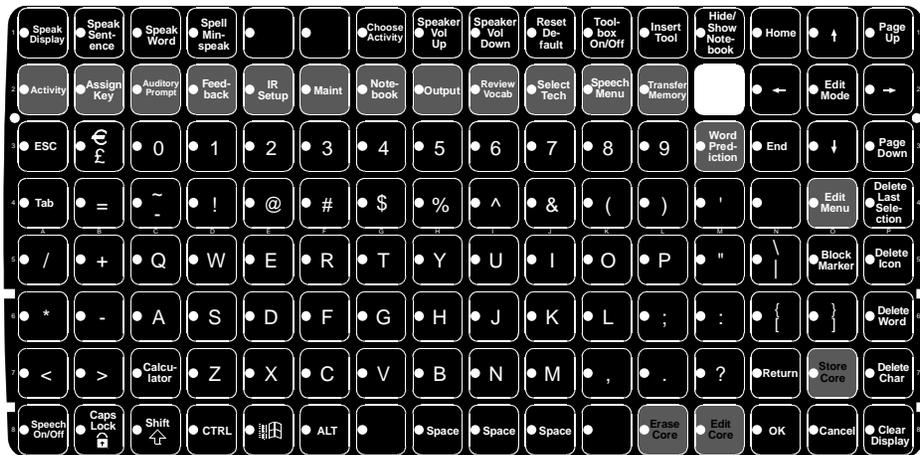

The Pathfinder Toolbox

About the Toolbox

Summary:

Remove the keyguard and overlay to see the Toolbox. To turn the Toolbox On, press the TOOLS key in the Status Area of the Pathfinder.


Remove the keyguard and overlay from the Pathfinder. What you see is the Toolbox. The Toolbox has 128 keys. Some of these keys contain tools that you can use to make the Pathfinder more useful to you. Some of the keys contain the letters and numbers that you use when you want to spell something. Some tools allow you to do specific things-like turn the speech on or off or erase the display. Other tools are multi-purpose-like the menu tools which give you a number of options to choose from.

Notice that some keys are black with white letters, some are gray with white letters and some are gray with black letters. **MENUS** are gray with white letters. Most of the Pathfinder menus are in the second row of keys. The WORD PREDICTION and EDIT menus are in different rows. The gray keys with black letters are the keys you use to store, erase and edit messages. All the rest of the keys are black with white letters and they are either Pathfinder tools or keyboard keys for the Pathfinder.

The Toolbox must be turned On before any of the keys will work. To turn the Toolbox On, press the TOOLS key on the front of the Pathfinder. When the Toolbox is On, the small red light on the TOOLS key is lit. This goes out when the Toolbox is turned Off.

The tools in the Toolbox are tools that you generally use most often. There are lots more tools available for use, but there aren't enough keys in the Toolbox for each individual tool. For a list of all the Pathfinder tools, see **Appendix A**.

Tools can be copied from the Toolbox to your overlay and stored under icon sequences or assigned to keys so that you can use them without having to go to the Toolbox. Pages 173-182 discuss these possibilities.

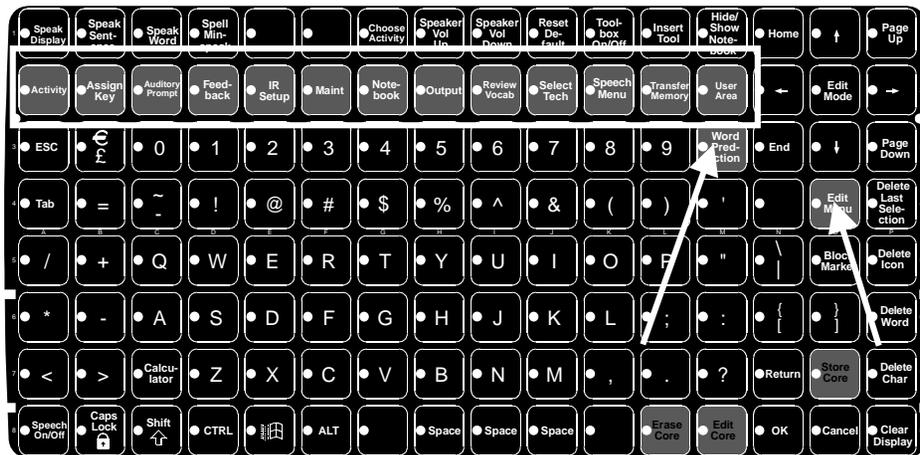

Menus in the Toolbox are Light Gray

Using Menus in the Toolbox

Remove the keyguard and overlay from the Pathfinder. Press the TOOLS button on the front of the case. When the light on the TOOLS button is lit, the Toolbox is On.

Using the menus in the Toolbox is easy. The menus are arranged alphabetically. Menu keys are light gray. The name of each menu is printed in white on the gray key.

Press the FEEDBACK menu. Let's look for a moment at the menu. The name of the menu and any instructions for using it are at the top of the menu. There are nine keys on the menu and each key is numbered. For example, the *Speaker Volume* key has the number "1" in the upper left corner of the key. *OK* and *Cancel* are number "8" and number "9". On any menu, OK and Cancel are #8 and #9. (You can also put *OK* and *Cancel* on your overlay if you want.)

To select an option you can press or scan to the number of the key you want in the Pathfinder's Toolbox, or you can press on or scan to the key area on the display itself. If you want to change the Activation Feedback from "Beep" to "light", you can activate the number "6" in the toolbox or you can activate the "Activation Feedback" key on the Pathfinder display.

Since there are more than 6 options to the FEEDBACK MENU, there is a *More Items* key in the menu (key #7). Activate "7" or press the *More Items* area on the display to see more options. Any *Yes* or *No* responses you may have when using a menu can always be answered by activating the numbers 8 or 9 on your overlay (the same as *OK* and *Cancel*).

Exit a menu by activating *OK*. (the #8 key in the Toolbox). The Pathfinder will save your changes and exit the menu.

To exit a menu without saving any changes, activate "Cancel". Do this before you activate the "OK" key. Once the "OK" key is activated any changes are saved.

Using the Arrow Keys to Navigate in Some Menus

Some menus will have screens with categories of icons on them, or Activity keys.



When you see screens like these, you can use the arrow keys on your overlay to navigate through the smaller squares. When you activate an arrow a square is highlighted. Activate it again and the next square in the direction of the arrow is highlighted. To select a square, press the key directly or activate the OK key on your overlay.

You can also use the up and down arrow keys to navigate through the options within menus. For example, if you are setting the date and you go past the number you want, activate the down arrow once then activate the option you want until you see the correct date in the box.

You Can Access Most Menus Without Removing the Keyguard and Overlay

There may be times when you would like to go to a menu without having to remove the keyguard and overlay and opening the Toolbox. You can do this by using the SYSTEM STATUS key on the top left front of the Pathfinder.



Activate the **SYSTEM STATUS** key, . Activate the **TOOLBOX** key on the display screen. Select the Toolbox menu you want. (Not all Toolbox tools are represented in the Status Screen Toolbox.)

Change Speaker and Earphone Volume

You can raise or lower the speaker volume in the Toolbox.

Remove the keyguard and overlay. Turn the Toolbox On by activating the TOOLS key on the front of the Pathfinder case. You can also activate the SYSTEM STATUS key and then select TOOLBOX.

Activate SPEAKER VOL. UP to raise the speaker volume.

Activate SPEAKER VOL. DOWN to lower the volume.

Press the display screen, or activate one of the SPEAK tools in the Toolbox to hear the volume changes.

➡ The FEEDBACK MENU in the Toolbox also contains the Volume Up and Volume Down functions for the Speaker and for an external Earphone.

➡ See **Appendix H** for speaker and earphone specifications.

Set Up a User Area and Selection Technique

USER AREA MENU - BETA UNITY Sequenced-Hit - Area 2			1
1	Switch User Area	4	7
		Erase a User Area	More Items
2	Rename a User Area	5	8
		Overlay Menu	OK
3	Copy a User Area	6	9
		Computer Access Menu	Cancel

USER AREA MENU - BETA UNITY Sequenced-Hit - Area 2			2
1	Copy Selected Activities		7
			More Items
2	Copy All Activities		8
			OK
			9
			Cancel

User Area Menu

User Area Menu

The USER AREA Menu allows you to select an area to use and then set it up with a keyboard and overlay.

There are six areas available for use in the Pathfinder. This means that six different people can use the Pathfinder and each person can have a separate area. If only one person is going to use the Pathfinder, then that user can set up different areas for different requirements.

In the USER AREA menu you can:

Switch between User Areas

Rename an Area

Copy an Area

Erase an Area

Select an overlay and vocabulary

Select a Computer for computer emulation (IBM or Mac)

Copy an Activity from one User Area to another

➡ See page 163 for information about the *Copy Activity* option.

➡ **Important Note:** Vocabulary, Activities, Notebooks and the Selection Technique are specific to each User Area and cannot be accessed from other areas.

➡ **Important Note:** Imported Icons, the Speech Dictionary, the Word Prediction Dictionary and any Infrared (IR) Commands can be accessed from all User Areas.

➡ The *Copy All Activities* and *Copy Selected Activities* options are discussed on pages 163.

The next few pages discuss the USER AREA menu options.

USER AREA MENU - UNITY Sequenced-Hit - Area 2			1
1	Switch User Area	4 Erase a User Area	7 More Items
2	Rename a User Area	5 Overlay Menu	8 OK
3	Copy a User Area	6 Computer Access Menu	9 Cancel

Select the User Area to Switch to.		
1	PF Unity 128 Single Hit	4 PF Unity 32 Sequenced
2	PF Unity 128 Sequenced	5 Area 5
3	PF Unity 32 Single Hit	6 Area 6
		9 Cancel

Select a User Area

Select a User Area

There are 6 possible user areas in the Pathfinder. Each user area can be set up differently.

If different people are using the Pathfinder, each person can have his or her own area complete with its own settings, overlays, notebooks and keyboards. If just one person is using the Pathfinder, he or she might like to set it up so that one user area can be accessed with a switch. Another area might be set up for use with the keyboard. Different overlays might be used in different areas. One area might be set up for school or work and another for when the user is at home. A clinician may choose to set up a single user area in order to teach vocabulary.

The default set-up is: User Area 1 contains Unity 1-Hit vocabulary; Area 2 contains Unity Sequenced-Hit Vocabulary; Area 3 contains Unity 32 Single Hit; Area 4 contains Unity 32 Sequenced; the other areas are blank. Both the Unity 1 and Unity Sequenced-Hit overlays use the 128-location overlay.

To Select a User Area

1. Remove the keyguard and overlay. Turn the toolbox On.
2. Open the USER AREA menu.
3. Activate the *Switch User Areas* option.
4. Select the area you want to use.
5. Activate OK to exit the menu, or select another option.

Choose a selection technique for your area from the SELECTION TECHNIQUE menu.

 The selection technique you choose will be active with any overlays in your user area. If you set up a new user area, the selection technique settings for another area do not transfer to the new area. You must select the technique and the settings you want for the new area.

About Overlays

There are two Unity-for-Pathfinder overlays with pre-stored vocabulary for 128-locations: Unity One-Hit in User Area 1 and Unity Sequenced-Hit in User Area 2. The default is User Area 2 with Unity Sequenced-Hit vocabulary. Areas 3 and 4 contain Unity 32 overlays. There are also blank overlays for 32, 64 and 128 locations

If you are going to use a 64 or 32 location overlay, place the appropriate keyguard on the blank overlay. This will define the key locations for you.

When you are using a 128-location overlay you cannot access any messages you may have stored using the 64 or 32 location overlays. To get to your 64 or 32 location vocabulary you must switch to the correct overlay (and User Area if you stored in a different area). 64 and 32-location overlays do not have pre-stored vocabulary.

To Load an Overlay:

1. Remove the keyguard and current overlay. Turn the Toolbox On.
2. Open the USER AREA menu.
3. Switch to the User Area you want (option #1).
4. Activate the *Overlay Menu* option.
5. Select the overlay you want from the menu. Activate YES.
6. Activate OK twice.
7. From the packet of overlays that you received, select the overlay you want. Place the overlay you have selected on the keyboard. Place the keyguard that goes with the overlay on top of the overlay.

➡ **About the Zorro Pattern:** When you change from one overlay to another in the same User Area, you are asked to enter the Zorro pattern in order to clear out the User Area first. To perform the Zorro pattern activate the following keys in the order they are given: A1, P1, A8, P8, A8, P1, A1. The pattern traces a "Z" across the keyboard.

☞ If you want to save the vocabulary in the area before you replace it, perform a Memory Transfer and put the vocabulary on a computer disk. *See your MTI manual for information about memory transfers.*

☞ All overlays can be accessed by direct selection (keyboard or head-pointer) or by using a switch.

Rename a User Area

The original names of the different areas are User Area 1, User Area 2, User Area 3, etc. You can rename the areas you are using. If you are one of a number of people using the Pathfinder you might want to name your area after yourself. Or, if you are the only person using the Pathfinder, you may have set up different User Areas for different parts of your life. You may have an area for computer use and an area for school or work, for example. You can give each area a separate name to help you remember which area is used for what.

Once you assign a name to an area, that area will be referred to by its name until you change the name or erase the area.

To Rename a User Area:

1. Remove the keyguard and overlay and open the Toolbox.
2. Open the USER AREA menu.
3. Activate the ***Rename User Area*** option.
4. Select the area you want to rename.
5. Type in the new name. Activate OK.
6. Activate OK to exit the menu.

☞ **Note:** Up to 19 letters/spaces can be used for the name.

Copy a User Area

You can copy a User Area and use it in another area. For example, if you are moving from a simple vocabulary to a more complex one, you might want to copy your original User Area into a new area. Then you can add vocabulary to the new area as your vocabulary increases. You will still have your original vocabulary to use until you have mastered the new vocabulary.

To Copy an Area:

1. Remove the keyguard and overlay. Turn the Toolbox On.
2. Open the USER AREA menu.
3. Activate the *Copy a User Area* option.
4. Select the area you want to copy.
5. Select the area where you want the copy to go.
6. Activate YES if you are sure you want to copy, NO if you change your mind.
7. Activate OK to exit the menu.

Note:

If the area that you are *copying to* already has an overlay loaded into it, you will lose the vocabulary that was in the area you are copying to by performing the Zorro pattern (see page 51). If you want to save the vocabulary in the second area, perform a memory transfer and save that area onto a computer disk before you copy another area into it. See the *MTI Manual* for complete instructions.

Erase a User Area

There may come a time when you will want to erase a whole User Area. Perhaps the area is no longer relevant. You have graduated from that course, or changed jobs, or added so significantly to your vocabulary that the old area is useless. The *Erase a User Area* option allows you to get rid of outdated areas and begin a new area if you want to.

To Erase a User Area:

1. Remove the keyguard and overlay. Turn the Toolbox On.
2. Open the USER AREA menu.
3. Activate the Erase User Area option.
4. Select the option to erase.
5. Answer YES if you want to erase; NO if you change your mind.
6. Activate OK to exit the menu.

About the Computer Access Menu

➤ See pages 75-124 in the *Power Tools* manual for complete Computer Access information.

The Pathfinder has a Computer mode that is already set up for computer work. To go to the Computer Mode you must put the GO TO COMPUTER tool on your overlay or under an icon sequence. (See pages 175-182 for information about moving and storing tools.) If you are using the Sequenced-Hit Unity overlay, GO TO COMPUTER is already stored under the icon sequence: TOOLS+COMPUTER. (The TOOLS icon has the picture of a Toolbox on it). When you are in Computer mode, the LED next to the Computer icon on the front of the case will light.

To access a computer from the Pathfinder, you must have a keyboard emulator and you must select the kind of computer you are using (IBM or Mac). You must also have the Computer Access overlay on your Pathfinder keyboard. Just place it over your regular overlay.

To Select a Computer from the USER AREA Menu:

1. Open the USER AREA menu and activate the *Computer Access Menu* (option #6)
2. If you are using an IBM computer, set the *Host Computer* option to IBM and click OK.
3. Macintosh computers are a little trickier. If you are using a Mac with monitor manufactured before March of 1993, set the *Host Computer* option to *Older Mac*. If you are using a Mac with a monitor manufactured after March 1993, set the *Host Computer* option to *Newer Mac*.
4. If you discover that your mouse cursor is not doing what you want it to do (e.g., instead of opening the File menu it is opening the Edit menu), return to the *Computer Access Menu* and change your current MAC setting to the other one. That is, if you originally chose *Newer Mac*, and your mouse did not work correctly, change the setting to *Older Mac* or vice-versa.

➤ When you are in Computer mode, you cannot retrieve sequences to speak messages.

Pathfinder Selection Techniques

Keyboard

Infrared Headpointer

Single Switch (1-Switch)

Dual Switch (2-Switch)

Joystick

Selection Techniques

The Selection Techniques

Selection Technique refers to the method you use to control the Pathfinder.

Direct Selection Techniques

Keyboard: The keyboard consists of 128 keys on the front of the Pathfinder. Press on the keys to activate them. *The keyboard always remains active, no matter which selection technique you choose.* This means someone can use the keyboard no matter what other selection technique may have been selected.

Headpointer: Pathfinder supports infrared headpointing. Select *Headpointer* from the SELECTION TECHNIQUE menu. Point the headpointer at a key. When the activation delay has passed, the key is activated. You can also plug a switch into either switch port (A or B) and use the switch to activate the key you have pointed to.


Scanning Techniques:

1-Switch: Use a single switch to scan to any key on the keyboard. Plug the switch into one of the switch ports on the top of the Pathfinder. You can scan by row/column or column/row.

Quarter Row/Column Scanning: This is a single switch scanning technique. The Pathfinder keyboard is divided into four quarters. You scan a particular quarter for a key you want.

Dual or 2-Switch: Use a dual switch or two single switches to scan left, right, up and down to any key on the keyboard. Plug a dual switch or two singles into the switch ports on the top of the Pathfinder.

Joystick: Use a joystick to scan up, down, left, right and diagonally. Plug a joystick into the Multi-Switch on the top of the Pathfinder.

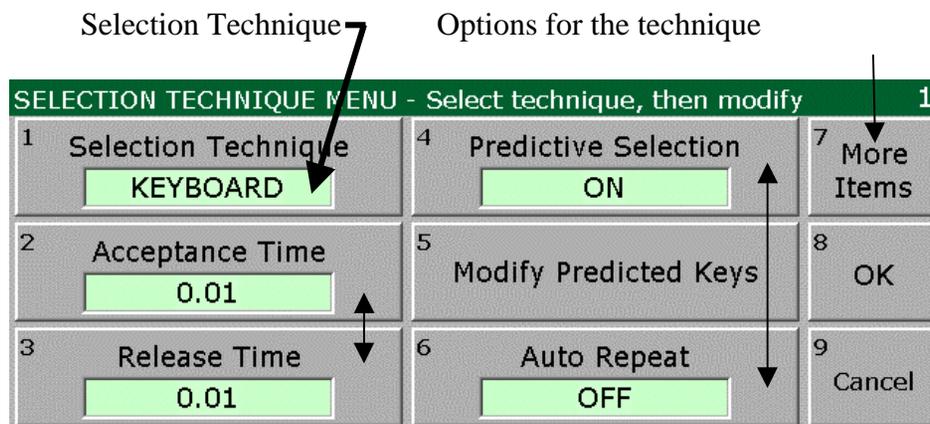

The SELECTION TECHNIQUE Menu

Using the SELECTION TECHNIQUE Menu

Summary:

Remove the keyguard and overlay. Turn the Toolbox On. Press the SELECTION TECHNIQUE key. Press the first option in the menu (#1) to step through the Pathfinder selection techniques. Each technique has its own optional settings.

➤ Before you choose a selection technique:

Make sure you have selected the User Area and overlay you want to set up for a particular selection technique.

If you are going to use a switch, headpointer or joystick, plug it into the correct jack on the top of the Pathfinder before you turn the Pathfinder on.

➤ The selection technique you choose to use with your overlay can also be used in the Toolbox. If you are scanning, you can scan in the Toolbox. ➤ Remember, the keyboard is always active and can be used, too.

Choose a Selection Technique

Remove the keyguard and overlay. Press the TOOLS button on the front of the Pathfinder case. The light on the TOOLS key will light when the Toolbox is on.

1. Press the SELECTION TECHNIQUE key in the Toolbox.
2. Press the "Selection Technique" option on the display screen or the #1 key on the Pathfinder to step through the different selection techniques. When you see the selection technique you want to use, stop pressing the key.
3. Select the set-up options you want to use.
4. When you have finished setting up the selection technique, press OK to save your changes and exit the menu. If you make a mistake or decide not to save changes, press "Cancel" to exit the menu.

SELECTION TECHNIQUE MENU - Select technique, then modify 1		
1 Selection Technique <input type="text" value="KEYBOARD"/>	4 Predictive Selection <input type="text" value="OFF"/>	7 More Items
2 Acceptance Time <input type="text" value="0.01"/>	5 Predictive Selection Maps	8 OK
3 Release Time <input type="text" value="0.01"/>	6 Auto Repeat <input type="text" value="OFF"/>	9 Cancel

SELECTION TECHNIQUE MENU - Select technique, then modify 2		
1 Activate Key When <input type="text" value="Pushed"/>		7 More Items
		8 OK
		9 Cancel

Keyboard Only

Keyboard Only


Summary:

If you select "Keyboard" for your selection technique you must be able to press the keys with your finger or some kind of "point and press" device. You may also press on the display screen to select appropriate options if that is easier for you. The keyboard is always available to use.



The keyboard and the headpointer are the quickest selection techniques. You make direct selections on the keyboard by pressing on the keys you want with your finger (or with the eraser end of a pencil or a "point and press" device). You can use the keyboard with any overlay.

The keyboard is the default selection technique for the Pathfinder. This means that when you receive your Pathfinder you can use it right away. Just turn it on and then press the keys on the keyboard.

The options for the keyboard selection technique are explained in the options section of this chapter on pages 61 and 75.

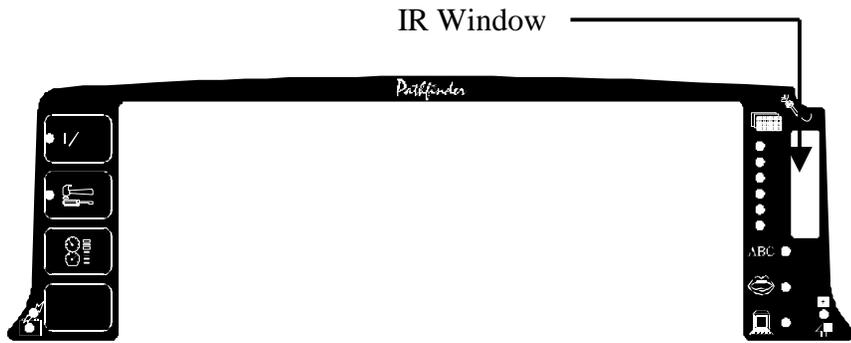

Headpointing Sensor Mounted on the Headset with the Body Pack

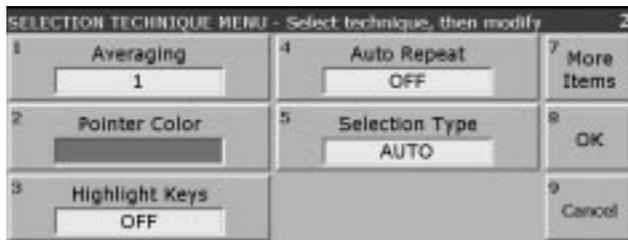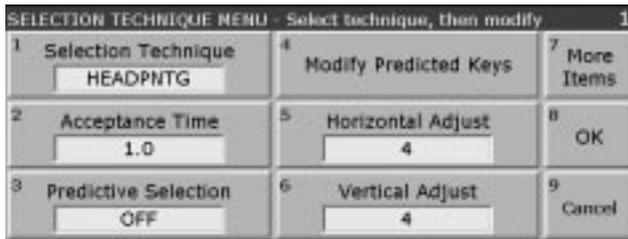

USER AREAS AND SELECTION TECHNIQUES

Pathfinder Manual: Basics

Infrared Headpointing

➡ **Always use the battery charger that came with your headpointer.**

Any other charger may damage your headpointing batteries. **Do not** use the Pathfinder charger to charge the headpointer body pack.


The Pathfinder infrared headpointing system consists of the infrared window on the top right front of the Pathfinder case, an infrared sensor (headpointer), a body pack and the battery charger for the body pack. The body pack contains the batteries that make the sensor work. (If the batteries run low, charge them by plugging the battery charger into the body pack.)

If you have been using an optical headpointer, you will find that the infrared headpointer works differently. You **do not** point the sensor at the keys on the keyboard to activate them. Instead, **you position the sensor so that it points at the infrared window on the front of the Pathfinder case.** (You can still look at the keyboard and move your head to find the keys you want, as usual.) The sensor must be able to "see" the IR window in order for the system to work correctly.

You can mount the sensor wherever it is most convenient for you. If you intend to mount the sensor on the PRC headset, headband or your own eyeglasses, ➡ please read the *Vanguard/Pathfinder Headpointer Mounting Instructions* that came with your headpointer.

➡ **If you want to use a switch** to activate selected keys, plug it into either switch port A or B on the top of the Pathfinder. Select **Manual** for the **Selection Type** option under **Headpointing** in the SELECTION TECHNIQUE menu.

The IR headpointing options in the SELECTION TECHNIQUE menu are discussed on pages 63-65, and 77-84.

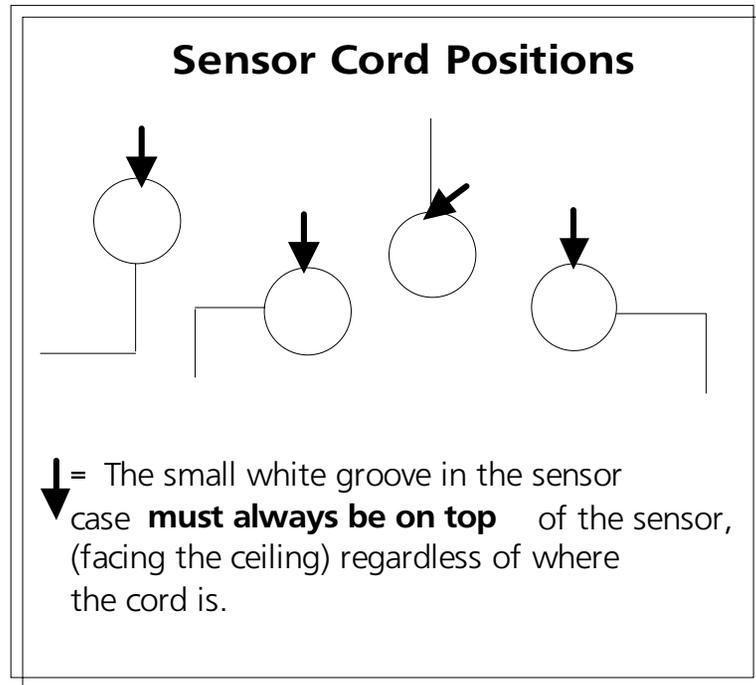

Position the Sensor Cord at 90-degree Angles

More About Your IR Headpointer

When the Body Pack Powers Down

The Body Pack will power down ("go to sleep") for two different reasons:

- If there is no sensor movement for 1 minute;
- If the sensor cannot see the Pathfinder for approximately 3 seconds.

"Wake up" the Body Pack by tilting your head up and down or by activating your switch.

About the Sensor Cord

If the sensor cord is in the user's way, you can rotate it at 90-degree angles around the sensor case. That is, you can move the cord from straight down to straight out to the left or right or straight up. See the graphic on the facing page. Instructions for rotating the cord are in the **Headpointer Mounting Instructions** that you received with the headpointer.

 **Warning! Do not pull on, or put undue stress on the sensor cord as you rotate it.**

The White Notch on the Sensor Case

Pick up the sensor and look at the round case. You can see a white notch on the sensor case. (This notch is **not** the cut-out area where the cord rests. It is a separate indentation on the main part of the sensor case.) No matter where you position the sensor cord, **this white notch must always be on top, facing upward (or toward the ceiling) when you mount the sensor.** If the white notch is not on the top of the sensor case, the sensor and headpointing system will not work.

 The sensor cord must be plugged into the **larger jack** on the top of the Body Pack before it will work.

SELECTION TECHNIQUE MENU - Select technique, then modify 1		
1 Selection Technique <input type="text" value="1 SWITCH"/>	4 Release Time <input type="text" value="0.01"/>	7 More Items
2 Scanning Speed <input type="text" value="1.0"/>	5 Predictive Selection <input type="text" value="ON"/>	8 OK
3 Acceptance Time <input type="text" value="0.01"/>	6 Modify Predicted Keys	9 Cancel

SELECTION TECHNIQUE MENU - Select technique, then modify 2		
1 Activation Delay <input type="text" value="OFF"/>	4 Scan Auto Restart <input type="text" value="ON"/>	7 More Items
2 Scan Mode <input type="text" value="MODE 1"/>	5 Auto Repeat <input type="text" value="OFF"/>	8 OK
3 Scan Type <input type="text" value="AUTO"/>	6 Scan Direction <input type="text" value="ROW/COL"/>	9 Cancel

SELECTION TECHNIQUE MENU - Select technique, then modify 3		
1 Number of Rescans <input type="text" value="3"/>	4 On-Screen Scanning <input type="text" value="ON"/>	7 More Items
2 Use Screen as Switch <input type="text" value="OFF"/>	5 Linear Scan Cutoff <input type="text" value="OFF"/>	8 OK
3 Quarter Row/Column <input type="text" value="OFF"/>		9 Cancel

1-Switch

1-Switch

Summary:

Plug a single switch or a dual switch (use only the SELECT side of a dual switch) into the Switch A or Switch B port on the top of the Pathfinder. The default scan is 'automatic'. The default scanning rate is '1 second'.

You can scan with a single switch, by using the SELECT side of a dual switch, or by using the display screen as a switch. You can use 1-switch scanning with the 128, 64 and 32-location overlays.

Plug the switch into the Switch A or Switch B port on the top of the Pathfinder.

Remove the keyguard and overlay from the Pathfinder and activate the TOOLS button on the front of the case.

Activate the SELECTION TECHNIQUE menu.

Press the "Selection Technique" key until you see "1-switch". The menu will then display the 1-switch scanning options. Pages 75, 89-100 describe these options.

Select the options you want. Activate OK to exit the menu when you are finished.

Note:

You will use your selection technique with any overlays in your area, (that is, whether you are using a 128, 64 or 32-location overlay).

SELECTION TECHNIQUE MENU - Select technique, then modify 1		
1 Selection Technique <input type="text" value="2 SWITCH"/>	4 Release Time <input type="text" value="0.01"/>	7 More Items
2 Scanning Speed <input type="text" value="1.0"/>	5 Predictive Selection <input type="text" value="ON"/>	8 OK
3 Acceptance Time <input type="text" value="0.01"/>	6 Modify Predicted Keys	9 Cancel

SELECTION TECHNIQUE MENU - Select technique, then modify 2		
1 Activation Delay <input type="text" value="OFF"/>	4 Scan Auto Restart <input type="text" value="ON"/>	7 More Items
2 Scan Mode <input type="text" value="MODE 1"/>	5 Auto Repeat <input type="text" value="OFF"/>	8 OK
3 Scan Type <input type="text" value="AUTO"/>	6 Scan Direction <input type="text" value="ROW/COL"/>	9 Cancel

SELECTION TECHNIQUE MENU - Select technique, then modify 3		
1 Number of Rescans <input type="text" value="3"/>	4 On-Screen Scanning <input type="text" value="ON"/>	7 More Items
2 Selection Type <input type="text" value="AUTO"/>	5 Linear Scan Cutoff <input type="text" value="OFF"/>	8 OK
3 Quarter Row/Column <input type="text" value="OFF"/>	6 Scan Start Position <input type="text" value="CENTER"/>	9 Cancel

2-Switch

2-Switch

Summary:

Use two single switches or a dual switch to scan up, down, left and right, or to do 2-switch row/column scanning. Plug the switch or switches into the Switch A and Switch B ports on the top of the Pathfinder. After each activation, the scan returns to the center of the keyboard.

With 2-switch scanning you can use either two single switches or one dual switch. If you are using two single switches, plug one switch into the Switch A port on the top of the Pathfinder. Plug the other switch into the Switch B port. *If you are using a Dual Switch, plug it into the Switch A port.*

Remove the keyguard and overlay.

Press the TOOLS key on the front of the Pathfinder case to turn the Toolbox on.

Activate the SELECTION TECHNIQUE menu.

Activate the "Selection Technique" option until you see "2-Switch".

Once you select 2-switch, you will see the 2-switch options. These options are described on pages 87-88 and 89-100. When you have finished with the SELECTION TECHNIQUE menu, activate the OK key to exit the menu.

Note:

You will use your selection technique with any overlays in your area, (that is, whether you are using a 128, 64 or 32- location overlay).

SELECTION TECHNIQUE MENU - Select technique, then modify

1 Selection Technique <input type="text" value="JOYSTICK"/>	4 Release Time <input type="text" value="0.01"/>	7 More Items
2 Scanning Speed <input type="text" value="1.0"/>	5 Predictive Selection <input type="text" value="ON"/>	8 OK
3 Acceptance Time <input type="text" value="0.01"/>	6 Modify Predicted Keys	9 Cancel

SELECTION TECHNIQUE MENU - Select technique, then modify 2

1 Activation Delay <input type="text" value="OFF"/>	4 Number of Rescans <input type="text" value="3"/>	7 More Items
2 Scan Type <input type="text" value="AUTO"/>	5 Selection Type <input type="text" value="AUTO"/>	8 OK
3 Auto Repeat <input type="text" value="OFF"/>	6 On-Screen Scanning <input type="text" value="ON"/>	9 Cancel

SELECTION TECHNIQUE MENU - Select technique, then modify 3

1 Scan Start Position <input type="text" value="CENTER"/>	7 More Items
	8 OK
	9 Cancel

Joystick

Joystick

Summary:

A joystick allows you to scan up, down, left, right and diagonally. Plug the joystick into the Multi-Switch port on the top of the Pathfinder.

You can push a joystick in any direction to scan up, down, left, right and diagonally. You can also use an Arm Slot control instead of a joystick to do the same things. An Arm Slot control allows joystick-type scanning for people who do not have enough fine motor skills to use a joystick.

Remove the keyguard and overlay.

Press the TOOLS key on the front of the Pathfinder case to turn the Toolbox on.

Activate the SELECTION TECHNIQUE menu until you see **Joystick**.

Once you select Joystick, you will see the Joystick options. These options are described on pages 71, 75, 90-100. When you have finished with the SELECTION TECHNIQUE menu, activate the OK key to exit the menu.


Move the joystick to scan to the location you want. The Pathfinder will automatically select the location you have chosen after the acceptance time has passed or you can activate a separate switch to select the location (*see Note below*).

Important Note:

If you choose "Manual" selection type from the menu, you must plug a separate switch into the Switch B connector on the top of the Pathfinder. Use the SELECT side of a dual switch or use a single switch). When you scan to the location you want, activate the separate switch to select your location.

 When you are scanning in the **Activity Row**, all joystick activations move to the right. You cannot go left.

SELECTION TECHNIQUE MENU - Select technique, then modify 3			
1	Number of Rescans <input type="text" value="3"/>	4 On-Screen Scanning <input type="text" value="ON"/>	7 More Items
2	Use Screen as Switch <input type="text" value="OFF"/>	5 Linear Scan Cutoff <input type="text" value="OFF"/>	8 OK
3	Quarter Row/Column <input type="text" value="OFF"/>		9 Cancel

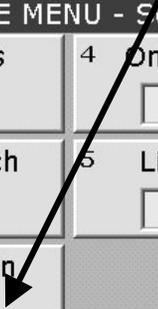

Quarter Row/Column Scan
OFF/ON

Quarter Row/Column Scan – Off/On

Quarter row/column scanning is a selection technique that can be used with 1-switch or 2-switch scanning. Turn this option On, in order to use it.

With Quarter Row/Column scanning the Pathfinder keyboard is divided into four quarters. Activate the switch to choose the quarter that contains the key you want. Once the quarter has been selected, the Pathfinder scans the quarter according to the other settings you have made, i.e., 1-Switch Row/Column Scan or 2-Switch (scan mode 2).

Quarter row/column scanning is somewhat faster than single row/column scanning but it does require more switch activations.

➡ Quarter Row/Column Scan **does not work** in 2-Switch Scan Mode 1 or in Linear Scan Directions.

SELECTION TECHNIQUE MENU - Select technique, then modify			2
1	Activation Delay <input type="button" value="OFF"/>	4	Scan Auto Restart <input type="button" value="ON"/>
2	Scan Mode <input type="button" value="MODE 1"/>	5	Auto Repeat <input type="button" value="OFF"/>
3	Scan Type <input type="button" value="AUTO"/>	6	Scan Direction <input type="button" value="ROW/COL"/>
		7	More Items
		8	OK
		9	Cancel

Some Scanning Options

Other Options in the SELECTION TECHNIQUE Menu

Below and on the next few pages you will find explanations of the options in the SELECTION TECHNIQUE menu.

Scan Type: Auto, Hold: *Auto* means the Pathfinder automatically scans the keyboard. *Hold* means you must press and hold your switch to scan. Step is an option if a Linear Scan Direction is chosen. Activate your switch each time you want to move to a new location.

Selection Type (2-Switch, Joystick, Headpointing): Auto/Manual: *Auto* means the Pathfinder will automatically select a location after the acceptance time has passed. *Manual* means that you must activate your switch to select a location.

Scan Auto Restart: On/Off: When this is *Off*, you must activate your switch to start a new scan. When this is *On*, the Pathfinder automatically starts a new scan after a message has been selected. If you are visually impaired, we suggest you set this to Off. This way the scan will always start at the same place (upper left corner of the keyboard).

Number of Rescans: 0-10: This allows you to select the number of times you want the Pathfinder to automatically re-scan the keyboard. You can select a number from 0-10. The default is 3 times.

On-Screen Scanning: Last/First/Off: **Last** scans core then display; **First** scans display then core; **Off** does not scan display.



Use Screen as Switch: this allows you to tap on the screen as you would a switch.

Activate Key When: (Keyboard) this allows you to tell the device to activate a key when you touch the key or when you release the key. The default is "Touched."

🔊 You can turn the Scanning Beep On or Off in the FEEDBACK MENU.

HEADPOINTING

Averaging

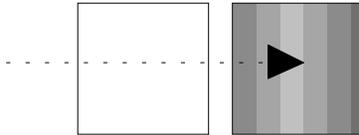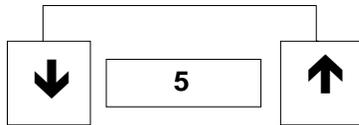

The **HIGHER** the Averaging number the farther you must move into the next key before it is activated.

↻ The default is "2"

Headpointing: Averaging


Summary:

Cut down on accidental key activations by selecting a number from 1-5. 2 is the default.

Averaging is a way of helping you cut down on accidental key activations caused by involuntary head movements. The number you select tells the pathfinder how far you must move out of one key location before a second key is activated. This generally means that when you cross the centerline between one key and the next key, the LED on the first key will go out and the LED on the second key will light.

If you select an averaging number **higher than "1"** you are telling the headpointer signal to "lag behind" any head movements you make. The **higher** the number you choose, the further into the key's area you must move before the LED on the second key will light.

Averaging numbers between 2 and 5 are recommended. The default is 2.

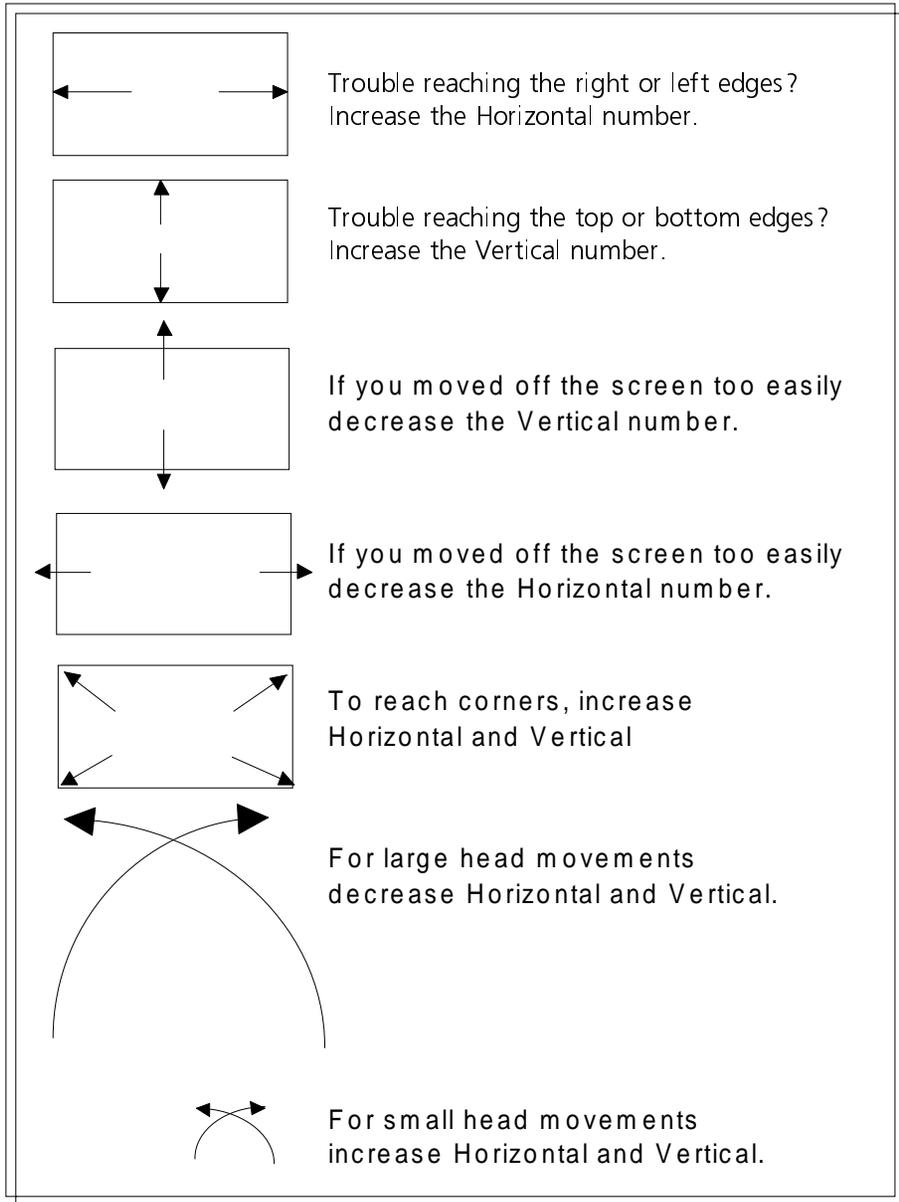

Horizontal Adjustment; Vertical Adjustment

IR Headpointing: Horizontal Adjustment and Vertical Adjustment

Summary:

If you have trouble selecting keys on the edges of the keyboard, or positioning the pointer in the corners of the keyboard, use the adjustments in this option.

Once your headpointer is connected and you are positioned properly, move your head and try to place the cursor in each corner of the screen—top right, top left, bottom right and bottom left. Move your head from left to right across the screen and from the top to the bottom of the screen. Did all the keys on the keyboard become highlighted?

If you had trouble reaching the right or left edges, **increase** the *Horizontal* adjustment number. If you had trouble reaching the top or bottom edges, **increase** the *Vertical* adjustment number.

If just the slightest head movement to the left or right moved you off the edges of the screen, **decrease** the *Horizontal* number. If the slightest head movement up or down moved you off the top or bottom edges, **decrease** the *Vertical* adjustment number.

If you had trouble positioning the cursor in the corners of the screen, **increase both** the *Horizontal* and *Vertical* numbers. If you make large head movements you may lose the cursor easily. Try **decreasing** the *Horizontal* and/or the *Vertical* numbers. If you make small head movements, try **increasing** the *Horizontal* and *Vertical* numbers.

Headpointing adjustments that you set when you are fresh and alert may not work as well for you later in the day when you become tired. If this happens, you can make new vertical and horizontal adjustments. You can also save the settings under an icon sequence on your overlay so that different settings are always available to you. See page 283.

SELECTION TECHNIQUE MENU - Select technique, then modify			2
1	Averaging <input type="text" value="2"/>	4	Highlight Keys <input type="text" value="OFF"/>
2	Select Window <input type="text" value="Right"/>	5	Auto Repeat <input type="text" value="OFF"/>
3	Pointer Color <input type="text" value="Red"/>		7 More Items
			8 OK
			9 Cancel

Select a Headpointer Color

IR Headpointing: Color Pointer

Summary:

You can change the color of the headpointer if you do not like the default color.

When you use Headpointing you will see a pointer in the shape of a hand on your display. The default color for the hand is gray. If this is hard for you to see, or if you would just like a different color, select the Pointer color option in the menu. Activate the option to see the different colors. When you see the color you want to use, activate the OK key.

IR Headpointing: Battery Information

When the headpointing battery is low, a red LED on top of the Body Pack lights to warn you to plug in the headpointing battery charger.

When you first receive your sensor and Body Pack, charge the batteries overnight.

Plug the small end of the charger that came with the headpointer into the battery charger connector on the side of the Body Pack. Plug the other end of the charger into a working A/C wall outlet. When the charger is plugged in, the red LED on the Body Pack lights until you unplug the charger. The next morning, unplug the charger from the Body Pack. Use the headpointer until the **red LED on the Body Pack** comes on, indicating that the batteries are low. This should take **3-4 days**, depending on how much you use the unit. When the light comes on, recharge the unit. If it took **4 days** before the red light came on, **re-charge the batteries once every 3 days**. If the low battery light came on after **3 days**, **re-charge the battery once every 2 days**. This re-charging schedule means you should be charging the batteries *before* the red light comes on.

Once a month, allow the batteries to go completely dead.

This means you should use the unit until it no longer works. When the batteries are dead, plug in the charger and charge the batteries overnight (8-10 hours). If you use the headpointer while the battery charger is plugged in, **the batteries will not charge**. Plan **not** to use the unit while it is charging. If you follow these guidelines the batteries should last for at least 2 years.

Replacing the Batteries

If the batteries quit working completely and will not re-charge, you will have to replace them. **You cannot use alkaline or NiCad** batteries in the Body Pack. If you think you need to replace your batteries, call the PRC Service Department for specific battery information.

Connecting the HeadMouse to the Pathfinder

➡ Before you begin **make sure the Pathfinder and the HeadMouse are turned Off.**

1. Find the cable that came with your HeadMouse. Plug the male end of the cable into **Port A** on the back of the HeadMouse Interface Unit (large box with a series of connectors on the back). Plug the female end of the cable into the **Multi** port  (picture of joystick) on the Pathfinder.
2. Use the enclosed Dual-Lock to attach the Optical Sensor (smaller box with black circle on it) in between the two headpointing windows on the Pathfinder.
3. Plug the Optical Sensor power cord into the appropriate connector on the front of the Interface Box.
4. **Place the thin reflective dot** on your forehead, the bridge of your glasses or any other moveable location on your body where the Optical Sensor can see it.
5. Turn the HeadMouse On (on the back of the Interface Unit), then turn the Pathfinder On. ➡ *Do not turn the Pathfinder on until you have turned the HeadMouse on first.*
6. Select the **Tools** key  (top right corner of the case).
7. Choose the **Selection Technique** Menu.
8. Choose **Headpointing**.
9. Adjust the **Acceptance Time** so it is appropriate for the person using the HeadMouse.
10. Adjust the mouse speed by using the **Horizontal Adjust** and **Vertical Adjust** options.
11. Activate OK to exit the menu when you are finished adjusting the HeadMouse.

Connecting Madentec's Tracker 2000 to PRC's Pathfinder

➡ Use these instructions in conjunction with the instructions in the Tracker 2000 manual.

➡ **Position the Tracker 2000 Tracker unit on the Pathfinder.** The round ball at the front of Tracker 2000 contains a sensor. This sensor must be positioned so that it points toward the user.

1. Plug the large connector directly into the back of the Tracker 2000 unit.
2. Plug the PS/2-to-serial-mouse-adapter-cable (yellow band) into one of the small round connectors of the PS/2 mouse cable (orange band)
3. Plug the 9-pin serial connector into the **Multi** port  (picture of joystick) on the top of the Pathfinder.
4. Plug the AC power adapter into the back of the Tracker 2000 unit. Plug the other end into a working wall socket or power bar.

Turn on Headpointing in the Selection Technique Menu:

1. Select the **Status** key twice  (left of the display)
2. Choose **Selection Technique** Menu (Select Tech)
3. Choose Selection Technique option #1 until it says **Headpointing**.
4. Select **OK** to exit the menu

Place the reflective dot on your forehead, the bridge of your glasses or any other moveable location on your body.

- If you see a **green light**, Tracker 2000 has detected the dot and is in the best position.
- If you see a **red light**, you need to center yourself or position yourself closer to the Tracker.

1-Switch-Scan Mode 1

Summary:

This is a single switch scan mode that requires 3 switch hits to make a selection when set to Row/Column scanning.

The default scan is *Automatic* and the scanning rate is *1 second*. You can also choose to use the display screen as a single switch.

To Scan in 1-Switch Scan Mode 1

Activate your switch to start the scan. If you are using **Row/Column** scan, the top row on the keyboard lights, then the second row, then the third, etc. Each time a row lights, you hear a beep. When the row you want is lit, activate your switch or tap the display screen to select that row. The scan then moves from left to right across the columns in the row. When the key you want lights up, activate your switch or tap the screen to select that key. If a message is stored in the location you have selected, and the speech is turned on, the Pathfinder will speak the message, then return to scanning down the rows.

If you do not make any activations after starting a scan, the Pathfinder automatically scans the keyboard 3 times and then stops scanning. To begin another scan, activate your switch or tap the display.

 You can select the amount of rescans you want in the "Number of Rescans" option.

1-Switch-Scan Mode 2

Summary

Scan Mode 2 involves 4 switch activations instead of 3.

Scan Mode 2 is another single switch technique. This technique involves 4 switch activations instead of 3 as in Scan Mode 1.

To Scan

Activate your switch. The top row on the keyboard lights. The scan automatically continues down the rows. When it reaches the row you want, activate and release your switch. The scan stops at the row you have chosen. Press the switch again and the scan proceeds automatically across the row. When the scan reaches the key you want, activate your switch to select the location.

➡ If you prefer Column/Row scan, columns will be chosen first, then the row.

2-Switch-Scan Mode 1

Summary:

Activate and hold switch A to start scanning the middle row. Activate switch A again to reverse scan. Activate Switch B to scan columns. Activate switch B again to reverse the scan. Tap the switch to select a location.

Scan Mode 1 is the default scan mode for 2 switches.

To Scan

2 Single Switches (Auto): "Auto" is the default. Each switch alternates the direction of the scan between up/down and left/right. The first activation of switch A starts a scan of the *middle row*. Press switch A again and the scan reverses. The first activation of switch B starts the scan across the columns. Activate switch B again and the scan reverses. Activate and *hold* the switch to start a scan. *Tap* a switch with a short hit to select a location that is lit.

Dual Switch: The principle is exactly the same as above. The SELECT side of your switch scans left/right and the ON/OFF side scans up/down.

You can choose **Auto** or **Hold** with 2-switch scanning.

Column/Row is not available with 2-switch Scan Mode 1.

➡ *In the Activity Row both switches will only move from left to right. You cannot move from right to left.*

2-Switch-Scan Mode 2

Summary:

This is a 2-switch row/column technique.

Switch 1 starts the scan. Switch 2 selects row or column. Switch 1 scans row or column. Switch 2 selects location.

You can use 2-Switch Scan Mode 2 with 128, 64 and 32 locations.

Press switch A to start the scan. The lights in the top row or the first column light. Hold down switch A to move the scan to the next row or column. You can also press and release the switch each time you want to move the scan to a new row or column.

When you reach the row or column you want, press switch B. Then press and hold switch B to scan across or down. You can also press and release the switch until you reach the location you want. When you have reached the location you want, press switch B again.

With this technique, switch 1 always moves you somewhere and switch 2 always chooses your locations.

Stepping or pressing and holding the switch are the only scanning options with this technique.

Scan Direction

Row/Column: This scan starts by lighting the top row (across the keyboard) of keys. It moves down to the second row, then the third row, etc. When the row that contains the key you want is highlighted, activate your switch. The scan then highlights the separate columns in the row you selected. The scan moves from left to right across the columns. When the key you want is highlighted, activate your switch.

Column/Row: This scan begins by scanning the columns first instead of the rows. The scan moves from left to right across the columns. When the scan reaches the column that contains the key you want, activate your switch. The scan then highlights each key in that row. When the key you want is highlighted, activate your switch.

☞ Row/Column and Column/Row are options for most of the scanning selection techniques.

Linear Scan: This allows you to step through keys on your overlay moving from **left to right** across the rows or from **top to bottom** down each column. If you select *Step* you must activate your switch each time you want to move to a new location. If you select *Auto*, the keys are scanned automatically. Activate your switch to start a new scan.

Scanning Speed

Summary:

Scanning speed is the amount of time it takes the Pathfinder to move from one key to the next in a scan.

➡ Scanning Speed is used with all scanning techniques.

Scanning speed refers to the time it takes the Pathfinder to move from one key in a scan to the next key. As each key is scanned, its light comes on and the Pathfinder beeps. The default scanning speed is 1 second. This means that the Pathfinder will move to a new key every second. You may want to use a slower speed to start with, perhaps 4 seconds between beeps. This means it will take 4 seconds from the time one key lights and the Pathfinder beeps to the time the next key lights and the pathfinder beeps. If you are used to scanning you may want to try a very fast scanning speed, such as .1 seconds

To Change the Scanning Speed

1. Remove the keyguard and overlay and turn the Toolbox On by pressing the TOOLS key on the front of the Pathfinder case.
2. Activate the SELECTION TECHNIQUE menu.
3. Activate the "Scanning Speed" option until you see the speed you want to use.
4. Activate OK to set the speed and to exit the menu.

Acceptance Time

Summary:

Scanning or Direct Selection: how long you must maintain contact with a switch or key before the Pathfinder accepts a location and makes a key activation. **Headpointing:** how long you must dwell at the lighted key.

➤ Acceptance Time is used with all selection techniques.

Acceptance time is the amount of time it takes for the Pathfinder to accept a location and activate the key. If you are scanning or pressing keys, the acceptance time is the amount of time you must **maintain contact** with your switch or the key before the key is accepted and activated. The default acceptance time is .01 second for scanning. If you are using a headpointer, acceptance time is the amount of time it takes between when the LED for the key you pointed to lights up and when the Pathfinder beeps that the location has been activated.

The default Acceptance Time is .01 of a second for scanning and 1.0 seconds for headpointing. If one-tenth of a second (.01) seems too fast or one second (1.0) seems too slow for you, you can change the time in the SELECTION TECHNIQUE menu in the Toolbox.

To Change the Acceptance Time:

1. Remove the keyguard and overlay and turn the Toolbox On by pressing the TOOLS key on the front of the Pathfinder case.
2. Activate the SELECTION TECHNIQUE menu.
3. Activate the "Acceptance Time" option until you see the time you want to use.
4. Activate OK to set the acceptance time and close the menu.

➤ If you select a high acceptance time, for example, 5 seconds, you **must** maintain contact with the key or switch until the key activates. This means you must hold down the key or switch for 5 seconds.

➤ If you selected *Scan Type = Hold* the Acceptance Time starts when you **release** your switch.

Release Time

Summary:

Release Time is the amount of time it takes for you to release a key or switch and be able to activate another key or start a new scan.

➡ Release Time is used with all selection techniques except Headpointing.

Release Time is how long you tell the Pathfinder to wait **after you have released a key or switch**, before it allows you to make another activation with the same key or switch. Release time allows you to create a "time window" after you release a key or switch. During this "window" the Pathfinder will not recognize repeated key or switch activations.

For example, if you set your release time for 1 second, you must wait 1 second after you release your switch or key before you can continue to scan or activate a new switch. This means that if you release your switch and then accidentally activate it again before the release time has passed, the Pathfinder will not make any activations. When the release time has passed you can activate the same key or the same switch again.

➡ If a user is accidentally making double activations, increasing the *Release Time* should help.

Activation Delay

Summary:

When this is On, you can scan to a key, select it and then change your mind and begin to scan again before the key is activated. Activation Delay is automatically set to your scanning speed.

Activation Delay gives you a few seconds to change your mind when you are scanning and have selected a key. When Activation Delay is turned On, you can scan to a key, select it, and then if you change your mind, you can activate your switch again to keep the key from being activated.

When this option is On, the amount of the delay is always the same as your scanning speed. If your scanning speed is set to one second (1.0), then once you select a key you have one second in which to "de-select" the key by activating your switch again.

The default for Activation Delay is Off. Any key you select is activated with no delay.

Auto Repeat Rate

Summary:

Select a key, hold it down and it will repeat at the rate you have set.

➡ Auto Repeat works with all selection techniques.

Auto Repeat Rate lets you choose how fast a key will repeat itself if you hold it down.

If you are using a switch you just hold the switch down. If you are using a headpointer, you must continue to point at the key to get it to repeat automatically. If you are using a 2-switch or joystick and selected ***Selection type: Manual***, hold your extra switch down when you reach the key you want to repeat.

You can choose to have a key repeat once every .10 seconds to once every 2 seconds. You can also turn this option Off.

➡ *Auto Repeat does not work if you are using 2-switch with Scan Mode 1.*

Predictive Selection

Summary:

When Predictive Selection and Icon Prediction are both turned On, only those keys where messages are stored are active.

➡ Predictive Selection works with *Icon Prediction*. Both options must be turned On for Predictive Selection to work correctly.

When *Icon Prediction* (in the FEEDBACK menu) is On along with *Predictive Selection* (in the SELECTION TECHNIQUE menu), the Pathfinder lights only those icons where messages are stored. These lighted icons are the only keys on the keyboard that will produce a response. This means that if you are scanning, the Pathfinder scans only the lighted keys. If you are using direct selection (keyboard or headpointer) you will only be able to activate lighted keys. If the light belonging to a key is not lit, it will not work. ➡ If you are using an optical headpointer, any options for LED intensity will not work, including this one.

➡ **Note:** Icon Prediction will work whether or not Predictive Selection is On. See page 29 for an explanation of Icon Prediction.

You can create your own predicted keys. These are like mini-applications. With your own predicted keys you can "mask out" any locations that are empty or that contain tools that you don't want the Pathfinder to scan. This greatly reduces the amount of time it takes to scan the complete keyboard. See the next page for details.

SELECTION TECHNIQUE MENU - Select technique, then modify			1
1 Selection Technique KEYBOARD	4 Predictive Selection ON	7 More Items	
2 Acceptance Time 0.01	5 Modify Predicted Keys	8 OK	
3 Release Time 0.01	6 Auto Repeat OFF	9 Cancel	

Modify Predicted Keys Menu		
1 Before 1st icon hit	4 Spell One Word mode	
2 After 1st icon hit	5 Computer Emulation Mode	8 OK
3 Spell mode		9 Cancel

Modify Predicted Keys

Modify Predicted Keys

Summary:

Keys with dim lights are inactive. Keys with bright lights are active. **Icon Prediction** must be **On**. Set **Icon Prediction** (in the FEEDBACK menu) to **Low** to better see the LEDs.

➤ Modified predictive keys work with all scanning techniques.

➤ If you are using an optical headpointer, Predictive Selection LED intensity (low, medium or high) in the FEEDBACK menu will not work.

When you choose the **Modify Predicted Keys** option, you can tell the Pathfinder exactly which keys you want it to scan and when.

1. Before 1st icon hit

These are keys that you want to be lit before you ever activate an icon. Consider them as "general" keys. For example, when you are in Minspeak, the time before you activate a key that begins a message and the time after you finish speaking the message are the times when you are most likely to decide to go to Spell mode, or to open a notebook or menu. Because of this, you will probably want most of the tools and menu keys to be active.

2. After 1st icon hit

These keys should be more specific because now you are committed to speaking and won't need so many options on your keyboard. For example, once you have selected the first icon in a sequence, you won't want all the tools and menus to be active while you are looking for your second icon. You can now make all the tools **inactive**, except those that are necessary for you while you are in Minspeak, such as DELETE ICON or CLEAR DISPLAY.

Spell Mode

You can make **inactive** all the tools in Spell mode that you seldom use, like DELETE ICON and any spelling characters that you know you aren't going to use.

Spell One Word

Make **active** any keys you use when you are in Spell One Word mode.

3. Computer Mode

At this time, we cannot think of any reason why you would want to inactivate any computer mode keys. If you find a use for this option, let us know.

To make a key inactive, select it **once**. Its LED dims to indicate that it is inactive. Select it **again** to make it active. Its LED will get bright.

Press #1, *Select All*, to set all the keys to **bright**.

Press #4, *Select None*, to set all the tools to **dim**.

Scan Start Position

 This option is available only with 2-Switch and Joystick scanning.

Scan Start Position tells the Pathfinder where to start each scan. *Center* indicates that each new scan will start from the middle of your keyboard. *Top Left* indicates that each new scan will start from the top left corner of the keyboard. *Scan Start Position* is used only with 2-Switch or Joystick scanning because these are the only scanning techniques where you can move up, down, left, right or diagonally from the center of the keyboard.

If you select *Center* with 2-Switch scanning, the scan will begin in the center of the keyboard. If you activate your "down" switch, the scan will continue down from the center. If you activate your "right" switch, the scan will continue to the right of the center. With a Joystick, the scan will continue in whichever direction you activate your switch.

About Linear Scan Cutoff

 **Linear Scan Cutoff** can be used with all scanning techniques.

This option allows a scan to change from row/column to key-to-key scanning.

When you select *Linear Scan Cutoff*, you are telling the Pathfinder that you want it to perform your regular scanning technique until there are only a certain number of selections left. At that point, you want the Pathfinder to stop scanning rows and columns or quarters and to directly scan only the keys that are left.

To Use This Option Effectively:

1. Set *Icon Prediction* to *Low* (this will not work if you are using an optical headpointer).
2. Turn *Predictive Selection* to *On*.
3. Use the *Modify Predicted Keys* option to eliminate all the keys that you don't want to have scanned at any time.
4. Activate *Linear Scan Cutoff*.
5. Enter the cutoff number you want to use. The cutoff number should be **equal to** or **less than** the number of predicted locations that are left on your overlay when you want to switch from regular scanning to key-to-key scanning.

For example, if you select "5" as your cutoff number, the Pathfinder will use your regular scanning technique to scan the keys on your overlay until there are only five predicted choices left. At that point, it will stop the regular scan and begin to scan those five keys with no regard to columns, rows or quarters.

Using Word Prediction



Use Minspeak for your core vocabulary

(all the words that you use frequently and can access easily with icon sequences).



Use Word Prediction

**for longer words that you don't use often
but want to have available.**

Word prediction saves you the time and energy of spelling these words out.

How to Blend Minspeak and Word Prediction

About Minspeak and Word Prediction

Minspeak Benefits

People who use PRC communication aids with Minspeak are known to be spontaneous, independent communicators. Using a single Minspeak overlay where icons are preset and basically unchanging, Minspeak users become automatic at recalling their stored vocabulary. When they are familiar with their vocabulary they no longer have to think about their communication device. Instead they can concentrate on actual communication—forming grammatical sentences, speaking to friends and colleagues, delivering speeches. A Minspeak user can learn to automatically recall and combine icon sequences to produce a very large vocabulary indeed.

Word Prediction Benefits

Word Prediction systems help improve a person's ability to communicate by cutting down on the amount of time spent in spelling out words.

Dynamic Word Prediction Systems

A Dynamic Word Prediction system predicts words based on how often (frequency) or how recently (recency) you have used the word. This means that the system is always changing and updating the words that it predicts. This makes it very difficult for you to automatically select words; you have to look at the screen to see where the word you want is slotted. With words that are less than six letters in length, using a dynamic system is often slower than actually spelling out the word. With words longer than six letters, significant time can be saved by finding the predicted word in the window.

**"We went to the museum
to see the Impressionist exhibit."**

Core vocabulary stored in Minspeak

we
went
to
the
museum
to see

**Special vocabulary to select from
Word Prediction**

Impressionist
exhibit

- ➔ Use Minspeak for core words.
- ➔ Use Word Prediction (Spell) for words that aren't part of your everyday vocabulary.

The Benefits of Minspeak and Word Prediction Together

The Benefits of Minspeak and Word Prediction Together

Minspeak vocabulary represents your **core** vocabulary—verbs, adjectives, nouns, interjections, phrases, etc. These are words that make up most sentences and turn sentences into conversation. But lots of sentences also contain words that are not exactly "core" words. They are more like "fringe" vocabulary; words that help to make the sentences interesting or to clarify something. These fringe words often are not stored in your Minspeak vocabulary because they are words that you only use occasionally. To use them, you must spell them out which takes time, cuts down on your ability to communicate rapidly and is labor-intensive. This is where Word Prediction helps.


Suppose you want to say, "We went to the museum to see the Impressionist exhibition." You can certainly recall "We went to the" and "to see the" using Minspeak. You might have "museum" stored in your device, too, in which case you can also recall it using Minspeak. But you probably don't have "Impressionist" or "exhibition" stored.

With Word Prediction turned On, you can recall the parts of the message that are in your Minspeak vocabulary, then begin to spell "Imp" or "exh" and select the word you want from the list of words on the screen.

Word Prediction allows you to more quickly select words that are not part of your core vocabulary so that you don't have to take the time to spell them out.

Your best bet for spontaneous, interactive communication is to become so familiar with your Minspeak vocabulary that recalling it is automatic. When situations arise where you need specialized vocabulary that you don't want to store in your system, use Word Prediction to select the words rather than taking the time to spell them. In many situations, a combination of Minspeak and Word Prediction will keep you talking!

WORD PREDICTION MENU		
1 Word Prediction On/Off <input type="checkbox"/> ON	4 Set Predictor Keys	7 More Items
2 Predict in Minspeak <input type="checkbox"/> OFF	5 Predictor Keys <input type="checkbox"/> ON	8 OK
3 Recency/Frequency	6 Flash Predictor Keys <input type="checkbox"/> ON	9 Cancel

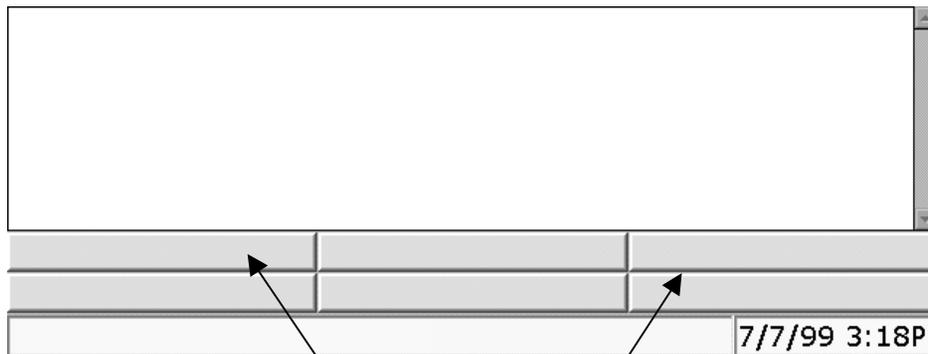

Word Prediction Keys

Turn Word Prediction On

Turn Word Prediction On or Off

The default for Word Prediction is Off.

To use it you must turn it On in the WORD PREDICTION menu.

➡ In order to see predicted words, you must be in SPELL mode.

When Word Prediction is On, the Word Prediction keys are visible in the bottom two rows on the Pathfinder display screen.

To Turn Word Prediction On or Off

1. Remove the keyguard and overlay. Press the TOOLS key on the front of the Pathfinder case.
2. Open the WORD PREDICTION menu.
3. Activate the first option in the menu until it says On
4. Activate OK.
5. Make sure you are in SPELL and begin to spell a word from the Toolbox overlay.

You can copy the tool, *Word Prediction On/Off*, to your communication overlay so that you do not have to go to the Toolbox to turn Word Prediction On or Off. See pages 175-182 for information about copying tools to your overlay.

imp

IMPRESSION	IMPACT	IMPORTANCE
IMPORTANT	IMPOSSIBLE	IMPRESSIVE

7/8/99 11:14A

Begin to spell a word. If the correct ending for the word you want is not listed, select the one that is most like the word you want.

impression

IMPRESSIONS	IMPRESSIONIST	IMPRESSIONISTIC
IMPRESSIONISTS		

7/8/99 11:15A

Now select the word that has the correct endings.

impressionist

7/8/99 11:15A

You are ready to begin a new word.

Spell and Select Words

Turn Word Prediction On in the Toolbox (WORD PREDICTION menu).
Make sure you are in SPELL mode.

Begin to spell a word. Notice that words begin to appear in the two rows of keys at the bottom of the Pathfinder display. These are the words that the Pathfinder is predicting from the letters you have spelled.

To select a word, press the word on the display that you want. The word is placed in the text area of the display and either the predicted words disappear, or you will see a list of words that are derived from the word you selected. You can select one of these words, or you can continue to spell and watch more words become visible in the Word Prediction area.

If you can't physically press on the display to select a word you can create predictor keys on the keyboard that correspond to the placement of the words in the Word Prediction area on the display. See p 113 for information on setting up predictor keys.

The words that are being predicted are from the Word Prediction dictionary. You may see a lot of unfamiliar words. Most of the Word Prediction dictionary words will be words that you do not use as often as you use your Minspeak vocabulary. Word Prediction words are meant to enhance your Minspeak vocabulary, not replace it.

➤ Important Note: If you are a scanner and want to scan the Word Prediction words on the screen, turn the *On-Screen Scanning* option *On* in the SELECTION TECHNIQUE menu.

WORD PREDICTION MENU		
1 Word Prediction On/Off <input type="checkbox"/> OFF	4 Set Predictor Keys <input type="checkbox"/> ON	7 More Items
2 Predict in Minspeak <input type="checkbox"/> OFF	5 Predictor Keys <input type="checkbox"/> ON	8 OK
3 Recency/Frequency	6 Flash Predictor Keys <input type="checkbox"/> ON	9 Cancel

WORD PREDICTION MENU		2
1 Word Prediction Dictionary		7 More Items
2 Computer command for "Backspace"		8 OK
3 Predicted Word Display		9 Cancel

WORD PREDICTION Menu

Set Up Word Prediction in the WORD PREDICTION Menu

Before you begin to use Word Prediction, you will want to set some options in the WORD PREDICTION menu so the Pathfinder knows what you want it to do.

Remove the keyguard and overlay. Press the TOOLS key on the front of the Pathfinder case to open the Toolbox. Open the WORD PREDICTION menu.

The next few pages explain the options in the WORD PREDICTION menu.

Predict in Minspeak

When this is On and you retrieve a Minspeak word that is also in the Word Prediction Dictionary, the word prediction system will show you all the derivations of the Minspeak word. For example, if you retrieved the Minspeak word, "Cool" (which is stored under KITCHEN THUMBS DOWN ADJ, the Word Prediction system would display words like, "cooled, cooling, coolest," etc.

 **Note:** You cannot use Activities if you turn *Predict in Minspeak* On.

Recency/Frequency

The Word Prediction system will keep track of and predict all the words that apply that you have used most recently. It also keeps track of the words you use most frequently. You can decide how many most recently used or frequently used words you would like to have predicted. If you select to see 6 most recently used words, all the Word Prediction keys will display the most recently used words that begin with the letters you have spelled. If you want to see 3 recently used words, three of the Word Prediction keys will display recently used words and the rest of the keys will display the three most frequently used words that begin with the letters you have spelled.

t		
TEMPER	TEST	TESTIMONY
TABLE	TAKE	TAKEN
		10/11/99 11:48

					F	G	H								
					1	2	3								
					4	5	6								

1 = TEMPER 2 = TEST 3 = TESTIMONY
 4 = TABLE 5 = TAKE 6 = TAKEN

 = Predictor Key

About Predictor Keys

You can select predicted words by pressing the word you want on the display screen, or, if you have trouble pressing on the screen, you can set up "predictor" keys on the Pathfinder that correspond to the Word Prediction keys on the display. This way you can press, or scan to, the appropriate key on the keyboard to select the word you want.

Predictor keys are only active when you are spelling a word. If you are in Minspeak mode, the predictor keys will act as icons and you can select them as part of a sequence. In Spell mode, if no words are currently in the Word Prediction list, the regular spelling assignment for each key is active.

Set your own Predictor Keys:

1. Turn the Toolbox On, open the WORD PREDICTION menu.
2. Activate *Set Predictor Keys*. The keys will say, *Unassigned*, or, if they are already set up, they will show the location where each is assigned.
3. Activate option #1.
4. Choose the key on your overlay that you want to use to pick the first word from the list of predicted words.
5. Choose option #2.
6. Choose the key on your overlay that you want to use to pick the second word from the list of predicted words.
7. Continue in this manner until you have assigned all the predictor keys.
8. When you have finished, activate OK.
9. Exit the menu.

Turn Predictor Keys Off and On

When Word Prediction is On, the Predictor Keys are On. When there is a word being predicted, even if you are in Minspeak mode, you will not be able to use the predictor key associated with the predicted word for anything other than selecting the word. If you want quick access to the predictor keys' regular key assignments, you will want to be able to turn Word Prediction Off easily.

We suggest that you store the tool, <**WORD PREDICTION MENU**>, #1 on your overlay. Either assign it to a key or store it in an icon sequence. This way you can turn the predictor keys Off when you need to and quickly turn them on again when you are back to typical spelling.

See pages 175-182 to learn how to store tools on your overlay.

Flash Predictor Keys

If you turn this On, the LEDs for the keys that you are using to select words from the Word Prediction display area will flash when there is a word predicted that you can select.

When this is Off, the LEDs for the predictor keys will light, but they will not flash when there is a word to be selected.

Word Prediction Dictionary

The Word Prediction Dictionary allows you to add, delete, edit and review words in the Word Prediction system. This dictionary is separate from your Speech Dictionary. The Word Prediction Dictionary contains the words that you want to have available for Word Prediction. It does not contain your Minspeak vocabulary. If you want to review your speech dictionary, you must go to the SPEECH menu.

The Word Prediction Dictionary can be accessed from all User Areas.

Add a Word to, or Delete a Word from, the Word Prediction Dictionary

To Add a Word:

Activate *Add a Word* and type in the root word of the word you want to add. (A root word has no endings to it, such as "es", "ing", etc.) For example, enter the word "abominate".

The display asks you to enter a word with similar word endings. This means you should enter another word that is used the same way in sentences and with the same kind of endings as the word you want to add has. For this example, type in "terminate".

The Pathfinder shows you a box containing the word you want to add with all its endings. If most of the endings are correct for the way you want to use the word, activate YES. The new word is entered into the Word Prediction Dictionary. If the endings aren't right for your word, activate NO, type in a new similar word and try again.

To Delete a Word

Activate *Delete a Word*. Type in the root word you want to delete. Activate OK. Answer YES if you want to delete the word. The root word and its derivations will be deleted. Answer NO if you change your mind. Nothing will be deleted.

Review Words in the Word Prediction Dictionary

You can review all the words in the Word Prediction Dictionary or you can see only those words that begin with specific letters. You can also delete or change word derivations.

Review Method

You can review by the *STEP* or *AUTO* method.

STEP means that you must activate your switch or touch a key or the display to move from one screen of words to the next.

AUTO means that the Pathfinder will automatically change screens for you after a pre-set amount of "looking" time has passed.

Review to Display

If you select this option, the words in the dictionary are shown on the Pathfinder display screen.

If you want to review words that begin with specific letters, type in the letters. The more letters you type in, the more specific the entries will be that the Pathfinder shows you. If you want to look at the whole dictionary, activate OK without typing anything.

 **Helpful Hint: Editing a Word and its Derivations:** You can edit a word by selecting the word on the display or by choosing the corresponding number from your keyboard. Then you can choose to change the endings for the word or to delete the word and its endings from the dictionary. If you want to keep the word but not have any endings associated with it, activate *Change Endings*, then activate OK without typing in anything. Activate YES when asked.

Review to Display and Speaker

This works the same way as *Review to Display*, except that the Pathfinder will also speak the entries on the screen.

Review to Output

This allows you to review the dictionary on a computer screen or print it out on a printer. You must be connected to a computer or printer and you must set the correct options in the OUTPUT menu for this option to work. ➤ See **Appendix C** for complete information.

About the Computer 'Backspace' Command

You will only use this option if you are using a computer emulation program that has a command for "backspace" that is different from the Pathfinder's command. The Pathfinder command is: ←backspace.

If the "backspace" command for your computer is different from the one that the Pathfinder uses, you must activate this option and type in the keyboard emulation command that your program uses. Activate OK when you have finished entering the command.

When Word Prediction is On and *Serial Output* is set to "Immediate" in the OUTPUT menu, any characters you type are sent to the Pathfinder display and to the computer monitor. When you select a word from the Word Prediction list, the characters that you spelled are deleted from the Pathfinder display and are replaced by the word you have selected. The characters you spelled must also be deleted from the computer monitor. The Pathfinder sends the "backspace" command for each character that is to be erased from your computer, so it is important that it knows the correct command.

For example, suppose you spell the word, "berry". In the Word Prediction list you see the word "berries" and you select it. The Pathfinder deletes the "y" in "berry" and adds "ies". On your screen you now see "berries."

➡ The computer "backspace" command must be used to delete one character at a time. The Pathfinder "Delete Character" and "Delete Word" keys will not work correctly.

Predicted Word Display

This option allows you to select a color for the Word Prediction keys and a color for the text on them that will help you to see them better. You can also choose the size of the text and whether the words are capitalized or not.

➤ If you choose "Large" text, the Pathfinder will use the large text for all words that fit easily on the key. If a word will not fit on the key in large text, the Pathfinder will make the text smaller and it may change from capital letters to lower case letters in order to make the word fit the key.

Background Color: This changes the color of the Word Prediction keys. Step through the color options and select the color you want.

Text Color

Step through color options for the text.

When you are finished, activate OK.

Storing Messages

Select type of message, then create it. Select OK when finished.

Message: <EMPTY>

Sequence: KNOW SUN

2	Type of Message Synthesized	5	Record Message	8	OK
3	Spell Text for Message	6	Listen to Recording	9	Cancel

Select type of message, then create it. Select OK when finished.

Message: <EMPTY>

Sequence: KNOW SUN

2	Type of Message Recorded	5	Record Message	8	OK
3	Spell Text for Message	6	Listen to Recording	9	Cancel

Synthesized and Recorded Speech

Store Digitized (Recorded) or Synthesized (Text-to-Speech) Messages

Summary:

Recorded = messages someone has stored by speaking into the Pathfinder's microphone.


Synthesized = messages that are typed into the Pathfinder which the Pathfinder then translates into its own speech.

The Pathfinder allows you to store messages in two different kinds of speech—digitized or synthesized. **Digitized** speech is speech that someone else has recorded into the Pathfinder, just like recording something into a tape recorder. **Synthesized** speech is speech that the Pathfinder produces from words that are typed into it. Each time you store a message you can select which kind of speech you want to use.

You will probably store most messages in synthesized speech. This is the best way to become familiar with Minspeak and to increase language abilities. But, it is fun to use recorded speech for things like clapping, laughter, singing and various animal or human sounds that don't sound natural when they are reproduced in synthesized speech.

The Pathfinder has enough digitized memory to record approximately 15 minutes of digitized speech.

➡ When you record messages, make sure the person doing the recording is a good representative of the person who uses the Pathfinder. If the Pathfinder user is an 8 year old girl, have an 8 year old girl do the recording. If the user is a middle-aged man, ask a middle-aged man to record the messages. We recommend that you do not use family members or close friends to do the recording. Remember, the Pathfinder is the *voice* of the person who uses it. You do not want to confuse it with someone else's voice.

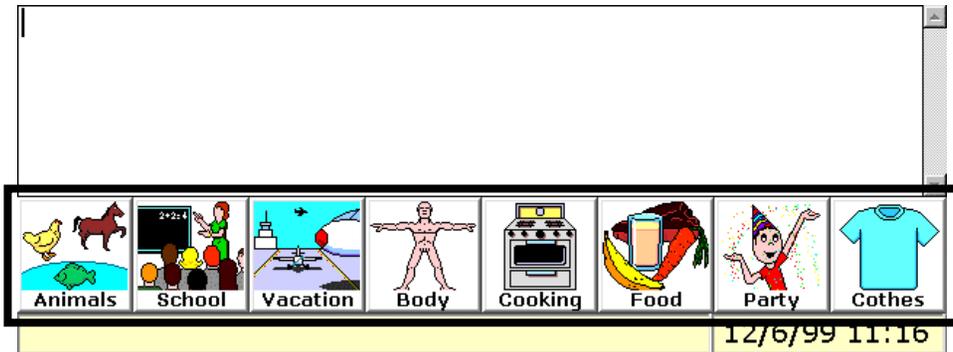

Activity Row

(The Activities in your Pathfinder will be different)

About Core Vocabulary and Activity Vocabulary

Summary:

Core Vocabulary is vocabulary that you use everyday.

Activity vocabulary can be vocabulary that you use in special circumstances.

The Unity overlays have two kinds of vocabularies stored in them, Core vocabulary and Activity vocabulary.

Core vocabulary is made up of words that you use everyday. Besides frequently used nouns and verbs, it also contains words like, "it", "and", "about", "the": all the words that help make complete sentences.

Core vocabulary is stored under icon sequences on your overlay. It usually takes one to three activations to retrieve a word or a message.

The **Activity Row** can contain vocabulary that you use in special circumstances. For example, vocabulary that you use mostly at your job or in a class at school, or when ordering in a restaurant could be stored under different activities in the Activity Row.

The Activity Row shows up on the bottom of your display screen when you have turned **Activities On** in the ACTIVITY menu.

The Unity overlays use the Activity Row *dynamically*. Dynamic icons are not visible in an Activity row. When you make an activity dynamic you can activate a key on your overlay and then activate an Activity Row key, speak a word or message and return to your overlay to continue speaking. You might have nouns or picture producing words stored in your Activity Row for this kind of use. For example, activate the APPLE icon on your overlay. You will see icons on your display that you can choose from to produce a message. See pages 157-159 to learn about Dynamic Activities.

About the SPELL/MINSPEAK, GO TO SPELL, GO TO MINSPEAK and GO TO COMPUTER Tools

Summary:

MINSPEAK allows you to retrieve messages by activating icons.

SPELL allows you to type messages into the Pathfinder.

When you are selecting icons to store messages or to speak the messages you have stored, you are in MINSPEAK mode. When you are putting text on the display, either by typing or scanning, you are in SPELL mode.

When you are connected to and using a computer you are in COMPUTER mode.

The SPELL/MINSPEAK tool in the Toolbox and on your overlay allows you to toggle between Spell mode and Minspeak mode. To use a computer, you must store or assign the tool on your overlay. (On the full Unity vocabulary, GO TO COMPUTER is already stored under the icon sequence: TOOLS+COMPUTER+COMPUTER.) You can also store or assign the tools, GO TO SPELL and GO TO MINSPEAK. (See pages 173-182 to learn about using INSERT TOOL and storing or assigning tools on your overlay.)

On the right front of the Pathfinder case there are LEDs next to the icons for SPELL, MINSPEAK and COMPUTER. When you are in one of these three modes, the icon next to the active mode will light. For example, if you are trying to spell, but every time you activate a key an icon appears on the display, check the front of the case. If the LED next to the Minspeak icon is lit, you will know that you are in Minspeak and not Spell mode. Activate the GO TO SPELL key until the light next to the ABC icon is lit.

 When *Auto Spell* is On (in the FEEDBACK menu), the Pathfinder automatically switches you from Minspeak to Spelling. You don't have to activate the GO TO SPELL key. See page 198 for more information about Auto Spell.

Before You Begin to Store, PLAN What You Want to Say

➡ If you want to store Activity Row Vocabulary, see pages 151-163.

You can store words, phrases and sentences in the Pathfinder's memory. You can store messages under an icon sequence of up to 8 icons, or you can store them under a single key on your overlay.

Before you begin to store, plan out the messages and the icon sequences you want to use. Do you want to store Core vocabulary or Activity vocabulary? Planning helps you to store the message faster, and, if you are recording, keeps you from holding the RECORD option down while you think about what you want to say. The longer you hold the RECORD option down, the more memory space you waste.

If you are using one of the Unity overlays, you might want to do a review of the overlay's vocabulary first to make sure that the message you are thinking of storing isn't already stored there. You can also use your Unity Vocabulary Sort to review the Unity vocabulary. (See page 165 for information about the REVIEW VOCABULARY menu.)

When you record a message, direct your voice toward the microphone on the top right corner on the front of the Pathfinder.

➡ Scanners can record digitized messages by holding their switch down on the Record Message button.

Select type of message, then create it. Select OK when finished.

Message:

Sequence:

2	Type of Message <input type="text" value="Recorded"/>	5	Record Message	8	OK
3	Spell Text to Display	6	Listen to Recording	9	Cancel



Store a Recorded Message

Store a Recorded Message

To Record and Store a Message

1. Decide what you want to say (or do) and the icon sequence you want to use.
2. Remove the keyguard and overlay. Open the Toolbox by pressing the TOOLS key on the front of the Pathfinder case.
3. Activate the STORE CORE key.
4. Activate the icon sequence you decided to use. Activate the OK key.
5. Activate the *Type of Message* key until you see *Recorded*.

☞ If you want to include written text that will be displayed while your message is spoken, activate *Spell Text for Message* and type in the text you want to have displayed. Activate OK when you have finished.

6. Press and hold the *Record Message* option. While you are holding the *Record Message* option, speak your message into the microphone at the top right corner of the Pathfinder case. Release the *Record Message* option as soon as you have finished recording.
7. Press *Listen to Message* to hear your recording.
8. If the message is fine, select OK.
9. If the message is too loud or too soft, press *Record Message* again. Activate the *Microphone Volume* button.
 - If the message was too loud select a lower number.
 - If the message was too soft, select a higher number.
10. Select *Record Message* and *re-record* your message.
11. Activate *Listen to Message* to see if you like the new recording.
12. If the recording is still too loud or too soft, *repeat steps 9-11*.
13. When the recording is the way you want it, activate OK.

☞ Once you have the **Microphone Volume** adjusted to your liking, you should not have to re-adjust it for recording any new messages.

☞ **If you are scanning**, select *Record Message* and hold your switch down until the message is recorded, then release the switch

Erase Recorded Messages

Eventually you will want to erase some recorded messages. The Pathfinder's memory may become full or some messages may just become out of date.

To Erase a Message:

1. Remove the keyguard and overlay and press the TOOLS button to turn the Toolbox On.
2. Activate the ERASE CORE key in the Toolbox.
3. Select the icon sequence you want to erase.
4. Answer YES if you want to erase. Answer NO if you change your mind.

 You cannot edit a recorded message. If you want to change a recorded message, erase it first and then re-record it.

Store a Synthesized Message

To Store a Synthesized Message

1. Remove the keyguard and overlay. Press the TOOLS key to open the Toolbox.
2. Activate the STORE CORE key (M8).
3. Activate the icon sequence you want to use.
4. Activate OK.
- ☞ Read the Instruction Area of the menu for directions.
5. Activate the *Type of Message* key until you see "Synthesized".
6. Activate the *Spell Text for Message* key.
7. Type in the message you want to store.
- ☞ Activate any *Speak* tool to hear the message at any time.
- ☞ Use the *arrow* keys and the *Delete* keys on the Toolbox overlay to edit your message.
8. When you have finished, activate the OK key.
9. If you see some more editing you want to do, activate the *Spell Text for Message* key again and then edit the message. Activate OK.
10. Activate OK again to exit the menu.

☞ **Note:** When storing, SPEAK DISPLAY speaks only the entered text. All other SPEAK tools speak the entered text and the menu prompts.

Edit a Synthesized Message

To Edit a Message

1. Remove the keyguard and overlay. Activate the TOOLS key on the front of the Pathfinder case to open the Toolbox.
2. Activate the STORE CORE key on the Toolbox overlay.
3. Activate the icon sequence of the message you want to edit.
4. Activate *Spell Text for Message*. Edit the message.
- Activate any *Speak* tool to listen to the message.
5. When you have finished, activate OK.
6. Activate OK to exit the menu.

Erase a Synthesized Message

To Erase a Message

1. Remove the keyguard and overlay. Press the TOOLS key on the front of the Pathfinder case to open the Toolbox.
2. Activate the ERASE CORE key on the Toolbox overlay.
3. Activate the icon sequence you want to erase.
4. Activate YES if you want to erase or NO if you change your mind.

Some Helpful Hints

When you store individual words, always store them with a space after the word. Always add one or two spaces after any punctuation.

To create a pause in a message, enclose an underscore (_) in right and left square brackets, i.e., [_] (do not leave spaces) where you want the Pathfinder to pause before continuing. You can specify the length of the pause by adding a number in angle brackets. For example, [_ <1000 >] equals a 1 second pause; [_ <2000 >] equals a 2 second pause, etc.

➡ See **Appendix F** for more DECTalk speech hints.

Storing and Retrieving Messages Using the 64 and 32-Location Keyboards

Storing and retrieving messages with a 64 or 32-location keyboard is exactly the same as with the 128 keyboard. However, you must make sure you have selected to use the 64 or 32 location keyboard in your user area.

If you are using a 64-location keyboard, you can only retrieve messages stored under this keyboard. The same is true for the 32-location keyboard. If you are using a 128-location keyboard you cannot retrieve messages stored under a 64 or 32-location keyboard.

The 64 and 32 location keyboards have blank overlays. This means that you must set up an overlay before you can store messages. See pages 173-182 to learn how to assign keys and tools to an overlay.

1 ICON THEME

 = **Colors**

 = green

 = red

 = yellow

2 ICON THEME

 = **Fruits**

 = kiwi

 = peach

 = pear

1 and 2 Icon Themes

About Themes

Summary:

"Theme" tells the Pathfinder to automatically "hold down" the same first icon or first two icons for any sequences you enter into the theme.

With Minspeak you can organize vocabulary by category, topic or situation. This helps you remember where messages are stored. You are also building strong language skills because you are learning to categorize, store and recall vocabulary.

For example, the icon  could represent the category, "colors".  and  might represent the color "green." The graphic on the facing page shows you a "Colors" category. To add new words to this category, you would store them using the  icon for your *first* icon. This helps you remember that each message is an interjection word. The second icon you use can help you remember the specific word or phrase.

The Pathfinder gives you Theme tools to choose from: PICK 1 (or PICK 2) ICON THEME or SET 1 (or SET 2) ICON THEME and EXIT THEME. If you store PICK 1 ICON THEME under a sequence, all you have to do is activate that sequence, then activate the theme icon you want the Pathfinder to hold down. Once you have selected the theme icon, you do not have to select it again as long as you remain in that theme.

If you store SET 1 ICON THEME under a sequence, you select the theme icon you want to use during the storing process. Every time you activate the sequence you stored SET 1 ICON THEME under, the Pathfinder will automatically hold down the icon you selected for your theme. You do not have to select it.

EXIT THEME allows you to exit a theme for a single out-of-theme message or to exit your theme altogether. Activate EXIT THEME **once** for a single out-of-theme message. Activate EXIT THEME **twice** to completely exit your theme.

When you Store Tools on your Overlay, use INSERT TOOL

In order to use Themes you must assign the different "Theme" tools to keys on your overlay or store them under icon sequences on your overlay. The THEME tools, SET or PICK 1 ICON THEME, SET or PICK 2 ICON THEME and EXIT THEME can be found in the All Tools List. When you 'move' a tool from the Toolbox or the All Tools List to your overlay you are really *copying* the tool and storing it on your overlay.

There are two places for tools in the Pathfinder: the Toolbox contains the most frequently used tools; the All Tools List contains the name of every tool available for use in the Pathfinder.

Tools actually tell the Pathfinder to *do* something. If you go to the Toolbox and activate SPEAK DISPLAY, the Pathfinder will speak whatever is on the display. Before you move a tool from the Toolbox or the All Tools List to your overlay, you must tell the Pathfinder that you don't want it to **do** what the tool says right now, you just want to move the tool.

To tell the Pathfinder that you want to move a tool to your overlay, **you must activate INSERT TOOL and then activate the tool you want to move.** You do this during the storing process. You can only use INSERT TOOL when you are assigning a tool to a key or storing it in an icon sequence. Once a tool is stored on your overlay it will perform its assignment just like it does in the Toolbox.

Assign PICK 1 ICON THEME to a Key or Store it under an Icon Sequence

When you assign PICK 1 ICON THEME to a key or store it under an icon sequence, you are telling the Pathfinder to put you in a 1 Icon Theme and that you will pick the theme icon.

Assign PICK 1 ICON THEME to a Key on your Overlay:

1. Remove keyguard and overlay. Turn Toolbox On.
2. Activate ASSIGN KEY MENU in the Toolbox.
3. Activate the key where you want to make the assignment.
4. Activate the *Minspeak Assignment* option.
5. Select *TOOL* in the *Assign Key As* option.
6. Activate *Enter Assignment*.
7. Activate the *Insert Tool* key on your storing screen.
8. Activate the first letter of the name of the tool you want, or activate *More Items* until you see PICK 1 ICON THEME, then activate it.
9. Activate OK until you exit the menu.

Now when you activate this key in Minspeak Mode, the Pathfinder will ask you to pick an icon for your theme. You will then be in your 1 Icon theme.

Store PICK 1 ICON THEME under an Icon Sequence

1. Remove the keyguard and overlay. Turn the Toolbox On.
2. Activate STORE in the Toolbox.
3. Select the icon sequence you want to store PICK 1 ICON THEME under.
4. Select *Synthesized* in the *Type of Message* option.
5. Activate *Spell Text for Message*.
6. Activate INSERT TOOL.
7. Activate the first letter of the name of the tool you want, or activate *More Items* until you see PICK 1 ICON THEME, then activate it.
8. Activate OK until you exit the menu.

Now when you activate this sequence, the Pathfinder will ask you to pick an icon for your theme. You will then be in a 1 Icon Theme.

A 1 Icon Theme Example:

Using the "Interjections" example from the previous page, you activate the key or icon sequence where you stored PICK 1 ICON THEME. Next you activate the "INTERJ" key. The Pathfinder will hold down this key. Now all you have to do is activate the individual interjection icons (PHONE, HAMMER, etc.).

About PICK 2 ICON THEME

Some vocabulary categories are represented by the same first *two* icons. For example, SENTENCE + EYE might be the icons for the vocabulary category, "Sentences About Me". When you have PICK 2 ICON THEME on your overlay you can activate it, then active SENTENCE and EYE for the Theme icons. The Pathfinder will hold these icons down while you select only the last icon in the sequence to produce a message.

Use the same procedure as with PICK 1 ICON THEME for storing or assigning PICK 2 ICON THEME.

Pathfinder Feedback when you are using a Theme:

When you are in a Theme, the status bar on the bottom of the display will show the Theme icon name followed by a colon and a right angle bracket, e.g., INTERJ:>

Whenever you see an icon name followed by a colon and a right angle bracket in the status bar you will know you are in a theme. When you select the next icon, its name will appear on the status bar in place of the right angle bracket. The colon will remain after the name of the Theme icon. For example, INTERJ: HAMMER.

SET 1 and SET 2 ICON THEME

The SET 1 or SET 2 ICON THEME tools allow you to store the Theme icon(s) at the same time you store the Theme tool. This means that when you activate the key or sequence you stored the Theme tool under, the Pathfinder will automatically hold down the icon (s) you stored with the tool. You don't have to select icons from your overlay as you do with PICK 1 or 2 ICON THEME.

For example, let's say you have a music synthesizer stored in your Pathfinder. It would be ideal if, every time you wanted to play with the synthesizer, you could just activate an icon sequence that would get you into your "Music synthesizer" theme. This sequence would tell the Pathfinder to hold the Theme icons down. Then you could activate single icons to compose a song.

To be able to do this, you store or assign SET 1 (or 2) ICON THEME along with the Theme icons. You might choose the MUSIC and COMPUTER icons for your theme icons.

To Assign SET 1 or 2 ICON THEME to a Key (2 Icon Theme Example)

1. Remove the keyguard and overlay. Turn the toolbox On.
2. Open the ASSIGN KEY menu.
3. Activate the key where you want to make the assignment.
4. Activate the *Minspeak Assignment* option.
5. Select **TOOL** in the *Assign Key As* option.
6. Activate *Enter Assignment*.
7. Activate *Insert Tool*.
8. Activate the first letter of the name of the tool you want, or activate *More Items* until you see SET 1(OR 2) ICON THEME, then activate it.

9. Activate the icons you want to for your Theme, e.g., MUSIC and COMPUTER
10. Activate OK until you exit the menu.

Now when you activate this key in Minspeak mode, the Pathfinder will automatically hold down MUSIC COMPUTER for you.

To Store SET 1 or 2 ICON THEME under an Icon Sequence (2 Icon Theme Example):

1. Remove the keyguard and overlay. Turn the Toolbox On.
2. Activate STORE.
3. Select the sequence you want to store the tool under (SET 1 (or SET 2) ICON THEME). In this example you could use. TOOLS+MUSICAL NOTE.
4. Select *Synthesized* in the *Type of Message* option.
5. Activate *Spell Text for Message*.
6. Activate *Insert Tool*.
7. Activate the first letter of the name of the tool you want, or activate *More Items* until you see PICK 1 ICON THEME, then activate it.
8. Activate the icon(s) you want to use for the theme (e.g., MUSIC COMPUTER).
9. Activate OK until you exit the menu.
10. Now when you activate TOOLS+MUSIC, the Pathfinder will automatically hold down the MUSIC+COMPUTER icons.

About EXIT THEME

Summary:

Activate EXIT THEME **once** to exit a theme temporarily. Activate it **twice** to permanently exit a theme.

➡ If you plan to use themes, you will have to assign the EXIT THEME tool (from the All Tools List) to a key or store it under an icon sequence on your overlay. See **Notes:** below.

Activate EXIT THEME **once** to temporarily exit your current theme for a one-sequence activation. When the message is completed, the Pathfinder automatically returns you to your theme.

If you want to leave a theme completely, activate EXIT THEME **twice**.

If you assign EXIT THEME to a key:

You can use it no matter what theme you are in.

If you store EXIT THEME under an icon sequence:

Store it using the same icons as your theme icons, that is, store it *within* your theme.

For Temporary Exits from the Theme:

If you know that, in general, you are going to leave your theme only to make a single sequence activation now and then, store EXIT THEME *once* under your theme.

For example, Let's say your theme is MUSIC COMPUTER and you enter it by activating TOOLS+MUSICAL NOTE.

1. Activate TOOLS MUSICAL NOTE.
2. Check the status bar at the bottom of the screen. If it shows MUSIC COMPUTER you are in your theme.
3. Go to the Toolbox and activate STORE.

4. Activate the icon you want to use to represent EXIT THEME (perhaps the LEFT ARROW icon).
5. Activate OK.
6. Activate *Spell Text for Message*, activate *Insert Tool*. Activate *More Items* until you see EXIT THEME, or type the letter "e" and then activate *More Items* until you see EXIT THEME. Select it.
7. Activate OK until you exit the menu.

Now when you are in your theme and you want to leave it to speak one outside message, just activate the Left Arrow icon on your overlay.

To Permanently Exit a Theme:

When you want to permanently exit your theme, you must have EXIT THEME stored *twice* in a single sequence.

For example:

1. Activate TOOLS MUSICAL NOTE.
2. Check the status bar at the bottom of the screen. If it shows MUSIC COMPUTER you are in your theme.
3. Go to the Toolbox and activate STORE.
4. Activate the icon you want to use to represent *permanently exiting the theme* (perhaps the NO icon).
5. Activate OK.
6. Activate *Spell Text for Message*, activate *Insert Tool*.
7. Activate *More Items* until you see EXIT THEME, or type the letter "e" and then activate *More Items* until you see EXIT THEME. Select it.
8. Select *Insert Tool* and EXIT THEME again.
9. Activate OK until you exit the menu.

About Single Meaning Icons

You may want words or phrases like, "Yes", "No", or "I need help!" stored under single icons. With these words or phrases, you can respond to a question while speaking on the phone, answer questions if you are unable to indicate the answer by gesturing or get help quickly in any situation. However, if messages like these are stored under single icons, and you are in a 1 or 2 ICON THEME, you can only access the non-theme messages by activating EXIT THEME and then activating the icon the message is stored under.

If you make "Yes", "No", or "I need help!" a Single Meaning icon, you can access the message when you are in a theme by simply activating the icon. You don't have to exit your theme. For example, if you are playing your music synthesizer (under the 2 ICON THEME: MUSIC COMPUTER), you can answer a question just by selecting the YES or NO key on the overlay. You don't need to exit the theme. This can be very valuable in an emergency when a quick response is needed.

A Single Meaning icon will **not** work if the icon you have chosen to be the Single Meaning icon is already stored under a sequence which could be used in a 1 or 2 icon theme.

For example, let's say you stored a musical note for your synthesizer under MUSIC COMPUTER MEDICAL. Let's also say that under the single icon, MEDICAL, you stored, "I need help!" You make, "I need help!" a Single Meaning icon (in the ASSIGN KEY menu) because it's a good phrase to have available at any time.

Now, imagine you are playing your synthesizer under the 2-icon theme: MUSIC COMPUTER. If you activate the MEDICAL icon, you will get the musical note instead of "I need help!" even though you made MEDICAL a Single Meaning icon. When you are in a theme and you activate a single icon, the Pathfinder will first look to see if that icon is used in your theme.

Activate ASSIGN KEY menu in the Toolbox.

ASSIGN KEY MENU - Choose assignments for key [E5]			1
1 Minspeak Assignment EAR	4 Swap Key	7 More Items	
2 Spell Assignment e	5 Clear Key Contents	8 OK	
3 Computer Assignment «SERIAL-MAR...	6 Select next key to define	9 Cancel	

Select *Minspeak Assignment*

Enter Minspeak Assignment for key [E5], then select OK		
Icon Name: EAR		
2 Assign Key As: ICON	5 Icon Type Single Meaning	8 OK
3 Rename Icon		

Select *ICON* and *Single Meaning*

Activate OK when finished.

Store a message under the icon.

Single Meaning Icons

Since MUSIC COMPUTER MEDICAL is a legitimate theme sequence, Pathfinder will select that sequence rather than "I need help!" A Single Meaning icon will work only if it is not connected in any other sequences to the theme you are using.

If you want an icon to be a Single Meaning icon no matter what theme you are in, you should consider selecting an isolated key on your overlay, such as we have done with the YES and NO keys.


To Create a Single Meaning Icon:

1. Remove the keyguard and overlay. Turn the Toolbox On.
2. Activate the ASSIGN KEY menu.
3. Select the Minspeak Assignment option.
4. Select ICON in the Assign Key As option.
5. Select Single Meaning in the Icon Type option.
6. Activate OK until you exit the menu.

You can now store the message you want under that icon. Use the regular storing process.

 Note: Text stored as "Single Meaning" in Minspeak Mode is also available in THEMES. For example, YES and NO in Unity are available in a Theme.

ACTIVITY MENU			1		
1	Activities <input type="checkbox"/> OFF	4	Move Activity	7	More Items
2	Create Activity	5	Modify Activity	8	OK
3	Dynamic Activities	6	Modify Item Within Activity	9	Cancel

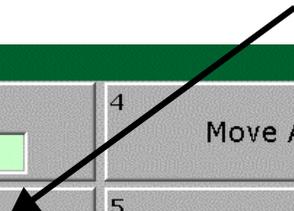

Create an Activity in the ACTIVITY Menu

About Creating an Activity

Regardless of how an activity is going to be accessed (through the Activity row list or from an icon), the first step is to create the activity.

To create an activity or to modify an already existing activity, open the ACTIVITY menu in the Toolbox. Use the *Create New Activity* option or the *Modify Item within Activity* option.

Perhaps you would like to create a "Baking Cookies" activity. The vocabulary for this activity can consist of words that are not part of your core vocabulary. You might need words like, "bowl", "flour", "teaspoon", "stir", "bake", etc. Each key that you define for this activity can contain one of these words. When you open this activity, you can use your core vocabulary to say, "I want" and then you can activate the "flour" key in the "Baking Cookies" activity. Or you might use your core vocabulary to say, "I need" and then activate "teaspoon" from the "Baking Cookies" activity.

By using these guidelines for creating activities, you can accomplish a number of goals. You will continue to use and learn your core vocabulary. You will not waste communication time by duplicating your core vocabulary in an activity. You will be able to easily access words that you don't use every day but which are necessary for specific activities. You will be able to participate in any communication that takes place during the activity.

Pathfinder icons are listed in categories and are already included in the device. You can select icons for activities from these categories. The Pathfinder takes you to the icon category page when you are creating or modifying activities. You can also import other icons, photos, etc. into the device. See the *MTI Manual* for more information about importing graphics.


Activate *Create Activity* from Screen 1




Spell name for Activity



Choose an icon: spell the first few letters of an icon name or use arrows to move through categories

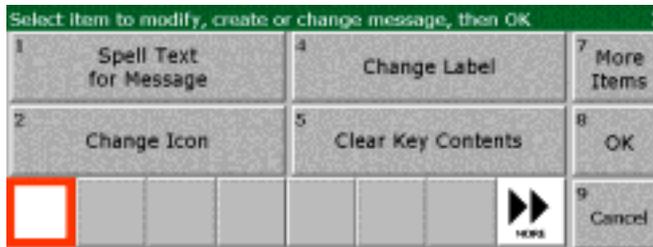

Add items to the activity beginning with the blank square

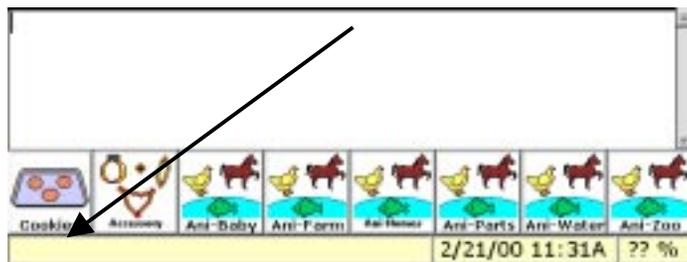

When finished, you will see this in your Activity Row

Create a New Activity: An Example

Baking Cookies

1. Go to the Toolbox and open the ACTIVITY menu.
2. Spell the name of your new activity. For example you might spell, "Cookies". Activate OK when you have finished spelling.
3. Choose an icon to represent the new activity from one of the Icon Categories. ➤ You can use the *Spell Icon To Find* key to spell the first few letters of an icon name. Activate OK. The Pathfinder will bring up all the icons that might have the name you are looking for. Use the arrow keys on your overlay to navigate the icons or categories. When the icon you want is highlighted, activate OK. You might want to open the AROUND THE HOME category and then open the KITCHEN sub-category. You could then select the *Cookies* icon, or you could just spell "coo" and select the COOKIE icon.
➤ This icon **will not** appear on the "Modify Item Within Activity" screen as you continue to create the items within the activity. It **will** appear in the Activity Row on your display when you have finished creating the activity and exited the menu.
4. Now you can create the items in your "Baking Cookies" activity. The first (blank) key in the activity is highlighted with a red box around it. This will be your first activity item.
5. Select the *Choose Icon* option and select an icon for the key. ➤ If the name of the icon you select represents what you want to say, you will not have to type in a message. The Pathfinder will speak the label on the icon.
6. **Options:** If you do not like the label for an icon, you can select the *Change Label* option. For example, you might want to use the "oven" icon, but you would like the Pathfinder to speak, "bake". You can change the label to "bake". If you want the Pathfinder to speak, "Bake at 350 degrees", you will also want to select the *Spell Text for Message* option, and type in the message.

ACTIVITY MENU			1
1	Activities <input type="checkbox"/> OFF	4	Move Activity
2	Create Activity	5	Modify Activity
3	Dynamic Activities	6	Modify Item Within Activity
		7	More Items
		8	OK
		9	Cancel

MODIFY ACTIVITY MENU			
1	Change Icon	4	Hide/Show Activity <input type="checkbox"/> Show
2	Change Label	5	Modify Item Within Activity
 Cookies		8	OK
		9	Cancel

More ACTIVITY Menu Options

If you want a blank label, select the **Change Label** option and activate a space instead of typing in a word for the label.

If you want to record a message to go with the icon, activate **More Items** twice and select **Recorded** in the **Type of Message** box. Hold down **Record Message** and speak your message. If you also want to include text, activate **More Items** once and then activate the #1 option: **Spell Text for Message**.

If you need an additional row in your activity, select **More Items** and then select the #1 option: **Add New Row**. When you add a row, the Pathfinder automatically creates a **More Items** key for that row. This allows you to move between the rows when communicating. While creating an activity with more than one row, use the **Previous Row** and **Next Row** options to move between rows. If you want to delete a row in your activity, make sure the row you want to delete is showing on the display. Activate the **Delete Current Row** option.

If you want to rearrange the order of the keys in your activity, make sure the key you want to move is highlighted. Select **More Items**, then select the **Change Position Move Left** or **Change Position Move Right** option.

➡ You can use the arrow keys on your overlay to step through the activities, one-by-one. Activate the OK key on your overlay to select a highlighted activity.

7. When you have finished with all the options that you want to use for this key, activate OK.
8. Activate the next key in the activity that you want to define and repeat these steps.

Delete an activity by selecting the activity and then activating the **Delete Activity** option.

➡ You can erase the contents of any key in the activity by selecting the **Clear Key Contents** option. Then you can just start over with the key.

ACTIVITY MENU			1
1	Activities <input type="checkbox"/> OFF	4	Move Activity
2	Create Activity	5	Modify Activity
3	Dynamic Activities	6	Modify Item Within Activity
		7	More Items
		8	OK
		9	Cancel

DYNAMIC ACTIVITY MENU			
1	Add Dynamic Activity	4	Dynamic Activities <input type="checkbox"/> ON
2	Delete Dynamic Activity	5	Activity Persistence <input type="checkbox"/> Temporary
3	Review Dynamic Activities	8	OK
		9	Cancel

The *Dynamic Activities* Sub-Menu

Dynamic Activities

A Dynamic Activity is an already created activity that you choose an icon or an icon sequence to represent. Once you make an activity dynamic, you can be in any other activity, activate the icon you selected for your dynamic activity, speak a message from it and return automatically to the original activity.

Let's make the "Cookies" Activity dynamic. You make an activity dynamic by selecting the *Dynamic Activities* option in the ACTIVITY menu and then by selecting an icon to represent the activity.

To Make an Activity Dynamic:

1. Go to the Toolbox and open the ACTIVITY menu.
 2. Activate the Dynamic Activities option.
 3. Make sure that the *Dynamic Activities* option is set to *On*.
 4. Make sure the *Activity Persistence* option is set to *Temporary*. This means that the activity will stay dynamic for one activation only, then take you back to where you were. If you want to remain in the new activity and not return to your original one, set this option to *Permanent*.
- ☞ Regular activities cannot be set to *Permanent* or *Temporary* unless they are made *Dynamic* first. Only Dynamic Activities can be *Permanent* or *Temporary*.
5. Activate the *Add Dynamic Activity* option.
 6. Select a sequence (a one-icon sequence is quickest) that will remind you of the activity. For example, for the "Cookies" activity you might select the KITCHEN icon which is the pan with a thermometer in it. Activate OK when you are finished. ☞ You will not mask any vocabulary when you select a sequence for an activity.

7. Now activate the activity you want to make dynamic (e.g., "Cookies"). You may have to activate *More Items* before you see the activity you want.
8. Activate OK. Activate OK to exit the menu.
9. Open any activity. Then activate the KITCHEN icon in your core vocabulary area.
10. The "Cookies" activity opens on your display. Activate an icon in the activity.
11. The word is spoken and you return to your original activity.

➡ Since a Dynamic Activity is intended to be selected from an icon or an icon sequence, you may not want to see a Dynamic Activity listed in the Activity list. Use the *Hide Activity* option in the ACTIVITY menu to keep a Dynamic Activity from being visible in the Activity list.

Erase a Dynamic Activity; Review All Dynamic Activities

To Erase a Dynamic Activity:

1. Go to the Toolbox and open the ACTIVITY menu. Select the *Dynamic Activities* option.
2. Activate the *Delete Dynamic Activity* option.
3. Activate the icon sequence for the Dynamic Activity you want to delete.
4. Select YES if you are ready to delete the activity.
5. Activate OK to exit the menu. Activate OK again to completely exit.

To Review Your Dynamic Activities

1. Go to the toolbox and open the ACTIVITY menu. Select the Dynamic Activities menu.
2. Activate the *Review Dynamic Activity* option.

Select *Step or Auto* as the method you want to use to review the

3. Activities.  If you want to *Review to Output*, see *Appendix C for information about connecting to a printer or computer*.
4. The Pathfinder shows you the sequences and the names of all the Dynamic Activities. If you are using *Step*, activate *More Items* if you have more than one page of Dynamic Activities.
5. Activate OK when you are finished and exit the ACTIVITY menu.

Modify an Activity

Once you have created an activity, you might like to change it in some way later. Use the **Modify Item Within Activity** option in the ACTIVITY menu to change an existing activity.

To Modify an Activity:

1. Go to the Toolbox and open the ACTIVITY menu.
2. Select the activity you want to modify. If necessary, activate **More Items** in the ACTIVITY menu until you see **Next Activities** or **Previous Activities**. Select the option you want.
3. You can change the icon or label of the main activity key. You can change its position in the Activity Row. You can hide it (or show it if you have already hidden it). You can also delete the entire activity.
4. If you want to change one of the items *within* the activity, select the **Modify Item within Activity** option. You can spell a new message for the key; change its icon and/or its label; clear the key's contents and start over; add a new row of keys to the activity or delete a whole row; change the position of the key in the row. You can use **Previous Row** and **Next Row** to navigate through more than one row in the activity.

Linking Activities/Creating Sub-Categories

Let's say you have created three separate activities: FRUIT names, MEAT names and DESSERT names. When you activate your FRUIT activity you see the names of fruits. You might like to link these similar activities as sub-categories under a dynamic activity called FOODS. You can choose the APPLE icon to represent the dynamic FOODS activity. This way, you can activate APPLE on your overlay and see the MEAT, FRUIT and DESSERT activities in your activity row.

To Link Activities and Create Sub-Categories

1. Create each activity that you will eventually want to link into sub-categories (e.g., MEATS, FRUITS, DESSERTS). Follow the steps for **Creating an Activity** on page 151-153.
2. Now, create a new activity and call it FOOD.
3. For the first location in the FOOD activity (this will be highlighted on your display), activate *Spell Text for Message*.
4. Activate INSERT TOOL.
5. Activate the letter, "S". Activate *More Items* until you see SET ACTIVITY. Activate SET ACTIVITY.
6. Choose the activity you want to set, for example, MEATS.
7. Activate OK
8. Use *Change Icon* and add an icon to this sub-category key.
9. Use *Change Label* and name the subcategory key (MEATS)
10. Highlight the next location in the activity row and repeat steps x-x for each of the sub-categories you want to add to your FOODS activity.
11. When you have, finished activate OK
12. Now, activate *More Items* until you see *Hide/Show Activity*. Set this to *Hide*.
13. Activate *More Items* until you see *Dynamic Activities*.
14. Activate *Dynamic Activities*.
15. Activate *Add Dynamic Activity*.
16. Select the icon on your overlay that you want to represent your FOOD activity, for example, select APPLE.

17. Activate OK.
18. On your display, select the FOOD activity from list of activities.
19. Make sure *Dynamic Activities* is set to *On*.
20. Activate OK twice.

Close the Toolbox.

Go to Minspeak mode.

Select the APPLE icon on your overlay.

You should see the MEATS, DESSERTS and FRUITS activities in your Activity Row.

Copy Activities from One User Area to Another

Eventually you may create rows of different activities in a User Area. If you set up a new User Area, you might like to have all the activities in your other User Area available in the new one. Instead of creating the activities all over again for the new area, the Pathfinder allows you to copy activities from one area to another.

➡ **A Note about Duplicate Activities:** If you are copying activities to an area that already has some activities in it, the Pathfinder will check for duplicate activities before copying. This means the Pathfinder will check to see if there are duplicate activity **names**. If you have an activity in one area that has the same name as an activity that you are copying, the Pathfinder will think they are the same activity whether they are or not. If you choose **Copy All Activities**, the Pathfinder will replace any duplicates that it finds with the activities you are copying. If you choose **Copy Selected Activities**, you can tell the Pathfinder to replace a duplicated activity or to keep the old activity and add the new one with same name.

To Copy Activities:

1. In the Toolbox, open the USER AREA menu.
2. Activate **More Items**.
3. If you want to copy only certain activities from one area to another, choose **Copy Selected Activities**. If you want to copy all the activities in one area to another area, choose **Copy All Activities**.
4. Select the area you want to copy the activities **to**.
5. ➡ To **Copy all Activities** in one area to another area:
Activate YES to copy or NO if you change your mind.. All the activities will be copied to the new area you selected. There will be no duplicates.

➡ **To Copy Selected Activities:**

If you are touching the display screen to select items:

Activate an activity. It turns gray to signify that it has been selected.

If you change your mind, activate the same activity again. It is no longer selected.

If you are using the arrow keys on the keyboard:

Activate an arrow key. The activity in the upper left corner is highlighted in a red box. To select this activity, activate the OK key on your **overlay**. The activity turns gray to signify that it has been selected. If you change your mind, activate the OK key again to de-select the activity. Activate the arrow key you want to highlight another activity. Activate OK on the overlay.

Select all the activities you want to copy. If you have more than one row of activities, the activities you have selected to copy in one row will remain selected while you move on to another row.

Activate the OK key on the display screen when you are ready for the Pathfinder to copy the activities you have selected.

If there are duplicates, the Pathfinder will ask you if you want to copy the activity. Activate YES or NO. If you activate YES, there will be two activities with that name in your new area.

6. Switch to the new User Area to check that the activities transferred.

Review Vocabulary

Summary:

You can review vocabulary stored under a specific icon, or you can see all the messages in the Pathfinder.

You can review the vocabulary in your Pathfinder. You can find out if there are messages you are no longer using and you can erase them. You can refresh your memory about what messages are already stored. If you want, you can print out the whole Pathfinder vocabulary.

To Review the Pathfinder Vocabulary

Remove the keyguard and overlay. Press the TOOLS button on the front of the Pathfinder case to open the Toolbox.

Activate the REVIEW VOCABULARY menu.

All Sequences or By First Icon

You can review all the sequences in the Pathfinder, or you can select a beginning icon and review all the sequences that begin with that icon.

Step or Auto

Review the vocabulary by activating a switch or key every time you want to move to a new sequence (step), or let the Pathfinder do it for you (auto).

Review to Display

The vocabulary is shown on the Pathfinder's display screen.

➡ You can edit a sequence or delete a sequence by selecting on the display screen the sequence that you want, or by selecting the corresponding number on your keyboard. Then select to edit the message or to delete it.

VOCABULARY REVIEW MENU		
1	Review All Sequences	4 Review to Display
2	Data Logging	5 Review to Display & Speaker 8 OK
3	Review Method STEP	6 Review to Output 9 Cancel

REVIEW VOCABULARY Menu

Review to Display and Speaker

You can see the vocabulary on the display and hear the Pathfinder speak it at the same time.

Review to Output (*See Appendix C for specific information*)

You can send the vocabulary to a printer or computer. You must be connected with the correct cable and the options in the OUTPUT menu must be set before the information will appear on a printer or computer.

Data Logging On/Off (*see Appendix B for specific information*)

When this option is On, any text (whether from Minspeak or from another program) that the Pathfinder user brings up on the screen is logged by date, time and a copy of the actual text. The data in this log can then be analyzed by clinicians. This is called Language Activity Monitoring (LAM).

To analyze language data, the Pathfinder must be connected to another computer with a serial cable. To upload the information, the computer must be running a terminal program. (If you do not have a terminal program such as "HyperTerminal" for Windows or "Comet" for Macintosh, you can find these programs on the Internet. All terminal programs work differently. After downloading, use your help screens to get started.) You must also set the appropriate settings in the Pathfinder's OUTPUT menu. When you are connected and ready, activate *Send Data Log to Computer*. (See below).

Send Data Log to Computer: sends any data that has accumulated in the log to a computer for editing and analysis.

Clear Data Log: tells the Pathfinder to delete all information in the data log.

Personalize Your Overlay

Menus to Help You Personalize Your Overlay

**ASSIGN
KEY
MENU**

**AUDITORY
PROMPT
MENU**

**FEED-
BACK
MENU**

**MAIN-
TENANCE
MENU**

**Personalize your Overlay using Options
from these Menus**

About Personalizing Your Overlay

In order to communicate rapidly and effectively, you will want to personalize, or customize, the overlay that you use the most. Whether you are using one of the Unity overlays or creating your own overlay from a blank one, you will want the tools, icons, spelling characters and computer commands that are most helpful to you on it.

There are many things you can do to personalize your overlay. If there are icons on your overlay that have no meaning for you, you can change them or rename them to something familiar. You can delete tools or spelling characters that you never use and add tools and spelling characters from the Toolbox that are more useful for you. You can swap key positions or disable keys so that they won't work. You can also assign keys, tools, spelling characters and instructions for their use to specific keys or under icon sequences.


To personalize your overlay you can use INSERT TOOL and the All Tools List. You will also use the ASSIGN KEY menu to perform many of the options mentioned above.

➡ If you want to personalize an area that is already set up, simply select it from the USER AREA menu.

➡ If you want to set up a totally new area, select an unused area from the USER AREA menu. (If you do not have an unused area, you may have to erase an already used area.) Select one of the pre-defined overlays or select a blank overlay from the *Overlay Menu* option. Finally, **Name the Area**. (See pages 47-55 for more complete USER AREA menu information.)



Activate ASSIGN KEY menu
Then activate a key to view or modify

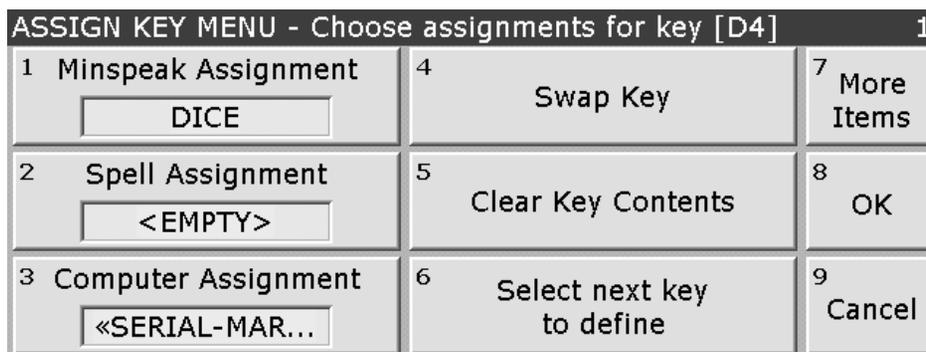


Once you have selected a key to view or modify,
you see this menu

The ASSIGN KEY Menu

You will use the ASSIGN KEY menu when you want to assign tools or icons to a key in Minspeak or Spell, when you want to assign spelling characters to a key in Minspeak or Spell or when you want to assign a computer command to a key. This menu also lets you swap keys, copy keys and hide keys.

Most keys have a Minspeak assignment and a Spell assignment. The Minspeak assignment is usually an icon and the spell assignment is usually an alphabet letter or a number or punctuation character.


Some keys also have a Computer assignment.

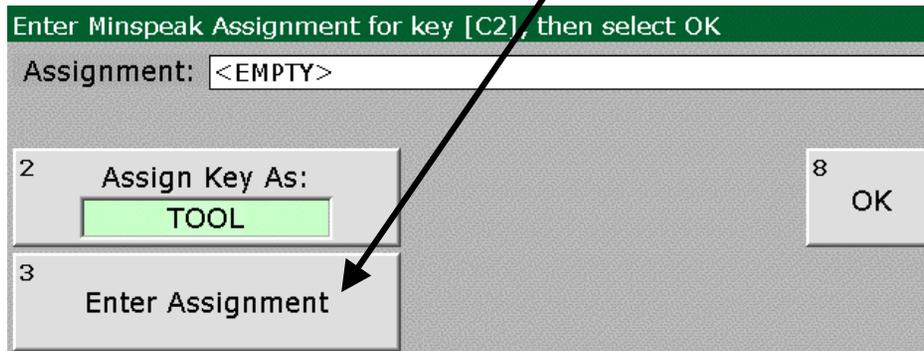

Activate *Enter Assignment*

Activate *Insert Tool*

Spell the first letter of the tool you want. Pathfinder will take you to the page of tools that begins with the letter you typed; or you can just activate *More Items* until you see the tool you want.

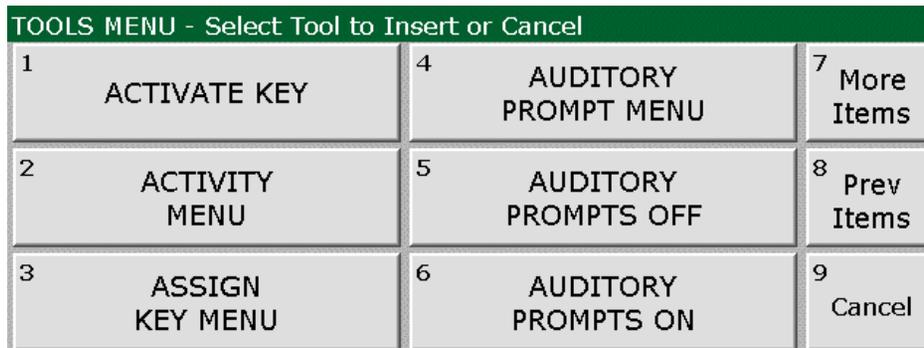

When you find it, activate the tool you want to insert.

About Copying Tools from the Toolbox or All Tools List to Your Overlay

If you use a certain tool frequently, you may want to store it on your communication overlay. When you store a tool under an icon sequence or assign it to a key on your overlay, it is always available to you. You don't have to open the Toolbox to get to it.

You can store tools under icon sequences, you can assign tool to keys on your overlay, and you can even combine tools and instructions together and store them in sequences or assign them to keys.

There are many tools available in the Pathfinder; so many, in fact, that they will not all fit in the Toolbox. The All Tools List is a list of all the tools that are available for use. You must activate INSERT TOOL before you activate the tool you want to assign or store.

The INSERT TOOL key allows you to copy tools from the All Tools List and assign them to keys or store them in sequences. INSERT TOOL only works when you are in storing mode or when you are in the ASSIGN KEY menu. When you activate INSERT TOOL, you are telling the Pathfinder that the next key you activate is going to be a tool that you want the Pathfinder to copy for you. You *must* activate INSERT TOOL before each tool that you want to copy.

Assign a Tool to a Key in Minspeak Mode

To assign tools to keys you must use the ASSIGN KEY menu and INSERT TOOL.

To Assign a Tool to a Key in Minspeak Mode:

1. Remove the keyguard and overlay. Press the TOOLS key on the front of the Pathfinder case to open the Toolbox.
2. Activate the ASSIGN KEY menu.
3. Select the key where you want to assign the tool.
4. You see the current key assignments, as well as some other options.
5. Activate the *Minspeak Assignment* option.  If you want the tool to be available in all modes activate *More Items* and turn option #2 *On*. Then return to the previous screen and select option #1,2 or 3.)
6. Activate the *Assign key as* option to read, *Tool*.
7. Activate *Enter Assignment*.
8. Activate the INSERT TOOL key on the Pathfinder.
9. From the list, select the tool you want to add to the key. Activate *More Items* to see more tools or begin to spell the name of the tool.
10. If you want to define more keys, select *Choose Next Key to Define*.
11. When finished, activate OK until you have exited the menu.

Note:

When You Assign a Tool to a Key in Minspeak:

The tool will always be available for use if you leave the *TOOL* option visible in the *Assign Key As* box when you exit the menu. If you leave the *ICON* option visible, the icon is available, not the tool. If you do not want the icon to be available at all, select *ICON* in the *Assign Key As* option. Select *Rename Icon*. Activate DELETE WORD. Activate OK. Exit the menu.

Assign a Minspeak Icon to a Key

You can also assign an icon to a key. If you are using a pre-defined overlay, all the icons on it are named (e.g., SUN, VOLT, DICE). If you are using a blank overlay all keys are named by their location code (e.g., E-5, A1, N3).

To Assign a Minspeak Icon to a Key in Minspeak Mode:

1. Remove the keyguard and overlay. Press the TOOLS key on the front of the Pathfinder case to open the Toolbox.
2. Activate the ASSIGN KEY menu.
3. Select the key where you want to assign the icon.
4. You see the current key assignments, as well as some other options.
5. Activate the *Minspeak Assignment* option. ➡ If you want the tool to be available in all modes activate *More Items* and turn option #1 *On*. Then return to the previous screen and select option #1,2 or 3.)
6. Activate the *Assign key as* option to read, *Icon*.
7. You see the icon name or key location on the screen. Select *Normal* as the *Icon Type*. (See p.145 for information about Single Meaning icons and *Appendix G* for information about Dependent icons.)
8. If you do not like the icon name, select *Rename Icon* and give the icon a new name. ➡ If you are setting up an overlay for the first time, you will probably want to name all the icons instead of using the pre-assigned location names
9. If you want to define more keys, select *Choose Next Key to Define*.
10. When you are finished, activate OK until you have exited the menu.

Using *Rename Icon*

When you activate this option you are asked to spell a new name for the icon. You can spell in the name of the original icon or you can spell in the name of a different icon. For example, if the original icon name is SUN, you might want to rename it to HAPPY. ➡ Your *Unity* manual tells you the original names of all the icons. **There is an 8 character limit for icon names.**

Assign a Tool to Spell Mode

To Assign a Tool to a Key in Spell Mode:

1. Remove the keyguard and overlay. Press the TOOLS key on the front of the Pathfinder case to open the Toolbox.
2. Activate the ASSIGN KEY menu.
3. Select the key where you want to assign the tool.
4. You see the current key assignments, as well as some other options.
5. Activate the *Spell Assignment* option.  If you want the tool to be available in all modes activate *More Items* and turn option #1 *On*. Then return to the previous screen and select *Spell Assignment*, option #2.
6. Delete the current assignment if there is one.
7. Activate the INSERT TOOL key on the Pathfinder.
8. From the list, select the tool you want to add to the key. Activate *More Items* to see more tools.
9. If you want to define more keys, select *Choose Next Key to Define*.
10. When you have finished assigning keys, activate OK until you have exited the menu.

Notes:


When you are choosing a tool from the all Tools List, you can type the first letter of the tool you want and the Pathfinder will display the first screen that contains tool names that begin with the letter you typed. This makes for quicker navigation.

You can access the ASSIGN KEY menu without removing the keyguard by activating the STATUS key on the left front of the Pathfinder case. Use the display to navigate.


Assign a Tool to Computer Mode

To Assign a Tool to a Key in Computer Mode

1. Remove the keyguard and overlay. Press the TOOLS key on the front of the Pathfinder case to open the Toolbox.
2. Activate the ASSIGN KEY menu.
3. Select the key where you want to assign the tool.
4. You see the current key assignments, as well as some other options.
5. Activate the *Computer Assignment* option.  If you want the tool to be available in all modes activate *More Items* and turn option #2 *On*. Then return to the previous screen and select option #3.)
6. Delete the current assignment if there is one.
7. Activate the INSERT TOOL key on the Pathfinder.
8. From the list, select the tool you want to add to the key. Activate *More Items* to see more tools.
9. If you want to define more keys, select *Choose Next Key to Define*.
10. When you have finished assigning keys, activate OK until you have exited the menu.

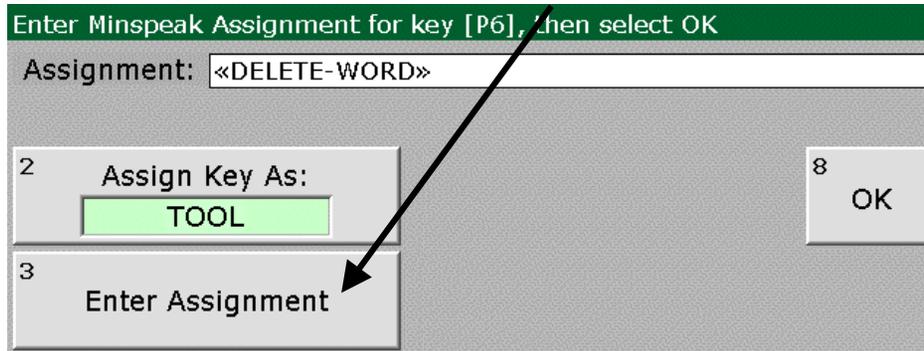

Activate *Enter Assignment*

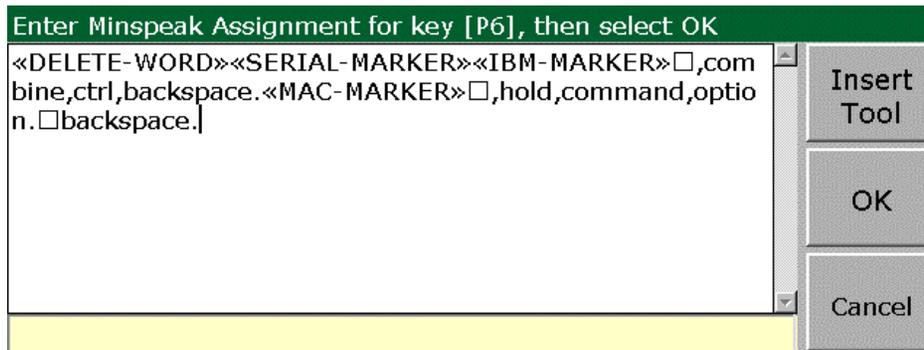

Activate DELETE CHARACTER or CLEAR DISPLAY on the Pathfinder to erase the tool.

Erase a Tool from a Key

To Erase a Tool from a Key:

1. Turn the Toolbox On and activate ASSIGN KEY menu.
2. Activate the key whose tool you want to erase.
3. Activate **TOOL** in the *Assign Key As* option.
4. Activate Enter Assignment.
5. Activate DELETE WORD or CLEAR DISPLAY on the Pathfinder to erase the tool.
6. Activate OK, or activate INSERT TOOL and assign a new tool, then activate OK.
7. Exit the menu.

Using *Clear Key Contents*

You can also erase a key by using the *Clear Key Contents* option. When you activate this option, it erases everything from the key and the key is blank. You can then assign new tools or spelling characters to the key.

Erase a Spelling Assignment from a Key

1. Turn the Toolbox On and activate the ASSIGN KEY menu.
2. Activate the key whose assignments you want to change.
3. Activate *Spell Assignment*.
4. Activate DELETE WORD to erase the assignment.
5. Activate OK, or type in a new assignment if you wish, then activate OK.
6. Exit the menu.

Copying a Tool and Storing it Under an Icon Sequence

Let's say you would like to copy the BEEP ON/OFF tool and store it under an icon sequence on your overlay.

To Copy a Tool and Store it Under an Icon Sequence:

1. Remove the keyguard and overlay. Press the TOOLS button on the front of the Pathfinder case to turn the Toolbox On.
2. Activate STORE CORE in the Toolbox.
3. Select the icon sequence you want to store the tool under and activate OK.
4. Select *Synthesized* in the *Type of Message* option.
5. Select the *Spell Text for Message* option.
6. Activate the INSERT TOOL key in the Toolbox.
7. You see the first screen of the All Tools List. Activate *More Items* until you see BEEP ON/OFF. Select it.
8. Activate OK.
9. Finish storing and exit the menu.
10. Make sure you are in Minspeak.
11. Activate the sequence you stored the tool under.
12. Activate a key. If the beep was On before, it should now be silent. If it was Off before, it should now be On.

To Erase the Sequence:

1. Turn the Toolbox On.
2. Activate ERASE CORE
3. Activate the icon sequence you stored the tool under.
4. Activate OK.
5. Activate YES to erase the sequence or NO if you change your mind.

About Macros

Summary:

Assign instructions that combine tools and text to a key or an icon sequence.

When you store tools and text together, you are making a macro. A macro is just a set of instructions that tell the Pathfinder to do something. Macros are often used as shortcuts for tasks that would normally require a lot of key activations.

When you copy and store a tool on your overlay you are creating a basic macro. You can also create and store macros that tell the Pathfinder to do a lot more than to just activate a tool. You can create macros that change your speaking voice in the middle of a sentence and return it to normal. You can create macros that involve changes in menu settings so that you don't have to go to the Toolbox to change settings.

You can assign macros to keys on your overlay or you can store them under icon sequences.

Assign a Macro to a Key

As an example, let's store the tool, BEEP ON/OFF, and add the words "Beep! Beep!" for the Pathfinder to speak.

To Assign Tools and Text to a Key:

1. Turn the Toolbox On and activate the ASSIGN KEY menu.
2. Select the key on your overlay where you want to store the macro.
3. Select *Minspeak Assignment*.
4. Select TOOL in the *Assign Key As* option.
5. Select *Enter Assignment*.
6. Activate INSERT TOOL.
7. Activate *More Items* until you see BEEP ON/OFF. Activate it.
8. Type Beep! Beep!
9. Activate OK.
10. Activate OK.
11. Go to your Minspeak overlay. Make sure your speech is On. Activate the key where you stored your macro. The Pathfinder should speak, "Beep! Beep!" and turn the activation beep Off (if it was already on. If the beep was Off to begin with, this should turn it On.)

Store a Macro Under an Icon Sequence

We will use the BEEP ON/OFF example from the previous page.

To Store a Macro Under an Icon Sequence:

1. Turn the Toolbox On and activate STORE CORE.
2. Select the icon sequence you want to use to store the macro.
3. Activate OK.
4. Activate *Spell Text for Message*.
5. Activate INSERT TOOL.
6. Activate *More Items* until you see BEEP ON/OFF. Activate it.
7. Type in Beep! Beep!
8. Activate OK.
9. Activate OK.
10. Go to your Minspeak overlay and activate the sequence you just stored the macro under. The Pathfinder should speak, "Beep! Beep!" and turn the feedback beeps Off. Activate the sequence again to turn them back On.

You can erase macros that you assigned to a key or a sequence the same way you erase tools from keys or sequences.

Some Common Macros

To switch between already-set-up user areas:

INSERT TOOL <USER AREA MENU>, select menu option #1, select the area you want to switch to, INSERT TOOL<OK>:

To Store the Notebook Menu so you can select from a list of menus

INSERT TOOL<NOTEBOOK MENU>, select menu option #3
(Select OK to exit the menu. You do not need to add OK to the macro.)

If you create messages in sentences and with punctuation, you can use the following macro to place a period at the end of a sentence, speak the sentence (if *Speech* is turned On) and set SHIFT for the beginning letter of a new sentence.

INSERT TOOL<DELETE CHARACTER> (space) INSERT
TOOL<CTRL>INSERT TOOL<SPEAK SENTENCE>INSERT TOOL<SHIFT>

Swap Key

After using your Unity overlay for awhile, you may find that you would like to swap some keys so that they would be in a better position for you to use. You can do this in the ASSIGN KEY menu.

To Swap Keys on your Overlay:

1. Turn on the Toolbox and activate the ASSIGN KEY menu.
2. Activate one of the keys you want to swap.
3. Activate the **Swap Key** option.
4. Activate the key you want to swap with.
5. Activate OK.
6. The keys are swapped.

➡ If you decide you don't want to swap the keys you have selected, activate **Cancel** and the keys will retain their original assignments.

➡ When you swap keys everything assigned to that key, including messages, is moved with the key.

Hide/Show Key

If you are just learning about your Unity overlay, you may want to use only a few keys on it to begin with. You can hide all the keys that you don't need.

When you hide a key, its LED goes out. If you activate the key, nothing will appear on the display. In order to use the key, you must select **Show** in the **Hide/Show Key** option.

If you choose to hide a key it will no longer work when you activate its location.

AUDITORY PROMPTS MENU		
1	Auditory Prompts <input type="button" value="OFF"/>	4 Key Prompts
2	Prompt Type <input type="button" value="Key"/>	5 Category Prompts
		8 OK
		9 Cancel
		6 Mid-Sequence Prompts

The AUDITORY PROMPT Menu

Auditory Prompts and the AUDITORY PROMPT Menu

Auditory prompts are prompts that you assign to individual keys, to a row or column of keys, or to a quadrant of keys. The Pathfinder will speak the prompts when *Auditory Prompts* is turned On. Spoken prompts help you to learn the location of different keys on the keyboard and what particular sequences contain. You can hear prompts through the speaker in the Pathfinder or you can plug an earphone into the earphone connector on the top of the case.

You can have *Key Prompts*, *Category Prompts* and *Mid-Sequence Prompts*. *Key Prompts* might be icon names, an alphabet letter, a color, etc. *Category Prompts, used for scanning only*, should help to identify what kind of messages a specific row, column or quadrant contains, such as "school" or "family". Or, they might simply identify a specific row, column, etc ("Row 1, Column 1"). *Mid-Sequence Prompts* allow you to change the name of a key in the middle of a sequence so that it relates to that particular sequence only, then reverts back to another (or original) name when the key is used in other sequences.

Prompt Type: Key or Content tells you what you will hear when you activate the last key in a sequence, either the key name or the stored message.

 It is a good idea to use a different voice for the auditory prompt voice. Change voices in the SPEECH menu.

Prompt Type: Key or Content

Prompt Type: Key or Content tells you what you will hear when you activate the last key in a sequence, either the key name or the stored message.

If you set this option to **Key**, the Pathfinder will speak the pre-assigned location code or icon name of each key. You can change pre-assigned names in the AUDITORY PROMPT menu.

If you set this option to **Content**, when you activate the last key in a sequence the Pathfinder will speak the contents of the message stored under that sequence.

➡ Once a prompt is spoken, activate the key *again* to tell the Pathfinder to select it.

⬇ Helpful Hint:

In Unity, some pre-stored sequences contain macros. If you are using **Content** prompts and you activate one of these sequences, the Pathfinder will say, "Macro". To store your own prompt, use the instructions for **Mid-Sequence Prompts** on page 193.

Key Prompts and Category Prompts

➡ **All prompts must be stored using synthesized speech.**

Key prompts are prompts that you assign to individual keys. You can assign a prompt to Minspeak, Spell or Computer.

You may not use this option much as Minspeak and Unity icons already have specific names that are used in pre-stored sequences. One area where a prompt might be helpful, though, is macros. If you activate a key where you have stored a macro (or a tool) the Pathfinder will say, "Macro". You might want to assign a prompt to that key that describes the macro, e.g., "Hide Clock" or "Beep On".

To Assign a Key Prompt:

1. Remove the keyguard and overlay. Turn the toolbox On.
2. Activate the AUDITORY PROMPTS menu.
3. Activate the **Key Prompt** option.
4. On your overlay, activate the key you want to assign the prompt to.
5. Select whether you want the prompt to be spoken when you are in Minspeak, Spell or Computer modes. (You can assign the prompt to all three modes by assigning the same thing to each mode, individually.) Activate OK when you have selected the mode you want.
6. Enter the prompt and activate OK. Continue to complete prompts for each mode. When you have finished, exit the menu.
7. The next time you activate this key, *if* you are in the mode you selected and *if Prompt Type* is set to **Key**, the Pathfinder will speak the prompt you assigned to the key.

Category Prompts are only used if you are scanning. These should be prompts that help you locate where you are on the keyboard, e.g., "Quadrant 1, Row 1 Column 1" or "Row 3 Column E". In some cases, you may have a number of similar messages stored in one row or column

or quadrant. In that case you might want the prompt to say something like, "Row 1, school".

You can assign category prompts to rows, columns and quadrants.

To Assign a Category Prompt:

1. Remove the keyguard and overlay. Turn the Toolbox On.
2. Activate the AUDITORY PROMPTS menu.
3. Activate the Category Prompt option.
4. Select the option you want.
5. Select a key in the appropriate quadrant or row.
6. Select whether you want the prompt to be spoken when you are in Minspeak, Spell or Computer modes. (You can assign the prompt to all three modes by assigning the same thing to each mode, individually.) Activate OK when you have selected the mode you want.
7. Enter the prompt and activate OK. Continue to complete prompts for each mode. When you have finished, exit the menu.
8. The next time you scan the row, column or quadrant, the Pathfinder will speak the prompt you assigned.

To Delete a Prompt:

1. Open the AUDITORY PROMPT menu.
2. Activate the kind of prompt you want to delete.
3. Activate the key where the prompt is stored.
4. Activate the option that has the prompt you want to delete.
5. Use the *Delete Character* key on your overlay to delete the prompt.
6. Activate OK until you exit the menu.

Mid-Sequence Prompts

The default name for the SCALES icon is "Scales". If you do not change the name and have **Prompt Type** set to **Key**, the Pathfinder will say "Scales" when you activate this key.

However, SCALES is also used to signify different categories of words on your Unity overlay. For example, to name different kinds of meat ("Meats" category), you activate the APPLE icon (Food), then the SCALES icon (meat), then a different icon for each type of meat you want to name—if you activate the DOG icon, the Pathfinder will say, "Hot dog". In this case, even when the SCALES icon means "Meat", the auditory prompt will still be "Scales".

The Pathfinder allows you to change the prompt of an icon in a sequence so that it signifies what the icon means in that particular sequence. When the icon is used in a different sequence, the Pathfinder will speak the icon's default name. This is called a **Mid-Sequence Prompt**. You can give any icon in a sequence a Sequence Prompt, *except for the first icon*. Generally, though, if you are using Unity, you will do this with the second icon in a three icon sequence.

To Assign a Mid-Sequence Prompt:

1. Remove the keyguard and overlay. Turn the toolbox On.
2. Activate the AUDITORY PROMPTS menu.
3. Activate the **Sequence Prompt** option.
4. Activate the **Add/Edit Dynamic** Prompt option.
5. Select the key in the sequence where you want to assign the sequence prompt. Activate OK (for example, SCALES).
6. Activate the key(s) leading up to the sequence prompt key. Activate OK (for example, APPLE).
7. Type in the prompt (for example, MEATS) and activate OK. Exit the menu.

You can also edit, delete and review any sequence prompts that you have assigned.

FEEDBACK MENU - Modify settings as needed, then select OK			1		
1	Speaker Volume 6	4	Brightness 5	7	More Items
2	Earphone Volume 10	5	Contrast 7	8	OK
3	Scan Beep ON	6	Activation Feedback Beep	9	Cancel

FEEDBACK MENU - Modify settings as needed, then select OK			2		
1	Auto Spell OFF	4	Text Size Small	7	More Items
2	Spell One Word ON	5	Icon Prediction Medium	8	OK
3	Speak When Spelling Words	6	Icon Tutor On-Screen	9	Cancel

FEEDBACK MENU - Modify settings as needed, then select OK			3		
1	Display Clock DATE and TIME	4	System Colors	7	More Items
2	Display Battery Status ON	5	Resize Text Area NO	8	OK
3	Show Status Bar ON	6	Scan LED Intensity Low	9	Cancel

Feedback Menu

Some Options in the FEEDBACK Menu

Speaker Volume: Change the speaker volume to a louder or softer level. The choices are from 0-20. 10 is the default volume.

Earphone Volume: If you have an earphone plugged into the Pathfinder you can change its volume. The choices are from 1-20. The default is 10.

Scan Beep: The options are On or Off. When this is off, you do not hear a beep as you scan.

Brightness: Choose a number from 1-10 to change the brightness of the display screen.

Contrast: Choose a number from 1-10 to change the contrast on your display.

Activation Feedback: This option lets you choose the kind of feedback you want to hear/see when you activate keys. You can choose to see a light, hear a beep, have both a light and a beep, or not have any feedback at all

Text Size: Select the size of text you want to see on the display, small, medium or large.  See the WORD PREDICTION menu to change the font size of the Word Prediction List.

Icon Prediction: This allows you to turn Icon Prediction On or Off. See page 29 for an explanation of Icon Prediction.

Display Clock: You can choose to display the time, the date, both date and time or to turn the clock display off. (Set up the clock in the MAINTENANCE menu.)

Display Battery Status: This will display battery information in the status bar at the bottom of the display.

Scan LED Intensity: You can change the brightness or dimness of the LEDs you are scanning by choosing a number from 0-7.

Show Status Bar: When this is On, you see a bar at the bottom of your display screen. If you selected to display your clock, it will show up here.

Speak When Spelling: When this is On and Speech is On, you can select to have the Pathfinder speak the letters, the words or the letters and words you type when you are in *Spell*. If you select **Words** you must spell the word and add a space or punctuation mark or quotation mark before the Pathfinder will speak the word.

Resize Text Area: When Word Prediction and the Activity Row are turned On, the display screen is a different size than when they are turned Off. Perhaps you don't use Word Prediction when you are spelling but you always use the Activity Row when you are in Minspeak. You will notice a difference in screen sizes when you move between Minspeak and Spell modes. This option allows you to tell the Pathfinder what to do with the screen size. If you select **No**, the screen size will remain constant whether Word Prediction and/or the Activity Row are turned On or Off. If you select **Yes**, the Pathfinder will resize the screen, always maintaining the largest text area possible.  If you select **Yes**, you will have to be able to re-orient yourself to the new screen size.

Spell One Word

When this option is On, you can go from Spell mode to Minspeak automatically.

If you are in Spell and you decide you want to go to your communication overlay, you can type a word, activate the SPACE key on your spelling overlay and the Pathfinder will automatically switch you to Minspeak.

The default for Spell One Word is Off.

Auto Spell

Summary:

Auto Spell ON allows you to transfer automatically from Minspeak to Spelling. The default is Off.

When you turn Auto-Spell mode On, and activate an empty sequence, the Pathfinder automatically switches you from Minspeak to Spell mode. It also displays the spelling characters that the icons you activated represent.

For example, Auto-Spell is On and you are in Minspeak. Let's say you have a message stored under a sequence that begins with HOUSE, but no message is stored under a sequence that begins with HOUSE APPLE. Let's also say that what you really want to do is spell the word, "Halloween".

Activate the HOUSE icon. The display status bar shows HOUSE>. Now activate the APPLE icon. Normally, the Pathfinder would display, HOUSE APPLE [EMPTY]. However, since Auto-Spell is On, the Pathfinder switches you to Spell mode and displays "ha" on the screen. You can now spell "Halloween".

When you are finished spelling, activate the GO TO MINSPEAK to return to Minspeak.

➡ If **Spell-One-Word** is *also* On, you can activate a space after typing your word and the Pathfinder will automatically switch you back to Minspeak mode. You do not have to activate the GO TO MINSPEAK.

➡ Auto-Spell Tips

Turn Icon Prediction On when you are using Auto Spell. When you activate an icon that is not lit, you know the Pathfinder will switch you to Spell mode. Auto-Spell will not work with Predictive Selection. You must use the GO TO SPELL key to go to Spell mode.

➡ Auto-Spell is very effective when you choose icon sequences wisely. For example, don't store a message under EAR EXIT if you want to use auto-spell to spell words that begin with "ex".

Icon Tutor

Summary:

The Icon Tutor helps you learn the sequences where you have stored specific words.

To use the Icon tutor, you must turn it On in the FEEDBACK menu. The default is Off. Select the number 1, 2, or 3 to tell the Pathfinder how many times to flash the light that will indicate where a word is stored.

Icon Tutor tells you where you have stored individual words.

An Icon Tutor Example:

Perhaps you want to say, "Please," and you can't remember if it is in your Unity vocabulary or not.

First, make sure you are in *Spell*. (If you are in Minspeak, press GO TO SPELL until you are in Spell.)

Spell "Please" on your keyboard. Activate a space key after you have spelled the word. If "please" is stored in your vocabulary, as soon as you activate the space key, the icon or icon sequence that "please" is stored under will be displayed on the screen and the lights for those icons will flash.

When Icon tutor is On, you can find out if any word is in your vocabulary and where it is stored by simply spelling the word and activating a space. If the word is in your vocabulary, the Pathfinder beeps, displays the icon name or sequence on the display and flashes the lights for the icons where the word is stored.

➡ Icon Tutor can only find individual words, not words in a message.

➡ When you *store* individual words, add a space at the end of the word.

System Colors

This option allows you to change the colors that are visible on parts of your display screen.

Changing System Colors:

In the first option, select the system component you want to change the color for. Use the second option to select the background color you want for that component. The third option allows you to select the color you want the text to be. You can see the selections you make in the display box on the right side of the menu. When you have finished, select OK.

Option 1: Change Color For

Menu Title: This is the bar across the top of the display where the menu name and instructions are shown.

Menu Parameter: Parameters are the different options that are listed in particular menu.

Notebook Status: If you have a notebook open, the bar at the bottom of the display will tell you information about it.

Workspace Status: When there is no notebook open, you see the Workspace Status bar at the bottom of the display.

Message Box: These are square "notice" boxes that you will see periodically when the Pathfinder has information to give to you or needs a response from you.

Option 2: Background Color

Step through the color options and select the color you want. You can see an example of your selection in the display box on the right side of the menu (where Option #4 would usually be).

Option 3: Text Color

Step through color options for the text. The Pathfinder will not allow you to select the same color for both background and text. It will automatically change the text color to a default that is visible. You can see an example of your selection in the display box on the right side of the menu (where Option #4 would usually be).

Option 4: Set Default Colors

This returns all the options to their default, or original, colors.

MAINTENANCE MENU			1
1	System Lock <input type="checkbox"/> OFF	4	Set Time/Date
2	Backlight Shut-off <input type="checkbox"/> 5 Minutes	5	Time/Date Format
3	Erase All Vocabulary	6	Delete Imported Icon
		7	More Items
		8	OK
		9	Cancel

MAINTENANCE MENU			2
1	Hardware Tests	4	Error Reporting <input type="checkbox"/> OFF
2	Enter Shipping Mode	5	License Message <input type="checkbox"/> OFF
3	Calibrate Touch Screen		
		7	More Items
		8	OK
		9	Cancel

MAINTENANCE Menu Screens

The MAINTENANCE Menu

The MAINTENANCE menu contains options that you will not use frequently, but which are still important. Use this menu to lock the Pathfinder, set the Backlight to On or Off, erase the Pathfinder's vocabulary and set the system clock. (You can find out the status of the batteries, your software version and other system information by activating the **SYSTEM** key on the left front of the Pathfinder.)

System Lock: When this is On, the use of the I/O and TOOLBOX keys on the front of the Pathfinder case are severely limited. For instance, if the system is locked and you activate the TOOLBOX key, you will see a message on the display that tells you the system is locked. You can press and hold the TOOLBOX key for 5 seconds to activate the Toolbox. This allows you to use *some* of the tools in the Toolbox, but not all of them. **To unlock the system:** open the MAINTENANCE Menu and change the *System Lock* option to *Off*.

Backlight Shut-Off: Set this for 5 or 10 minutes, or Off. When the time you have selected passes and no key activations are made, the Pathfinder will automatically power down into "sleep" mode. This conserves your batteries. To wake up the Pathfinder, activate a key or switch. If you are Headpointing and cannot touch the display or activate a switch to turn the Pathfinder back on, select Off. This means the Pathfinder will never power down until you activate the I/O button on the front of the case. If the Pathfinder stays On all day, remember to charge your batteries regularly.

Erase All Vocabulary: When you select this option, the system erases all its vocabulary and resets itself to the original defaults.

 **Important Note:** When you select the *Erase All Vocabulary* option, the Speech and Word Prediction dictionaries return to their default entries; any customization is lost. You will lose any icons you have imported and all IR commands will be erased.

Set Time/Date: Use this option to set the current time and date for the system clock.

Time/Date Format: Use this option to set the format for the clock readout. You can select either a 12 hour or 24 hour format and a Month/Day/Year or Day/Month/Year format.

Delete Imported Icon: This allows you to delete any imported icons that you no longer want.

Hardware Tests: If you have called the PRC Service Department because something isn't working right with the Pathfinder, a service technician may ask you to open this sub-menu. The technician will walk you through it.

Enter Shipping Mode: if you need to ship the Pathfinder to a different destination, **always** put the device in shipping mode before mailing it.

 Activate *Shipping mode*. Wait ten seconds. After 10 seconds, shipping mode will be in effect and no keys will work.

Calibrate Touch Screen: Use this option if you are pressing keys on the keyboard but not seeing the correct response on the display. Use a stylus or pencil point and follow the directions that you see on the display. When you have calibrated the screen correctly, you will see a screen that tells you to tap the display and wait 30 seconds while the Pathfinder saves the calibrations. If you do not see this screen, repeat the initial calibration.

Error Reporting: When this is On, the Pathfinder reports any system errors that occur. You will see a box on the display that asks you write down what happened and report it to PRC. You are given some options for continuing. When this is Off, the Pathfinder does not report any system errors.

License Message: This option is for PRC's Technical Service Department's use only.

Reset Menu Options to their Default Settings

The RESET DEFAULT tool in the Toolbox allows you to reset defaults (original settings) in menus, option by option.

RESET DEFAULTS is most helpful when writing macros. Tools that toggle between two or more settings are difficult to include in macros without using **RESET DEFAULTS**.

For example, if you want to return any options you have changed in the SPEECH menu to their original settings, you would open the SPEECH menu, activate the RESET DEFAULT key in the Toolbox and then in the SPEECH menu activate the option you want to reset. Repeat these steps for each option you want to reset, in any menu.

To Reset Defaults in Menus:

1. Open the menu that contains the options you want to reset.
2. Activate the RESET DEFAULT tool in the Toolbox.
3. Activate the option you want to reset in the menu that is open.
4. Repeat steps #2 and #3 for each option that you want to return to its original settings.

Macro Examples:

For Loud Speech:

```
<FEEDBACK MENU><RESET DEFAULTS>111118
```

For Quiet Speech:

```
<FEEDBACK MENU><RESET DEFAULTS>1<DOWN>11118
```


Appendices

Appendix A: All Tools List

Summary:

When you activate INSERT TOOL during the storing procedure, the Pathfinder shows you the screens that contain all the tools that you can assign to a key or store under an icon sequence. The tools are listed alphabetically. Explanations of the tools follow the list.

Activate Key	Disable Headpointer
Activity Menu	Down
Alt	Ear Volume Down
Assign Core Key Menu	Ear Volume Up
Auditory Prompt Menu	Edit Core
Auditory Prompts Off	Edit Menu
Auditory Prompts On	Edit Mode
Auditory Prompts On/Off	End
Backlight Intensity Down	Erase Core
Backlight Intensity Up	Exit Theme
Beep Off	Feedback Menu
Beep On	GoTo Computer Mode
Beep On/Off	GoTo Minspeak
Block Marker	GoTo Spell
Calculator	Hide/Show Notebook
Call	Home
Cancel	IBM Marker
Caps Lock	Icon
Caps Off	Insert Tool
Caps On	IR Setup Menu
Choose Activity	IR Command
Clear Display	Left
Close Notebook	MAC Marker
Close Toolbox	Maintenance Menu
Copy Block	Minsert
Ctrl	More Items
Cut Block	Mouse Move
Date Display	Mouse Set
Delete Character	Notebook Menu
Delete Icon	OK
Delete Last Selection	Open Notebook
Delete Word	Open Toolbox
	Output Menu

APPENDIX A: ALL TOOLS LIST

Pathfinder Manual: Basics

Page Down
Page Up
Paste Block
Pick 1 Icon Theme
Pick 2 Icon Theme
Plural
Power Off
Reset Default
Review Vocabulary Menu
Right
Save Changes
Save Notebook
Search/Replace
Selection Technique Menu
Serial Marker
Serial Off
Serial On
Serial On/Off
Set 1 icon theme
Set 2 icon theme
Set Permanent Activity
Set Temporary Activity
Settings Menu
Settings 1
Settings 2
Settings 3
Shift
Sing
Speak All
Speak Date
Speak Display
Speak Last Selection
Speak Next Sentence
Speak Next Word
Speak Notebook
Speak Sentence
Speak Time
Speak Word
Speaker Volume Down
Speaker Volume Up
Speech Marker
Speech Menu
Speech Off

Speech On
Speech On/Off
Spell/Minspeak
Store Core
System Function Menu
Text Marker
Text/Speech Marker
Themes Menu
Transfer Memory Menu
Up
User Area Menu
Voice
Wait Any Key or Switch
Word Prediction Menu
CTRL Any Speak Function
CTRL Block Mark
CTRL Edit Mode
CTRL Serial Marker
CTRL Serial On
CTRL Serial Off
CTRL Speech Off
CTRL Speech On

APPENDIX A: ALL TOOLS LIST

Pathfinder Manual: Basics

Explanations of the Tools in the All Tools List

ACTIVATE KEY

Use this tool to tell the Pathfinder that you want to activate a particular key as part of your macro.

ACTIVITY MENU

Allows you to assign this menu to a key or sequence on your overlay.

ALT (not available at this time; for future use)

When used in conjunction with another tool, and possibly with **CONTROL**, this tool modifies the next character or function.

ALT SPACE

During storing, insert this tool if you want the Pathfinder to speak the alternative pronunciation of the word.

ASSIGN CORE KEY

Allows you to assign a tool or tools and text to a core vocabulary key on your communication overlay. You can also swap keys, hide/show keys, assign tools and/or text to Spell, Computer or Minspeak modes with this tool.

AUDITORY PROMPT MENU

Assign the **AUDITORY PROMPTS MENU** to a key or icon sequence.

AUDITORY PROMPTS OFF

AUDITORY PROMPTS ON

AUDITORY PROMPTS ON/OFF

Turns auditory prompts Off or On or toggles between the two actions.

BACKLIGHT INTENSITY DOWN

BACKLIGHT INTENSITY UP

Turns the display backlighting up or down.

BEEP OFF
BEEP ON
BEEP ON/OFF

Turns the feedback beep on or off or toggles between the two actions.

BLOCK MARKER

Allows you to highlight a block of text.

CALCULATOR

Allows you to store the calculator to a key or sequence.

CALL

When activated, this tool emits a series of tones that will alert people that you need something immediately.

CANCEL

Allows you to cancel an action. You would most likely use this when you are navigating through menus.

CAPS LOCK

Allows you to type using all capital letters. Activate it again to return to lowercase letters. When CAPS LOCK is On, the LED lights next to the picture of the padlock on the upper right front part of the Pathfinder case.

CAPS OFF

Activate to switch to lower case letters.

CAPS ON

Activate to type with all capital letters.

CHOOSE ACTIVITY

Allows you to select an activity from the Activity Row.

CLEAR DISPLAY

Removes all text and prompts from the display and text buffer.

CLOSE NOTEBOOK

Allows you to close an open Notebook from your overlay.

CLOSE TOOLBOX

Allows you to close the Toolbox with a macro.

CONTROL (CTRL)

Activating this key can change what the next key that you activate does. The behavior of some keys will be altered, or the next character typed will generate a "control character" which can be sent to a computer, for example. *See the CTRL commands at the end of this list.*

COPY

Allows you to copy text in one application and paste it into another application. For example, if you have created some text in the regular Workspace that you decide is something that you would like to save in a Notebook, you can highlight the text using BLOCK MARK and copy it to a Notebook.

CUT

Allows you to cut text from a document and paste it somewhere else.

DELETE CHAR

DELETE CHARACTER removes the character to the *left* of the cursor.

DELETE ICON

This tool deletes the last icon that you selected. You can select a new icon or re-select the deleted one if you want.

DELETE LAST SELECT

DELETE LAST SELECTION removes any text that you brought up on the display when you activated your last (previous) selection. If the previous selection was a tool, DELETE LAST SELECTION does nothing.

DELETE WORD

This tool removes the word or group of letters to the *left* of the cursor up to the next space or punctuation mark. If the character to the left of the cursor is a space, DELETE WORD removes the space.

DISABLE HEADPOINTER

Activate this tool to disable your headpointer. When the headpointer is

disabled you will see a box on the Pathfinder screen. When you are ready to re-enable the headpointer, click OK in the box. When the headpointer is disabled you will only be able to activate OK in the box. The headpointer will not activate any other keys.

DOWN

Moves the cursor down one line, one screen, to the end of the workspace, depending on the EDIT MODE setting.

EAR VOLUME DOWN**EAR VOLUME UP**

Turns the earphone volume up or down.

EDIT CORE

Allows you to edit core vocabulary sequences.

EDIT MENU

Allows you to cut, copy edit and replace text from your overlay.

EDIT MODE

Allows you to move through text by sentence, word, or character.

END

Moves the cursor to the end of the current line.

ERASE CORE

Erases sequences stored in your core vocabulary.

EXIT THEME

Allows you to exit a theme from your overlay. Activate once to exit for one sequence; activate twice to permanently exit theme.

FEEDBACK MENU

Allows you to change the feedback the Pathfinder gives you.

GO TO COMPUTER MODE

Allows you to go from the mode you are currently using to Computer mode.

GO TO MINSPEAK

Use this tool to move from the mode you are currently using to Minspeak mode.

GO TO SPELL

This takes you from the mode you are currently using to Spell mode.

HIDE/SHOW NOTEBOOK

This allows you to hide the Notebook you are working in and to retrieve it later.

HOME

Moves the cursor to the beginning of the current line.

IBM MARKER

Use this tool in conjunction with the MAC MARKER tool. Perhaps you use an IBM at work and a Mac at home. Often a Mac and an IBM have a key that does the same thing but it is called something different on each computer (e.g. Mac <DELETE> command or IBM <BACKSPACE> command). When you store computer emulation commands, you can store INSERT TOOL, IBM MARKER, and the specific IBM key command; then add INSERT TOOL, MAC MARKER and the specific Mac key command. Store both these commands under the same icon sequence. When you activate this sequence, the Pathfinder will select the correct command for the computer you are using, (Select the computer you are going to use for computer access in the *Computer Access Menu* of the USER AREA MENU.)

ICON

Places a selected icon into the icon buffer.

INSERT TOOL

Allows you to store tools or text and tools under icon sequences or assign them to keys. Activate this tool before each tool that you want to store. Then activate the tool you want to store.

IR SETUP MENU

Allows you to store infrared commands for remote controls.

IR COMMAND

Use this tool when you are storing remote control commands.

LEFT

Moves the cursor to the left one character, word or sentence, depending on the setting of EDIT MODE.

MAC MARKER

See IBM MARKER.

MAINTENANCE MENU

This menu allows you to load software, erase the Pathfinder, lock the system, set the time and date, set backlight shut-off , etc.

MINSERT

Allows you to put blanks in a message that you can fill when you speak the message.

MORE ITEMS

This allows you to access more items in a menu or activity.

MOUSE MOVE

This allows you to store mouse movements from a list of movements which are relative to your current mouse position, e.g., large movement, very large movement, small movement, etc. The mouse moves in increments.

MOUSE SET

This allows you to store mouse movements from a list of specific movements, e.g., upper left corner, bottom right corner, middle top, etc. The mouse moves to a specific location.

NOTEBOOK MENU

Allows you to place this menu on your overlay.

OK

This is used in creating macros that navigate the menus. This tells the Pathfinder to save what you have done and return to the previous page. The OK key also acts as a "repeat" key in Minspeak mode. After you

have retrieved an icon sequence, you can activate OK to repeat the sequence.

OPEN NOTEBOOK

Allows you to open a specific Notebook from your overlay.

OPEN TOOLBOX

Allows you to open the Toolbox from your overlay.

OUTPUT MENU

This menu allows you to set up the Pathfinder so that it will work with a computer or printer, remote control or memory transfer.

PASTE

Allows you to paste cut or copied text into a new place.

PAGE DOWN

Moves the text on the display down a page.

PAGE UP

Moves the text on the display up a page.

PICK 1 ICON THEME

Switches you to a 1 icon theme and allows you to select the theme icon.

PICK 2 ICON THEME

Switches you to a 2 icon theme and allows you to select the first two icons.

PLURAL

When activated, the word to the left of the cursor is made plural.

POWER OFF

When this tool is activated the Pathfinder 'goes to sleep.' Activate a key or switch to wake it back up.

RESET DEFAULT

This resets all the Pathfinder defaults.

REVIEW VOCABULARY MENU

Allows you to review all the vocabulary in the Pathfinder. You can print it out or see it on the display.

RIGHT

Moves the cursor right one character, word or sentence depending on the setting of EDIT MODE.

SAVE CHANGES

Normally, when you turn the Pathfinder off, you see "Saving changes..." on the screen. This tool allows you to save changes without turning the Pathfinder off. You might want to store this tool if you are doing a lot of programming and want to periodically save your work.

SAVE NOTEBOOK

Allows you to save the current Notebook from your overlay.

SEARCH AND REPLACE

Allows you to search for a word and replace it.

SELECTION TECHNIQUE MENU

Choose a selection technique and the options for it that you want from this menu.

SERIAL MARKER

Any text that follows a SERIAL MARKER is sent out the serial port but not displayed or spoken.

SERIAL OFF

Turns the serial output Off.

SERIAL ON

Sets the serial output to its last active setting, either Immediate or Delayed. Text is automatically sent out the serial port when this is turned On.

SERIAL ON/OFF

Toggles between Serial On and Serial Off.

SET 1 ICON THEME

Switches you to a 1 icon theme you stored previously.

SET 2 ICON THEME

Switches you to a 2 icon theme you stored previously.

SET ACTIVITY

Allows you to go to a specific activity. Use the arrow keys and OK key on your overlay to navigate through your activities if you don't want to use the display screen.

SETTINGS MENU

Allows you to save settings and change between settings without navigating each of the menus. Use this with SETTINGS 1, 2 and 3 below.

SETTINGS 1**SETTINGS 2****SETTINGS 3**

Allows you to store and then access 3 different groups of settings.

SHIFT

When this is set, the next character you activate will be an UPPERCASE letter. The second character you activate will be in lowercase. When you activate SHIFT, the "Shift" LED flashes on the upper right front of the Pathfinder case.

SING

This tool shows you a list of all the songs in your Pathfinder and allows you to pick a song.

SPEAK ALL

When this tool is activated, the Pathfinder speaks all of the text in the text area of the display and in the buffer. (When you type, letters and words will move off the screen in order to make room for new letters and words. Words that are no longer on the screen have moved to the workspace. You can use the arrow keys to bring them back onto the display screen.) Clearing the display (CLEAR DISPLAY tool) or turning the device Off clears the buffer. This tool works whether speech is On or Off.

SPEAK DATE

When you store this tool and then activate it, the Pathfinder will speak the current date.

SPEAK DISPLAY

When this tool is activated, the Pathfinder speaks everything that is visible in the text display area. This tool works whether speech is On or Off.

SPEAK LAST SELECTION

When you activate this tool, the Pathfinder speaks the last stored message you activated or the last spelling characters you put on the display. This is helpful if you have to repeat yourself so that someone listening can understand you better, or to get someone's attention. This tool works whether speech is On or Off.

SPEAK NEXT SENTENCE

For use with Notebooks. If you want to speak a Notebook sentence by sentence, use this tool. Activate it each time you are ready to move to the next sentence.

SPEAK NEXT WORD

For use with Notebooks. If you want to speak a Notebook word by word, use this tool. Activate it each time you are ready to go on to the next word.

SPEAK NOTEBOOK

Tells the Pathfinder to speak the Notebook you stored with the tool.

SPEAK SENTENCE

When you activate this key, the Pathfinder speaks the message that contains the cursor. This is helpful if you want to hear how something you are working on sounds before you store it. It is also useful if you are going to give a speech or presentation. People seem to understand DECtalk Access32 speech better when they hear it one sentence at a time until they are used to listening to it. This tool works whether speech is On or Off.

SPEAK TIME

When you store this tool and then activate it, the Pathfinder will speak the current time.

SPEAK WORD

When this tool is activated, the Pathfinder speaks the last word you typed on your spelling overlay, or the last word you brought up in Minspeak mode. If you are typing in Spell mode, you can turn speech OFF (so you don't have to listen to the Pathfinder speak each letter) and type the word. Then you can activate this key and the word you typed will be spoken. See also the FEEDBACK menu, *Speak When Spelling*.

SPEAKER VOLUME DOWN**SPEAKER VOLUME UP**

Turns the speaking volume up or down.

SPEECH MARKER

This marks the beginning of text in a sequence that is to be spoken but not displayed.

SPEECH MENU

This allows you to insert the SPEECH MENU in a macro.

SPEECH OFF**SPEECH ON****SPEECH ON/OFF**

When this tool is On, all the sequences you select or letters you type are automatically spoken. When this tool is Off, nothing is spoken unless you activate one of the SPEAK tools. SPEECH ON/OFF toggles between the two actions. (If *Speak When Spelling* is On in the FEEDBACK menu, the Pathfinder will speak regardless of where this tool is set.)

SPELL/MINSPEAK

Toggles between Spell mode and Minspeak mode.

START SCAN

Add this to a macro to force the Pathfinder to start scanning each time, regardless of how long a switch is held down.

STORE CORE

Allows you to store core vocabulary messages under icon sequences.

SYSTEM FUNCTION MENU

This single menu gives you access to all other menus. With this menu, you do not have to go to the Toolbox to access menus, nor do you have to store each individual menu in an icon sequence or on your overlay.

TEXT MARKER

This marks the beginning of text in a sequence that is to be displayed but not spoken.

TEXT AND SPEECH MARKER

This marks the end of the text in a sequence where a SPEECH MARKER and/or TEXT MARKER has been used. Any text after the TEXT AND SPEECH MARKER will be both displayed and spoken.

THEMES MENU

Allows you to have this menu on your overlay.

TRANSFER MEMORY MENU

This allows you to save the Pathfinder's memory to a computer disk, transfer memory from one Pathfinder to another and import graphics, and Save/Load Notebooks to and from a computer. (See the *MTI Manual* for more information.)

UP

Moves the cursor up one line, paragraph, to the end of the workspace, depending on how EDIT MODE is set.

USER AREA MENU

Allow you to select user areas and overlays.

VOICE

Allows you to use different voices in a message, without having to store the whole SPEECH menu in a macro. This gives you only the voice options from the menu.

WAIT ANY KEY OR SWITCH

Causes a pause in a macro until a key is activated.

WORD PREDICTION MENU

This menu allows you to turn Word Prediction On and Off, add words to the word prediction dictionary, delete words, etc.

CTRL EDIT MODE

Clears the current edit mode.

CTRL BLOCK MARK

Clears the current block marker.

CTRL SERIAL MARKER

This tells the Pathfinder to send data out the serial port only if the Serial option is turned On.

CTRL SERIAL ON

Turns the serial output back on in whatever the last serial mode was, for the remainder of the macro in which it is found.

CTRL SERIAL OFF

Turns the serial output back off in whatever the last active serial mode was, for the remainder of the macro in which this is found.

CTRL SPEAK ANY FUNCTION

Performs any regular speech function (SPEAK ALL, SPEAK DISPLAY, SPEAK NOTEBOOK, etc.) only if speech is turned on.

CTRL SPEECH ON

Turns speech on for the remainder of the macro in which it is found.

CTRL SPEECH OFF

Turns the speech off for the remainder of the macro in which this is found.

CTRL START-SCAN

This only starts scanning if *Scan Auto Restart* is set to ON (in the SELECTION TECHNIQUE MENU)

Appendix B: About Data Logging and Language Activity Monitoring Using the LAM

A significant milestone in the field of augmentative and alternative communication (AAC) is the development of the notion of language activity monitoring. LAM is the intentional monitoring of the language activity of people who rely on AAC systems. LAM has implications in three primary areas of AAC: clinical intervention, outcomes measurement, and research. (Romich and Hill, 1999)

LAM is being implemented in high performance AAC systems. (For older AAC systems, an external LAM device, available from PRC, can be connected to the serial port.) The LAM records each language event (one or more letters or words) and adds a time stamp (24-hour format with a resolution of one second). This data provides a record of what was said and when it was said. Time information allows determination of the language representation method(s) being used and the communication rate being achieved.

The LAM data can be uploaded periodically to a computer for editing and analysis. The use of a standardized reporting protocol provides for universal compatibility with the computer-based tools. (Hill and Romich, 1999)

The practice of any discipline can be either an art or a science. Science requires data. LAM data is changing AAC from an art to a science.

References

- Romich, B. and Hill, K., "A Language Activity Monitor for AAC and Writing Systems: Clinical Intervention, Outcomes Measurement, and Research", *Proceedings of RESNA '99*, pp. 19-21, 1999.
- Hill, K., and Romich, B., "A Proposed Standard for AAC and Writing System Data Logging for Clinical Intervention, Outcomes Measurement, and Research", *Proceedings of RESNA '99*, pp. 22-24, 1999.

Instructions for Using the LAM

For the instruction on the following pages you will need: an IBM PC or clone running Windows 95 or 98 with the LAM Terminal program; a Pathfinder with LAM software installed, and a null modem cable to connect the communication device to the computer.

➡ The LAM Terminal Program is on a floppy disk in an envelope at the back of the LAM *Procedural Manual*.

These LAM instructions will teach you to:

- Setup the LAM terminal program on your Windows Desktop
- Setup the Pathfinder for the upload
- Perform a LAM upload so that you can do it easily at any time.

⬇ Before You Begin:

In order to accumulate language data to be analyzed, you must turn *Data Logging On* in the REVIEW VOCABULARY MENU **prior** to the time you want to analyze the data. The LAM software will begin to record all language activity data from the time it is turned On. You will upload this data to a computer running the LAM Terminal program at a later date.

To check if there is any data to be transferred and analyzed, activate the STATUS key on the left front of the Pathfinder. Check the *Data in Log* line. If it says *Data Logging = xx bytes*, there is information you can upload. If it says *Data Logging = OFF 0 bytes* there is no information to be uploaded. You must turn *Data Logging On* (in the REVIEW VOCABULARY MENU) and wait to perform an upload until the Pathfinder has been used enough to accumulate language data.

Set Up Instructions

On both the Computer and the Pathfinder

1. Connect your computer and the Pathfinder with a null modem serial cable. You can use the MTI cable that came with your Pathfinder.
2. Plug one end of the cable into COM 1 on your computer.
3. Plug the other end of the cable into the serial port on the Pathfinder. The serial port is on the top of the case and has a picture of a computer above it. If your computer has more than one port, and COM 1 is being used for something else, plug the cable into another port.

On the Computer

 There is a floppy disk in an envelope at the back of the *LAM Procedural Manual*.

1. Remove the disk from the envelope and put it in the floppy drive on your computer.
2. On your computer, go to the drive where you put the floppy disk and open the disk.
3. Click on and drag the **Lamterm.exe** icon to your desktop.
4. Double click on the **Lamterm.exe** icon on your desktop.
5. You see the **Lamterm** screen.

The **Status** box indicates the baud rate and the computer COM port that are currently selected. Initially these are set to the factory defaults, using COM port 1 and a 9600 baud rate.

If you want to change COM ports or baud rates:

Click on the *Options* menu.

Click on *Port Settings*.

Click in the box labeled *Use Defaults*. This removes the check mark in the box.

Click on the *Port* and/or the *Baud* rate you want to use. (Make sure the LAM box is connected to the port you have selected.)

Click on *OK*.

On the Pathfinder

1. Open the Toolbox.
2. Open the OUTPUT menu.
 - Set *Baud Rate* to **9600**. The OUTPUT menu baud rate and the *Lamterm* baud rate must be the same.
 - Set *Serial* to *On*.
 - Set *Auto Line Feed* to *On*.
 - Activate *OK* to exit.

To Confirm that Everything is Properly Set Up between your Computer and the Pathfinder:

On the Pathfinder, open the REVIEW VOCABULARY MENU in the Toolbox.

Activate *Send Data Log to Computer*.

Look at the *Number of Bytes Received* line in the Lamterm **Status** box on your computer. You will see the data transferring in byte numbers.

When the data has stopped transferring:

You see a *Save As* screen on your computer.

You see a suggested format for saving the data

YYMMDD_n_ls_name

- **YYMMDD = the year, month and date;**
- n = the number of the language sample for this person on this date
- ls = the Language Sample context (2 characters)
(for example, pd=picture description; in=interview;
co=conversation; etc.)
- name = the first four letters of the person's last name.

In the ***Save In*** box at the top of the screen, select where you want to save the upload.

 The upload will be saved as a Text file (.txt). You will have to re-save the file as a Microsoft Word file or a WordPad file.

To Save the Upload file in MSW or WordPad:

1. Open Microsoft Word or WordPad.
2. Open the file you recently uploaded.
3. Save the upload again as a MSW or Word Pad file. Delete the original .txt file.

Using the LAM Upload

Follow these instructions each time you want to perform an upload.

1. Connect the Pathfinder to a computer using a null modem cable (the '95 MTI cable will work).
2. Double-click on the *Lamterm.exe* icon on your computer desktop.
3. On the Pathfinder, open the REVIEW VOCABULARY MENU in the Toolbox. Activate *Send Data Log to Computer*.

When the data has all been transferred, save it on your computer using the Lamterm *Save As* screen that appears on your computer screen

1. Name the file.
2. Select where you want the file to be saved.
3. Select **Save**.

When the file is saved:

1. Open your word processing application (MSW or WordPad).
2. Open the upload file you just saved.
3. Save it as a file in MSW or WordPad.
4. Delete the original .txt file

On the Pathfinder

1. Open the REVIEW VOCABULARY MENU in the Toolbox.
2. Activate *Clear Data Log*.
3. Activate OK to exit the menu.

 If you want to accumulate more date for a later upload, make sure *Data Logging* is *On* in the REVIEW VOCABULARY MENU.

About the LAM Data

When you look at the raw LAM data that you saved in your computer file or folder, the first thing you will see is the Time Stamp. This is the first item in the first line. The time is presented in a 24-hour format.

Then you see the SOURCE of the language event. The following SOURCE codes are currently valid: ACT (Activity Row); SEM (Semantic compaction = Minspeak); SPE (spelling); WPR (word prediction); DWP (Derivative from Word Prediction). This code is in capital letters.

The character string that follows the Source code is the language event that occurred at that time.

Any characters spaced more than 0.2 seconds apart are considered to be separate language events.

Below is an example of how the raw LAM data looks:

```
20:37:00 SEM "I need "  
20:37:05"*[VOLUME UP]*"  
20:37:06"*[VOLUME UP]*"  
20:37:07"*[VOLUME UP]*"  
20:37:14 SEM "something "  
20:37:16 SEM "to drink "  
20:37:19 SPE "i"  
20:37:20 SPE "m"  
20:37:24 SPE "m"  
20:37:28 WPR "ediately "
```

From this example you can see that whole words and phrases were accessed, in this case using the semantic compaction language representation method. The "[VOLUME UP]" event indicates that the system user is capable of identifying the need to perform that function

and can implement it. The word "immediately" was selected using spelling and word prediction. The utterance contains six words. Including the volume adjustment, the time from the first event to the last event was 28 seconds.

LAM data can be analyzed in several different ways. Using a word processor, you can paste the LAM data into a new document and save it in its raw form. Some word processors, such as Word for Windows, allow you to create a macro that can be used for some editing, such as removal of the Time Stamp. However, simply viewing the raw data can provide information about the vocabulary that is being selected and how the AAC system is being used. For example, if a particular word was introduced or practiced at a previous therapy session, a search for that word in the data can indicate how many times the word was used in later conversations.

Other performance monitoring tools soon to be available from Prentke Romich Company include an editor, language analysis program(s) and implementation programs.

See the *Procedure Manual for the LAM* for more information.

Appendix C: How to Use the *Review to Output* Option in Various Pathfinder Menus

You can review the Pathfinder vocabulary, the Word Prediction Dictionary or the Speech Dictionary by looking at the review on a computer or printing it out on a printer.

To Review to an IBM or clone Computer:

1. Connect your computer and the Pathfinder with a null modem serial cable. You can use the 95 MTI cable that came with your Pathfinder. Plug one end of the cable into COM 1 on your computer. If your computer has more than one COM port and COM 1 is being used for something else, plug the cable into another port. Plug the other end of the cable into the serial port with the computer icon on the top of the Pathfinder.
2. Select the computer *Start* menu.
 1. Select *Programs*.
 2. Select *Accessories*.
 3. Select *HyperTerminal*.
 4. In the HyperTerminal window, select *Hyperterm.exe*.
 5. In the *New Connection* dialogue box:
 - Type, "Review" as the name.
 - Select an icon.
 - Press *OK*.
 6. In the *Phone Number* dialogue box:
 - At *Connect Using*, select *COM 1*. (If you connected your cable to a different COM port select the correct port.)
 - Press *OK*.
 7. In the *COM 1 Properties* dialogue box:
 - Either select the *Defaults* option, or set the following:
 - Data Bits = 8*
 - Parity = None*

Stop Bits = 1

At ***Bits per Second***, select ***9600*** if this is not already set there.

At ***Flow Control***, select ***None***.

8. The "Review" window should now be open.
11. Open the ***Transfer*** menu in the menu bar. Select ***Capture Text***. Type in the ***drive*** and the ***file or folder*** where you want the text to go. Click on ***Start***.

On the Pathfinder

1. Open the Toolbox on the Pathfinder.
2. Open the OUTPUT menu.
 - Set ***Output Method*** to ***RS232***.
 - Set ***Output*** to ***On***.
 - Set ***Auto Line Feed*** to ***On***.
 - Set ***Output Mode*** to ***Immediate***.
 - Set ***Baud Rate*** = ***9600***. (The HyperTerminal and Pathfinder baud rates must be the same.)
 - Activate ***OK*** to exit.
3. From the menu you are using in the Pathfinder, activate ***Review to Output***.

When the review has been transferred to your computer you can ***Copy*** it into a word processing program, ***Save*** it in a file or folder and then use the computer ***Print*** function to print the review.

On a Macintosh Computer

1. Connect your computer and the Pathfinder with the MTI Mac cable that came with your Pathfinder. Plug one end of the cable into either the Printer port or the Modem port of your computer, depending on what your Terminal program calls for. (Some Terminal programs allow you to connect to either port, some only allow you to connect to the modem port.)
2. Open the Terminal program on your computer.

3. Select the port that you are connected to.
4. Set the baud rate to match the Pathfinder's baud rate. (The Pathfinder baud rate is in the OUTPUT menu.)

On the Pathfinder

1. Open the OUTPUT menu.
 - Set **Output Method** to **RS232**.
 - Set **Output** to **On**.
 - Set **Baud Rate** = **9600**. (You can select a higher baud rate, but the terminal program and Pathfinder baud rates must be the same.)
 - Set **Auto Line Feed** to **On**.
 - Activate **OK** to exit.
3. From the menu you are using in the Pathfinder, activate **Review to Output..**

When the review has been transferred to your computer you can **Copy** it into a word processing program, **Save** it in a file or folder and then use the computer **Print** function to print the review.

To Print Directly to a Serial Printer

1. Attach the printer cable to the **serial port** with the picture of a computer on the top of the Pathfinder.
2. Set the printer baud rate to match the baud rate in the Pathfinder's OUTPUT menu.
3. The printer should be set for **8 Data Bits**, **2 Stop Bits** and **Parity** should be set to **None**.
4. Open the OUTPUT menu in the Pathfinder Toolbox.
 - Set **Output Method** to **RS232**.
 - Set **Output** to **On**.
 - Set **Baud Rate** = **9600**. (You can select a higher baud rate, but the printer and Pathfinder baud rates must be the same.)
 - Set **Output Mode** to **Immediate**.

5. Make sure the printer is on-line and ready to print.
6. From the menu you are using in the Pathfinder, activate ***Review to Output***.

To Print to a Parallel Printer

If you want to print to a parallel printer, you will need a parallel-to-serial adapter. You can usually get these adapters from the place you got your printer.

1. Plug the printer cable into the adapter and plug the adapter into the serial port on the top of the Pathfinder.
2. Set the baud rate on the printer to match the Pathfinder's baud rate.
3. The printer should be set for ***8 Data Bits***, ***2 Stop Bits*** and ***Parity*** should be set to ***None***.
4. In the Pathfinder OUTPUT Menu
Set ***Output Method*** to ***RS232***.
Set ***Output*** to ***On***.
Set ***Baud Rate*** = ***9600***. (You can select a higher baud rate, but the printer and Pathfinder baud rates must be the same.)
Set ***Output Mode*** to ***Immediate***.
5. Make sure the printer is on-line and ready to print.
6. From the menu you are using in the Pathfinder that contains the ***Review to Output*** option, activate ***Printer***.

Appendix D: The Pre-Stored Songs in the Pathfinder

The following songs are pre-stored in the Pathfinder for your use.

Amazing Grace
Joy to the World
O Come Little Children
Silent Night
Away in a Manger
My Country 'Tis of Thee
Old MacDonald
Mary Had A Little Lamb
Humpty Dumpty
Twinkle Twinkle Little Star
London Bridge
I'm a Little Teapot
Hickory Dickory Dock
Jack and Jill
Jingle Bells
For He's a Jolly good Fellow
Michael, Row the Boat Ashore
Ring Around the Mulberry Bush
Polly Put the Kettle On
Waltzing Matilda
Alouette
Auld Lang Syne
My Bonnie Lies Over the Ocean
The Farmer in the Dell

Appendix D: To Copy a Pre-Stored Song to Your Overlay

The songs on the previous page are already stored in the Pathfinder. In order to sing them, you have to move them from the Toolbox to your Minspeak overlay or into an Activity.

To Move a Song to Your Overlay:

1. Go to the Toolbox.
2. Activate STORE CORE.
3. Select the icon sequence you want to use.
4. Select "Spell text for message."
4. Activate INSERT TOOL.
5. Activate SING.
6. Select the song you want to copy.
7. Activate OK twice.

The next few pages explain how to program songs that you want the Pathfinder to sing.

Appendix E: How to Program and Store Your Own Songs in the Pathfinder

 Remember: You can always record songs rather than programming them

Programming a song in DECTalk speech is different from programming regular text because you must tell the Pathfinder to *sing* the song, not speak it.

To do this you must type each letter of text phonetically and add time and pitch values, which correspond, to musical notes. At the end of this appendix you will find a list of PITCH and TIME values, a chart of musical notes and a list of vowel and consonant phonemes to use.

Programming songs is not difficult but it *is* time consuming, so plan accordingly.



Highly Recommended

Write down the programming of your song *before* you begin to program it into the Pathfinder.

This means writing down the words to the song in phonemes and adding the time and pitch values exactly the way you want to program them. This will save you a great deal of time when you are ready to enter the song into the device.



Helpful Hint

You will be doing quite a bit of cross-referencing with the lists at the end of this appendix. We encourage you to use the sheet music for any song you want to program. This way you can match the notes on the sheet music with the notes in the pitch and time tables. You don't have to be able to read music.

When programming a song you must change the letters in your song to phonemes (the written form of sounds). Find the sound in the Vowel and Consonant Phoneme List for each letter in your song and write down the phoneme.

Once you have entered the correct phoneme for a letter, you must tell the Pathfinder how *long* (time) to hold which *note* (pitch). The TIME and PITCH values are *always* enclosed in angle brackets < > immediately following the phoneme they represent. The TIME value *always* comes before the PITCH value inside the brackets. The two values are *always* separated by a comma, i.e., <300,24>.

You always start to program a song by activating the left square bracket [. This tells the Pathfinder that you are programming a song. Program all songs in STORE mode for DECTalk speech.

Break the Song into Phrases or Lines

Most songs, unless they are very, very short (one line) must be broken down into separate phrases or single lines to be stored. After each phrase (or line) you must activate the *right* square bracket,] and then the **CTRL** key and the "k" key on your keyboard. You will see the right bracket and two small characters beside it on your display.

The right square bracket] usually tells the Pathfinder that it has come to the end of the song. When it is combined with the **CTRLk** activation, it tells the Pathfinder to keep on singing the song until it comes to a right square bracket that is **not** followed by a CTRL k character. The] *without* a CTRL k following it signifies the end of the song.

➡ Before you Begin to Store:

Write down the song you want to store *exactly* the way you will enter it into the Pathfinder.

1. Start with the left square bracket [.
2. Use the **phoneme lists** to change each letter of the song into the correct phoneme.
3. Use the **TIME and PITCH tables** to match the note values in your song.
4. Write down the TIME and PITCH values after each phoneme.
Enclose the TIME and PITCH values in angle brackets < > and separate the two values with a comma.
5. Break the song into phrases or lines by adding a **right square bracket] at the end of the line or phrase**. Use a reminder of some kind to tell you to activate the **CTRL and "k" keys** when you are actually programming.
6. Start the next phrase or line with the left square bracket [.
7. Continue until you have completed the song.

About TIME Values

The TIME value refers to the *amount of time in milliseconds* that a note is to be held.

1,000 milliseconds = 1 second;

500 milliseconds = ½ second.

You can use the above numbers as a reference point.

You choose the amount of time you want a note to be held according to the tempo of the song. For example, a *quarter note* is held half as long as a *half note*, so you would adjust the time value for a phoneme accordingly. Use the **Note Time Value Chart** at the end of the appendix.

After you have entered a phoneme, type the left angle bracket < and then enter the TIME value (how long you want that phoneme to be held).

About PITCH Values

PITCH values correspond to notes as you see them on sheet music or hear them on a piano. The PITCH value of notes for the Pathfinder range from 1 to 37. 1 corresponds to C2 (or the second C below middle C on a piano). 25 corresponds to the note C4 (middle C on a piano). 37 corresponds to C5 (the first C above middle C on a piano).

The two tables, **PITCH NOTES AND VALUES** and **PITCH VALUE CODE**, at the end of the appendix will help you to figure out PITCH values.


Enter the PITCH value that corresponds to the note you want after you have entered the TIME value and a comma. Then type a right angle bracket > . Do not put a space before or after the comma.

An Example of Programming a Song

Mary Had a Little Lamb



Ma - ry had a lit - tle lamb lit-tle lamb lit-tle lamb

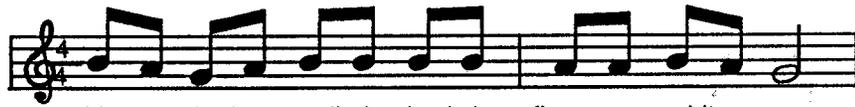

Ma - ry had a lit-tle lamb its fleece was white as snow

#1 Change the words of the song into phonemes. Use the Vowel and Consonant Phoneme Tables to find the correct phonemes.

Mar y ha d a li t tle la mb
MER IY HXAE D AX LIH T TUHL LEYAE M

li t tle la mb li t tle la mb
LIH T TUHL LEYAE M LIH T TUHL LEYAE M

mar y ha d a li t tle la mb
MER IY HXAE D AX LIH T TUHL LEYAE M

I t s fle e ce wa s whi te a s snow
IH T S FLIY S WAX Z WAY T AE Z SNOW.

If you are new to using phonemes, it will take a little practice to become familiar with the different sounds. Sometimes a beginning consonant will blend into a vowel, but other times it may have to stand alone.

Most ending consonants will stand alone but will be held for a very short time. Use the phoneme tables at the end of the appendix.

Practice by spelling words into the Pathfinder using the phonemes you think might work. Then listen to them. Soon you will develop an idea about which phonemes are going to correspond to which sounds.

#2 Find the correct **PITCH** value in the **PITCH Values** charts.

Below is the regular notation for *Mary had a Little Lamb*.

The image shows two staves of musical notation for the song 'Mary had a Little Lamb'. The first staff is in 4/4 time and contains the melody for the first line of the song. The notes are: Ma (B), ry (A), had (G), a (A), lit- (B), tle (A), la (A), mb (B), lit- (A), tle (B), lamb (D). The second staff is also in 4/4 time and contains the melody for the second line of the song. The notes are: Ma (B), ry (A), had (G), a (A), lit- (B), tle (A), lam (A), b (A), its (A), fleece (A), was (A), white (A), as (A), snow (A). The pitch letters are placed above the notes.



How You Can Make the Song

Sound Better:

Most songs that have notes above "37" (see the PITCH Value Code Chart) are scored too high for any Pathfinder voice to sing well. These higher songs will sound much better if you give the notes their value in an octave lower. For example, the first note in *Mary Had a Little Lamb* is B above middle C. This B has a PITCH value of 36. To assign the correct value in an octave lower, go to the PITCH Value Chart and find

the first B *below* middle C. Notice that its value is 24. Assign 24 to your B.

Also, some songs sung in a male Pathfinder voice may sound better if you program them in an octave lower than they are written.

Below, the PITCH values for *Mary had a Little Lamb* are shown an octave lower than written in the music.

The image shows two staves of musical notation in 4/4 time. The first staff contains the melody for the first line of the song, with pitch values written above the notes: 24, 22, 20, 22, 24, 22, 24, 27. The second staff contains the melody for the second line, with pitch values written below the notes: 24, 22, 20, 22, 24, 22, 24, 22, 20. The lyrics are written below the notes.

Ma - ry had a lit - tle lamb lit-tle lamb lit-tle lamb
 Ma - ry had a lit-tle lamb its fleece was white as snow

#3 Find the correct **TIME** value for each note and write it down.

1. Using the secondhand on a watch, time the amount of time it takes you to sing your song.
2. On the sheet music, look at the tempo of the song. In our example it is 4/4 time.
3. The top "4" means that there are four beats to a measure. There are four measures in our song. The vertical lines between groups of notes define a measure.
4. The bottom "4" means that each $\frac{1}{4}$ note gets one beat.
5. Suppose it takes 10 seconds to sing *Mary had a Little Lamb*.

To find out how many seconds each measure gets, divide 4 into 10.
In this case each measure gets 2.5 seconds.

6. To find out how many milliseconds each $\frac{1}{4}$ note (one beat) gets, divide 4 into 2.5. In this case each $\frac{1}{4}$ note gets 625 milliseconds per beat. (You can round these numbers up or down so that they are easier to work with; e.g., 600.)
7. An $\frac{1}{8}$ note is held for one-half the time of a $\frac{1}{4}$ quarter note, or 312 milliseconds (624 divided by 2). A $\frac{1}{16}$ note is held for one-half as long as an $\frac{1}{8}$ note, or 156 milliseconds (312 divided by 2).
8. A $\frac{1}{2}$ note is held twice as long as a $\frac{1}{4}$ note, or 1250 milliseconds (625 times 2).
9. A whole note ($\frac{4}{4}$) is held twice as long as a $\frac{1}{2}$ note, or 2500 milliseconds (1250 times 2).

This gives you a rough formula for figuring out the TIME values.

10. In our example, "Mary" (two $\frac{1}{8}$ notes in the music) makes up one beat and a beat is one $\frac{1}{4}$ note. But "Mary" is broken down into three phonemes: m er iy (the "a" and "r" in Mary is one vowel sound). The "m" "ea" "r" and "iy" phonemes are equivalent to one $\frac{1}{8}$ note. According to our calculations above, an $\frac{1}{8}$ note gets 312 milliseconds.



Shortcut for Finding a TIME Value

If you are not concerned with achieving the exact tempo of a song, there is a quicker way to figure out the TIME values.

Hum a song. Look at a watch or have someone count to calculate the rough timing of the song. Decide how many seconds you want a whole note to be held. Then divide the $\frac{1}{2}$, $\frac{1}{4}$, and $\frac{1}{8}$ notes accordingly.

For example: If you give a whole note the time of 2 seconds (or 2000 milliseconds), then a $\frac{1}{2}$ note gets the value of 1000 milliseconds; a $\frac{1}{4}$ note gets the value of 500 milliseconds; an $\frac{1}{8}$ note gets the value of 250 milliseconds.

By practicing giving whole notes different time values and listening to the results on the Pathfinder, you can approximate the tempo of a song fairly accurately.

#4 Assign TIME and PITCH values to the phonemes.

Now let's figure out how many milliseconds from among the 312 (or 625 for a $\frac{1}{4}$ note) to give to each phoneme.

Vowels are always held for a longer amount of time than consonants. Consonants that stand alone at the beginning and end of words are usually held for short amounts of time.

Now that we know this we can begin to put the phonemes, TIME value and PITCH value together.

The "m" and "er" of "Mary" can probably be sung together—"MER". Together, these two phonemes make up an $\frac{1}{8}$ note. We discovered that

in this song an 1/8 note uses 312 milliseconds, so let's assign a TIME value of 300 to "MER". (300 is easier to work with than 312 and there is no noticeable difference between the length of 300 and 312.)

When you compare the notes on your sheet music with the notes on the PITCH Charts, you see that the first note (MER) has a PITCH value of 24. So "MER" will be sung with a TIME value of 300 and a PITCH value of 24.

The "IY" phoneme is held for about the same amount of time as "MER" when you sing it so let's give it a TIME value of 300, too. When we look at the PITCH charts, we see that its PITCH value is 22.

The first word of our song will look like this when you enter it into the Pathfinder:

```
[MER<300,24>IY<300,22>
```

The second word, "HAD" uses an 1/8 note and contains the phonemes "HXAE D". The "D" is an ending consonant and stands alone. We don't want it to sound too pronounced when the word is sung, so it should get a lesser TIME value than "HXAE".

Now we have:

```
[MER<300,24>IY<300,22>HXAE<220,20>D<80,20>
```

Let's give "HXAE" a TIME value of 220. Its PITCH value is 20. And let's give "D" a TIME value of 80. Its PITCH value is 20. (220 plus 80 equals 300, or an 1/8 note.)

"A" (phonemes "AX") is easy. It's a one-phoneme 1/8 note. It gets a TIME value of 300 and the assigned PITCH of 22.

"LITTLE" must be broken up into "LIH T TUHL". Because we want to hear the "T" pronounced in "little" we should give it its own TIME value. "LIH" and "T" make up the first 1/8 note, so let's give LIH the TIME value of 250 along with its PITCH value of 24. Since an 1/8 note has 300 milliseconds and we have used 250 for "LIH", that leaves us with a TIME value of 50 for "T" (300 minus 250 equals 50). Its PITCH value is 24. "TUHL" uses an 1/8 note so it gets a TIME value of 300 and its PITCH value of 24.

Now we have

```
[MER<300,24>IY<300,22>HXAE<220,20>D<80,20>AX<300,22>
LIH<250,24>T<50,24>TUHL<300,24>"]
```

LAMB" ("LEYAE M") uses a ¼ note or about 600 milliseconds. It has an ending consonant, "M" (the "b" is silent) which must stand alone so we can hear it. Let's give "LEYAE" a TIME value of 500. Its PITCH is 24. ("LEYAE" sounds better when sung than "LAE".) "M" then gets a TIME value of 100. Its PITCH value is 24.

When "LAMB" is repeated in the second line of the song, it uses an 1/8 note (300 milliseconds) so give "LEYAE" the TIME value of 200 and "M" the TIME value of 100. The pitch is 24.

Now we have

```
[MER<300,24>IY<300,22>HXAE<220,20>D<80,20>AX<
300,22>LIH<250,24>T<50,24>TUHL<300,24>LEYAE<
500,24>M<100,24>"]
```

MARY HAD A LITTLE LAMB

The song is written down exactly as it is to be entered into the Pathfinder

[MER<300,24>IY<300,22>HXAE<220,20>D<80,20>AX<300,22>

LIH<250,24>T<50,24>TUHL<300,24>LEYAE<500,24>M<100,24>
[LIH<250,22>T<50,22>TUHL<300,22>LEYAE<500,22>M<100,22>

[LIH<250,24>T<50,24>TUHL<300,27>LEYAE<500,27>M<100,27>

[MER<300,24>IY<300,22>HXAE<220,20>D<80,20>AX<300,22>

LIH<250,24>T<50,24>TUHL<300,24>LEYAE<500,24>M<100,24

IH<170,24>T<50,24>S<80,24>FLIY<200,22>S<100,22>WAX<200,2>Z<100,22>

WAY<220,24>T<80,24>AE<200,22>Z<100,22>SNOW<600,20>]

Helpful Hint:

The Pathfinder will power down if no activations are made after a pre-set number of minutes pass. Press any key or switch to continue. We also *highly recommend* that you have the song written down on a piece of paper exactly as you are going to program it before you begin to enter it. This will save you a lot of time.

⤷ Note: Use of the] and CTRLk in very long songs:

If you are programming a really long song, you may need to use the] and CTRLk to tell the Pathfinder to continue singing. For most songs, this will not be needed.

Select the] key on your keyboard and then select the CTRL key and the "k" key. (You will see a symbol for CTRLk () on your display, not the letters.) When you end a phrase or line of a song with a] (right square bracket) and the **CTRLk** characters, you *must* start the next phrase or line with a [(left square bracket).

If you do not do this, the Pathfinder will not know it is supposed to continue singing the song.

At the end of the entire song you *must* enter a] (right square bracket).

Storing the Programmed Song

If you do not have all the extra song programming characters on your Minspeak overlay (including TEXT and SPEECH MARKERS if you plan to use them), store the song from the Toolbox.

Storing Using TEXT and SPEECH MARKERS

When you store a song under an icon sequence, each time you activate the sequence to sing the song, the display will show the song just as you programmed it. Using TEXT and SPEECH MARKERS tells the Pathfinder to display only the title of the song. See below.

1. Activate STORE.
2. Select the icon sequence you want to use.
3. Activate INSERT TOOL.
4. Activate TEXT MARKER.
5. Type the title of the song.
6. Activate INSERT TOOL.
7. Activate SPEECH MARKER.
8. Activate the left square bracket [.
9. Type in the song exactly as you wrote it out for programming.
You can listen to your song at any time by activating SPEAK ALL.
You can use the arrow keys and DELETE keys to edit the song and correct any mistakes.
1. Activate INSERT TOOL.
12. Activate TEXT AND SPEECH MARKER.
13. Select OK.
14. Select OK again.

If you do not want to use TEXT and SPEECH MARKERS, skip steps #3-7 and #10-11.

**The following pages contain some musical charts
which you may find helpful when programming
songs into the Pathfinder.**



English Vowel Phonemes

Phoneme Symbol	Example	Phoneme Symbol	Example
Vowels			
aa	Bo b	ih	bi t
ae	ba t	ix	kiss es
ah	bu t	iy	bea t
ao	bou ght	ow	boa t
aw	bou t	oy	boy
ax	a bout	rr	bir d
ay	bi te	uh	bo ok
eh	be t	uw	lu te
ey	ba ke	yu	cu te

R-Colored Diphthongs

ar	b ar
er	b ear
ir	b eer
or	b or e
ur	p oor

Allophones

dx	rid er
lx	wi ll
q	we e at
rx	or ation
tx	Lat in

Silence _ (underscore)

Chart copied with permission from Digital Equipment Corporation

English Consonant Phonemes

Phoneme

Symbol	Example
--------	---------

Consonants

b	b et
ch	ch in
d	d ebt
dh	th is
el	bottle
en	button
f	f in
g	g uess
hx	h ead
jh	g in
k	K en
l	l et
m	m et
n	n et
nx	sing
p	p et
r	r ed
s	s it
sh	sh in
t	t est
th	th in
v	v est
w	w est
yx	y et
z	z oo
zh	az ure

Chart copied with permission from Digital Equipment Corporation

Pitch Notes and Pitch Values

PITCH NOTES AND PITCH VALUES

BASS

A bass clef staff with seven notes: C, C#, D, D#, E, F, and F#. Each note is a quarter note with a stem pointing down. The notes are positioned on the first line (C), first space (C#), second line (D), second space (D#), third line (E), third space (F), and fourth line (F#).

C	C#	D	D#	E	F	F#
1	2	3	4	5	6	7

A bass clef staff with seven notes: G, G#, A, A#, B, C3, and C#. Each note is a quarter note with a stem pointing down. The notes are positioned on the second space (G), second line (G#), third line (A), third space (A#), fourth line (B), fourth space (C3), and fifth line (C#).

G	G#	A	A#	B	C3	C#
8	9	10	11	12	13	14

A bass clef staff with seven notes: D, D#, E, F, F#, G, and G#. Each note is a quarter note with a stem pointing down. The notes are positioned on the first space (D), first line (D#), second line (E), second space (F), second space (F#), third line (G), and third space (G#).

D	D#	E	F	F#	G	G#
15	16	17	18	19	20	21

TREBLE

A treble clef staff with eight notes: G, G#, A, #A, B, C4, C#, and D. Each note is a quarter note with a stem pointing up. The notes are positioned on the first space (G), first space (G#), second line (A), second space (A#), third line (B), third space (C4), third space (C#), and fourth line (D).

G	G#	A	#A	B	C4	C#	D
20	21	22	32	24	25	26	27

PITCH NOTES AND VALUES (cont'd)

G G# A A# B C4 C# D
 20 21 22 32 24 25 26 27

D# E F #F G G# A A#
 28 29 30 31 32 33 34 35

B C5 C# D D# E F F#
 36 37 554.5 587.6 622.4 659.6 697.8 740.2

G G# A A# B C6 C#
 284.4 831.0 880.0 932.4 988.0 1046.6 1109.0

PITCH Value Code

PITCH NOTE	PITCH VALUE	OCTAVE RANGE
C4	25	C4 = Middle C
B	24	
A#	23	
A	22	ONE
G#	21	
G	20	
F#	19	OCTAVE
F	18	
E	17	
D#	16	BELOW
D	15	
C#	14	
C3	13	MIDDLE C
B	12	
A#	11	TWO
A	10	
G#	9	
G	8	OCTAVES
F#	7	
F	6	
E	5	BELOW
D#	4	
D	3	
C#	2	MIDDLE C
C2	1	

PITCH Value Code (cont'd)

PITCH NOTE	PITCH VALUE	OCTAVE RANGE
C4	25	C4 = Middle
C#	26	-
D	27	ONE
D#	28	
E	29	
F	30	OCTAVE
F#	31	
G	32	
G#	33	ABOVE
A	34	
A#	35	
B	36	MIDDLE C
C5	37	
C#	555.4*	TWO
D	587.6	
D#	622.4	
E	659.6	OCTAVES
F	697.8	
F#	740.2	
G	784.4	
G#	831.0	ABOVE
A	880.0	
A#	932.4	
B	988.0	MIDDLE C
C6	1046.8	

* For two octaves **above** middle C, you must enter PITCH values in Herz. This range is quite high for DECtalk speech voices. If you don't like the way a song sounds in this range, try entering the same notes with PITCH values that correspond to one octave **lower**.

TIME Values of Notes

NOTE TIME VALUE CHART

1/16 note 1/8 note 1/4 note 1/2 note dotted 1/2 note whole note

1/8 note triplet 1/4 note triplet | pair of 1/8 notes pair of 1/16 notes
 (= 1/4 note) (= 1/2 note)

(each note above gets part of the time value)

1/8 rest 1/4 rest half rest Whole

Appendix F: Helpful DECtalk Speech Tips

Alternate Pronunciations of Words

Sometimes a word has two different pronunciations. For instance, the word "record" can be pronounced REcOrd or reCORD.

On the next page there is a list of words with alternate pronunciations which the Pathfinder knows. You can use **CTRL<space>** to tell the Pathfinder which pronunciation to use.

For example, if you type "record" it is usually pronounced "REcOrd." If you want the *alternate* pronunciation for the word, insert the CTRL<space> tool just prior to typing in the word, **record**. Speech must be On.

If you frequently use the alternate pronunciation of a word that is not on the word list, you can store it into the Speech Dictionary, using CTRL<space>. Once it is in the Dictionary, the Pathfinder will always use the alternate pronunciation. If there comes a time when you do not want to use the alternate pronunciation, use CTRL<space> and then spell the word in your message.

We recommend that you **do not** store into the Speech Dictionary any alternate pronunciations of words that are on the list that follows, unless you always want the Pathfinder to use the alternate pronunciation.

The next few pages contain a list of words that have alternate pronunciations which DECtalk speech already knows. There is also a list of all the abbreviations that you can use which DECtalk speech knows how to pronounce.

The last part of this appendix explains how the Pathfinder and the DECtalk speech board handles numbers.

Words with Alternate Pronunciations

abstract	coordinate
abuse	decrease
addict	defect
advocate	delegate
affix	deliberate
ally	desert
alternate	desolate
animate	diffuse
annex	digest
appropriate	discharge
arithmetic	discount
articulate	dove
associate	duplicate
attribute	elaborate
august	estimate
bass	excerpt
baton	excuse
close	expatriate
combat	exploit
combine	export
compact	extract
complex	ferment
compound	frequent
compress	geminate
concert	graduate
conduct	impact
confederate	implant
confine	import
conflict	imprint
conglomerate	incense
console	incline
construct	increase
contest	insert
contract	insult
contrast	interchange
converse	intimate
convert	invalid
convict	just

lead	refresh
live	refund
minute	refuse
miscount	reject
misprint	relapse
misuse	relay
moderate	remake
object	rerun
overrun	research
perfect	resume
permit	retake
pervert	rewrite
polish	segment
postulate	separate
predicate	sow
predominate	subject
present	sublet
proceed	subordinate
produce	survey
progress	suspect
project	syndicate
protest	tear
read	torment
reading	transform
rebel	transplant
recall	transport
recap	upset
recess	use
record	wind
recount	wound
refill	

Chart copied with permission of Digital Equipment Corporation

Abbreviations Recognized by DECtalk Speech

Abbreviation	Word
all-in-one	All in one
Apr.	April
Assoc.	Associates
Aug.	August
Av.	Avenue
Ave.	Avenue
bldg.	Building
Blvd.	Boulevard
CH.	Chapter
cm.	centimeters
cms.	centimeters
Co.	Company
COD	"see-oh-dee"
cont.	continued
Corp.	corporation
Ctr.	center
CTRL	control
cu.	cubic
dec	DEC (as in Digital)
Dec.	December
Dist.	District
deg.	degrees
Dept.	Department
doz.	dozen
Dr.	Doctor (when followed by a name, i.e., Doctor West)
Dr.	Drive (when preceded by a name , i.e., West Drive)
e.g.	"ee""gee" (not 'for example')
ESC	escape
esp.	especially

Abbreviation	Word
est.	established ("eee ess tee")
etc.	et cetera
ext.	extension
Feb.	February
Flt.	flight
FOB	freight on board
fn.	Footnote
Fr.	Father
Fri.	Friday
Ft.	Fort (not 'Foot')
Gen.	General
gm.	grams
Gov.	Governor
hrs.	hours
i.e.	"eye ee" (not 'that is')
Inc.	Incorporated
Intl.	international
Jan.	January
Jr.	Junior
Jul.	July
Jun.	June
kg.	kilograms
kgs.	kilograms
km.	kilometers
lb.	pounds
lbs.	pounds
Ltd.	Limited
Mar.	March
mg.	milligrams
mgs.	milligrams
misc.	miscellaneous
ml.	milliliters
Mon.	Monday
Mr.	mister
Mrs.	missus

Abbreviation	Word
msde.	merchandise
msec.	milliseconds
msecs.	milliseconds
mss.	manuscripts
Mt.	Mount
Nov.	November
Oct.	October
oz.	ounces
ozs.	ounces
Pl.	Place
pp.	pages
ppd.	post paid
Pres.	President
Rd.	Road
recd.	received
Rep.	Representative
Rev.	Reverend
rsts	"ris-tis"
Rte.	route
Sen.	Senator
Sep.	September
Sept.	September
sq.	square
Sr.	Senior
St.	Saint (when followed by a name)
St.	Street (when preceded by a name)
Sun.	Sunday
Thu.	Thursday
Thurs.	Thursday
Tue.	Tuesday
Tues.	Tuesday
Univ.	University
USA	"you-ess-aye"

Abbreviation	Word
---------------------	-------------

Vol.	Volume
vs.	versus
yds.	yards

Exceptions*

nt.wt.	net weight
op.cit.	op cit
p.p.d.	post paid
pat.pend.	patent pending.
Wednesday	

- These abbreviations will **not** work at this time (12/90)

List copied with permission of Digital Equipment Corporation

Producing a Long "A" in the Middle of a Sentence

Usually an "A" by itself in the middle of a sentence is pronounced as "uh," but there may be times when you want it to be pronounced as a long "A" ("ayee"). To achieve this, type in the capital letter "A" and then add a punctuation mark immediately following it:

A, A? A! (A A;

This will change the inflection of the sentence somewhat, but it will produce a long "A".

To Create a Pause within a Word or between Words

Enclose an underscore, `_`, and a number designating the length of the pause within square brackets. For example, `[_<1000>]` tells the Pathfinder to pause for 1 second; `[_<2000>]` creates a 2 second pause. Changing the number changes the amount of seconds the Pathfinder pauses.

Creating a Pause after a Period or a Comma

Generally, the Pathfinder pauses for ½ second after a period, before beginning the next sentence. It will pause for 1/6 of a second after a comma.

If you want to create a longer pause after a period or a comma, add the following indicator to your message:

[:pp4500]

This adds a 4.5 second pause to the regular period pause. You can type in the number you want to use.

[:cp4800]

This creates a 4.8 second pause after the regular comma pause. You can type in the number you want to use.

If you want to return the pauses to their default settings, enter these indicators:

For a period default:

[:pp0]

For a comma default:

[:cp0]

How DECTalk Speech Handles Numbers

1998 is a year. 1,998 is a big number. 1998 (no comma) is spoken as "nineteen-ninety-eight." 1,998 (with a comma) is spoken as "one-thousand-nine-hundred-and-ninety-eight." *Any four-digit number is spoken as a year unless you add a comma after the first digit.* Commas are not necessary if a number contains more than four digits (e.g., 10,375). Numbers containing up to 9 digits are pronounced correctly.

Examples of Entering Math Problems:

When you enter math problems, they must be spaced correctly for DECTalk speech and the Pathfinder to pronounce them correctly.

For example: 5 - 3 = 2

Enter this problem in the following manner:

5 (space) -3(no space between "-" and "3" but add a space after "3")
=(space) 2.

The problem looks like this when it is entered:

5_ -3_=_2 ("_" indicates a space)

Examples of Spacing for Other Math Problems:

2*8_=_16; 1*2_=_2 ("*" indicates "times" but the Pathfinder speaks "asterisk".)

2_+8_=_10

1 ½ must be entered as 1&_1/2

(1/2 can be entered by activating "1" "slash" "2".)

At this time DECTalk speech and the Pathfinder cannot speak division problems correctly.

Examples of Entering Dates:

DECTalk speech does not recognize dates if you enter them the way they are normally written, e.g.,

6-8-1990 or 11/20/95

Enter dates in the following manner for them to be spoken correctly:

6-Aug-1990 or 20-Nov-1995

These dates are spoken as "August 6, 1990" and "November 20, 1995".

Example of Entering a Phone Number:

Phone numbers must be entered by placing spaces and commas as shown below:

3 3 0, 2 6 2, 1 9 8 4 or 2 6 2, 1 9 3 3

Appendix G: About Dependent Icons

☞ If you are using one of the Unity overlays, dependent icons are already programmed into the vocabulary. If you are creating your own overlay from scratch, you will want to create dependent icons to make programming the overlay and the retrieval of words from the overlay easier.

Let's say you have an icon that represents NOUN and an icon that represents NOUN PL (plural nouns). Let's also say that you have stored your most commonly used nouns under icon sequences that end with the icon: NOUN. For example, HOUSE+NOUN = home. If you want the word "homes" on your overlay, traditionally you would have to store it under another sequence that ended with the icon: NOUN PL. For example, HOUSE+NOUN PL = homes. This means you would have to store under icon sequences, all the nouns you want for your vocabulary. Then you would have to store all those same nouns again as plurals!

With Pathfinder you can make NOUN PL a **dependent** icon. This means that an icon sequence ending with NOUN PL will be dependent on an identical sequence which ends in NOUN. For example, if you store NOUN PL as a dependent icon and then activate the sequence, HOUSE+NOUN PL, the Pathfinder will look to see what is stored under HOUSE+NOUN. It will see "home" and automatically add an "s" to it. You would see "homes" on your display. Dependent icons save you hours of storing. You can do the same thing with verbs by making icons such as VERB+ed and VERB+ing dependent on sequences that contain the icon: VERB. You store all of the regular verbs you want to use and the Pathfinder will put the correct ending on them.

The next page explains how to create dependent icons.

Creating Dependent Icons

➡ this example assumes that you have added a space after any words that you have stored in your vocabulary.

In this example we are going to assign the icons NOUN and NOUN PL to separate keys. We are also going to make NOUN PL dependent on the NOUN icon. Then we will type in instructions for the Pathfinder to follow whenever the NOUN PL key is activated. These instructions will tell the Pathfinder to find the word stored under a sequence that ends in NOUN, to delete the space after the word that is stored there and to add an "s" and a space to the word.

Assign NOUN to a Key:

1. Go to the Toolbox and open the ASSIGN KEY MENU.
2. Select the key you want to assign NOUN to.
3. Activate the *Minspeak Assignment* option (#1).
4. Make sure the *Assign Key As:* box is set to *ICON*.
5. Activate *Rename Icon*.
6. Type in the name: NOUN.
7. Activate OK twice.

Assign NOUN PL to a Key:

8. Activate the *Select Next Key to Define* option (#6).
9. Select the key you want to assign NOUN PL to.
10. Activate the *Minspeak Assignment* option (#1).
11. Make sure the *Assign Key As:* box is set to *ICON*.
12. Activate *Rename Icon*.
13. Type in the name: NOUN PL.

Make NOUN PL Dependent on NOUN and Type in Directions:

14. Activate OK.
15. Activate the *Icon Type* option (#5) twice.
16. You see a box that asks if you want to make the icon NOUN PL a dependent icon. Answer YES.

17. Select the icon you want NOUN PL to depend on. In this case select the key where you stored NOUN.
18. Now you must tell the Pathfinder what you want it to do when it sees NOUN PL.
 1. Activate INSERT TOOL.
 2. Type "d" and select <DELETE CHARACTER> or activate *More Items* until you see <DELETE CHARACTER> and then select it.
Type the letter "s" and add a space.
19. Activate OK until you exit the ASSIGN KEY menu.

NOUN PL is now a dependent icon. It will add an "s" to any word stored under a sequence that ends with the NOUN icon.

Some Notes and Exceptions for Dependent Icons

You can create dependent icons for many different circumstances. You could have a dependent icon that automatically adds "the" before any noun (e.g., the HOME) or one that adds "es" to nouns (e.g., HOUSES).

However, some irregular nouns and verbs will have to be programmed separately. For example, suppose you have the noun WOMAN stored under the sequence WOMAN NOUN. Since you have made NOUN PL a dependent icon, if you activate WOMAN and NOUN PL, the Pathfinder will display "Womans" instead of Women. You would have a hard time writing a macro to make a dependent icon work with all irregular nouns or verbs because they do not follow a standard format. Because of this, you will have to store irregular nouns and verbs under separate sequences.

For example, using the regular storage process, you can store the noun "women" under WOMAN+NOUN PL. In this case NOUN PL is not used as a dependent icon. It is simply part of an icon sequence for the word, "women".

To Return a Dependent Icon to a Normal Icon

If, for some reason, you decide that you no longer want an icon to be dependent, you can return it to a normal icon.

1. Go to the Toolbox and open the ASSIGN KEY menu.
2. Select the key which you want to return to normal.
3. Activate the *Minspeak Assignment* option (#1)
4. Make sure the *Assign Key As* option is set to **ICON**.
5. Activate the *Icon Type* option (#5) until you see **NORMAL**.
6. Activate OK twice to exit the menu.

Your dependent icon is no longer dependent, and will now act as a regular icon.

Appendix H: External Speaker and Earphone

You can connect an external speaker and earphone to the Pathfinder. (External earphones should be used if you are also using Auditory Prompts.) Their connectors must fit into the following jack sizes on the Pathfinder's top panel:

External Speaker Jack—3.5 millimeter mono; use an 8 ohm speaker

Earphone Jack—2.5 millimeter; mono



Warning!

Do not plug an earphone into the **EXTERNAL SPEAKER** jack.
This could cause damage to your hearing.

Appendix I: Creating Macros That Contain Menu Options

Many menus have options that are not listed in the All Tools List as insertable tools. For example, you cannot insert the tool "1-Switch" into a macro. But you can create a macro that includes part of the SELECTION TECHNIQUE MENU that allows you to change your selection technique to 1-Switch Scanning.

Let's say that you normally use the keyboard for your selection technique, but when you get tired you like to change to 1-switch scanning. You might like to create a macro that allows you to change from keyboard to 1-switch scanning. You can store this macro on your overlay and activate it whenever you want to switch selection techniques.

In order to create a macro that includes part of a menu you have to do some homework first.

- You need to know and write down (or remember) the keys that you have to activate in the menu.
- For our 'Keyboard to Scanning' example, you should set up the Keyboard and 1-Switch selection techniques **before** you begin to create the macro so that all the options are set the way you want them when you use a switch.

Homework:

Let's find out what key activations you have to make in the SELECTION TECHNIQUE MENU to change from Keyboard to 1-Switch Scanning.

1. Go to the Toolbox and open the SELECTION TECHNIQUE MENU.
2. The current selection technique should be "Keyboard."
1. Activate option #1.
2. Activate OK.

To Return to the Keyboard

1. Activate TOOLBOX.
2. Activate SELECTION TECHNIQUE MENU
3. Activate Option #1 *four* times.
4. Activate OK.

To Create the "Keyboard to Scanning" Macro:

1. In the Toolbox, activate STORE CORE.
2. Select the icon sequence you want to store the macro under, for example, you might use the pointing finger (REMEMBER) and the medical icon (MEDICAL).
3. Activate OK.
4. Activate "Spell Text for Message" (we will add a message to be spoken at the end of the macro.)
5. Activate INSERT TOOL. Activate "S" on the keyboard.
6. Activate the SELECTION TECHNIQUE menu
7. Type #1.
8. Activate INSERT TOOL
9. Activate OK.
12. Now type a message if you want one, perhaps: "Please plug in my switch."
13. Activate OK when you have completed your message. Activate OK again.
14. Close the toolbox. Go to Minspeak mode.
15. Activate the icon sequence you chose (REMEMBER MEDICAL).
16. The Pathfinder should speak, "Please plug in my switch." Activate your switch once. If you are using Row/Column scanning, the top row of the display screen should be highlighted in red, ready for you to scan.

 If you want to **return** to using the keyboard (that is, if you want to toggle between two selection techniques) you must create and store a second macro that takes you back to the first technique.

➡ Besides storing these macros under icon sequences, you can store them in Activity keys or assign them to Core keys.

Other Macros That Might Be Used:

You could modify other selection technique settings for different times of day without going to the Toolbox.

( **Remember**, you can also do this by saving your settings and storing them on your overlay, see pages 283.)

You could turn scan beeps on and off using menu selection #3 from the FEEDBACK menu.

Appendix J: Computer Commands

for IBM Keyboard and Keypad Keys

➡ Create the <esc> command that you see below by activating **CTRL** [(Control left bracket) in your Toolbox. This sequence (CTRL[]) is shown as <esc> in all commands below. On your Pathfinder display you will see a small left arrow (←).

Keyboard Keys

	<esc>f10.
	<esc>f11.
<esc>backspace.	<esc>f12.
<esc>capslock.	<esc>home.
<esc>delete.	<esc>insert.
<esc>down.	<esc>left.
<esc>end.	<esc>pagedown.
<esc>enter.	<esc>pageup.
<esc>esc.	<esc>pause.
<esc>f1.	<esc>print.
<esc>f2.	<esc>right.
<esc>f3.	<esc>scroll.
<esc>f4.	<esc>tab.
<esc>f5.	<esc>up.
<esc>f6.	
<esc>f7.	<esc>,hold,alt.
<esc>f8.	<esc>,hold,control.
<esc>f9.	<esc>,hold,shift.

Specify a right or left key by placing an "l" or "r" in front of the key name.

For example, to hold the *left shift key*: <esc>,hold,lshift.

Keypad Keys

<esc>kp0.
<esc>kp1.
<esc>kp2.
<esc>kp3.
<esc>kp4.
<esc>kp5.
<esc>kp7.
<esc>kp8.
<esc>kp9.
<esc>kpslash.
<esc>kpstar.
<esc>kpup.
<esc>numlock.
<esc>sysreq.

<esc>kpdel.
<esc>kpdown.
<esc>kpend.
<esc>kpenter.
<esc>kphome.
<esc>kpins.
<esc>kpleft.
<esc>kpminus.
<esc>kppgdn.
<esc>kpplus.
<esc>kpright.

Computer Commands (cont'd)

For Apple/Macintosh Computer Keyboard and Keypad

➡ Create the <esc> command that you see below by activating **CTRL[** (Control left bracket) in your Toolbox. This sequence (CTRL[) is shown as <esc> in all commands below. On your Pathfinder display you will see a small left arrow (←).

Keyboard Keys

<esc>backspace.	<esc>f9.
<esc>capslock.	<esc>help.
<esc>delete.	<esc>home.
<esc>end.	<esc>insert.
<esc>esc.	<esc>left.
<esc>f1.	<esc>pagedown.
<esc>f10.	<esc>pageup.
<esc>f11.	<esc>pause.
<esc>f12.	<esc>print.
<esc>f13.	<esc>return.
<esc>f14.	<esc>right.
<esc>f15.	<esc>scroll.
<esc>f2.	<esc>tab.
<esc>f3.	<esc>up.
<esc>f4.	<esc>,hold,command.
<esc>f5.	<esc>,hold,control.
<esc>f6.	<esc>,hold,option.
<esc>f7.	<esc>,hold,shift.
<esc>f8.	<esc>hold,alt.

Specify a right or left key by placing an "l" or "r" in front of the key name.

For example, to hold the *left shift key*: <esc>,hold,lshift.

Keypad Keys

<esc>kp0.
<esc>kp1.
<esc>kp2.
<esc>kp3.
<esc>kp4.
<esc>kp5.
<esc>kp6.
<esc>kp7.
<esc>kp8.
<esc>kp9.
<esc>clear.
<esc>kpdel.
<esc>kpdown.
<esc>kpend.
<esc>kpenter.
<esc>kpequal.
<esc>kphome.
<esc>kpins.
<esc>kpleft.
<esc>kpminus.
<esc>kppd.
<esc>kppgdn.
<esc>kppgup.
<esc>kpplus.
<esc>kpright.
<esc>kpslash.
<esc>kpstar.
<esc>kpup.
<esc>numlock.
<esc>reset.

Appendix K: Save Your Setup Settings

Summary:

The insertable tool, SETTINGS MENU, allows you to save from 1-3 different groups of settings.

You may want to set up more than one selection technique and its options so that you can use different techniques during different times of the day. To do this, you need to assign the SETTINGS MENU key to your overlay or to store it under an icon sequence. The Pathfinder will save anything you have set in the SELECTION TECHNIQUE menu, AUDITORY PROMPT menu, FEEDBACK menu, SPEECH menu and OUTPUT menu.

To store the SAVE SETTINGS MENU under an icon sequence:

1. Activate STORE CORE.
2. Select the sequence you want to use.
3. Select *Synthesized* and *Spell Text for Message*.
4. Activate INSERT TOOL.
5. Activate "s" and *More Items* until you see SETTINGS MENU. Activate it.
6. Activate OK.

To Save your settings into the Menu

1. Make sure you are currently using the first set of settings you want to save.
2. Go to Minspeak and activate the sequence you saved the SETTINGS MENU under.
3. Save your current settings to #1 by activating *Save #1 Settings*.
4. Now set up your second technique and save those settings to *Save #2 Settings*.
5. Activate OK to until you exit the menu.

When you are using one group of settings and want to change to another, activate the SETTINGS MENU and then activate the "Use #2 Settings option, or whichever option you have saved settings to.

Once you have saved the settings you want in the SETTINGS MENU, you can store the tools, SETTINGS 1, SETTINGS 2 and SETTINGS 3 under icon sequences on your overlay. Then you can delete the SETTINGS menu if you want. To go to different settings, you just activate the sequence you stored the settings under.

➡ You can assign all these tools to keys on your overlay if you have any keys available.

Store the SETTINGS 1, SETTINGS 2 and SETTINGS 3 tools:

Repeat the Menu storing steps on the previous page to store the SETTINGS 1, SETTINGS 2 and SETTINGS 3 tools under sequences. Replace SETTINGS MENU in step #5 with SETTINGS 1,2 or 3.

To Change Your Settings

1. Activate the sequence where you stored the SETTINGS 1 tool.
2. You see a box that says: "Are you sure you want to replace your current settings with Settings #1?" Answer YES to store the settings.

Now when you are tired of using your headpointer and want to start using a switch or the keyboard, you can activate the sequence where you stored your switch or keyboard settings and the Pathfinder will automatically switch you to the new settings. You can delete the SETTINGS MENU if you are not going to set new settings any time soon.

Appendix L: Troubleshooting



Warning!

DO NOT open the Pathfinder case!

If you encounter a problem with your unit,
call the PRC Service Department at 1-800-262-1990.
If you are outside the U.S., call your local distributor.

Do not attempt to open the case!

Battery Charger Warning:



Warning!

Use *only* the battery charger that came with your Pathfinder. Other chargers may damage the batteries or the Pathfinder. This also applies to your headpointer battery charger.

The Battery Status display window on the Pathfinder screen will tell you the battery status at any given time. You must turn the display window **ON** in the FEEDBACK menu, in order to see the status of the batteries.



Remember! You cannot STORE when the batteries are low.

➡ **If the battery charger is plugged in and the battery LED is solid red:**

This indicates that the battery pack is not communicating with the Pathfinder. Allow one hour of charging. If the battery charger LED continues to remain red after the hour of charging, call the PRC Service Department.

➡ **If the batteries in your Pathfinder or Headpointer need to be replaced, dispose of the old batteries properly.** Follow the regulations of your community, state or country for the disposal of batteries.

➡ **Important Note!**

When you turn on the Pathfinder On you must wait a few minutes for the software to open up. The computer screen will remain dark while the software loads. **Do not press the On/Off button again.** Doing so will lock up your device.

If you have locked up the Pathfinder:

Try pressing and holding the I/O and TOOLS keys at the same time.

If the above suggestions does not help:

Press all four keys on the top left side of the Pathfinder case **at the same time**. This re-boots the Pathfinder. All unsaved work will be lost.

➡ **The Pathfinder does not seem to be speaking:**

Go to the toolbox and make sure the speech is On.

Go to the **FEEDBACK MENU** in the Toolbox and set the "Volume" option to it highest setting. Check to see if *Speak when Spelling* is On or Off.

If it is On, check to see if it is set to Words or Letters. If it is set to Words, the Pathfinder will not speak until a whole word is spelled.

➤ **If you are having trouble storing or retrieving messages:**

You may need to exit a theme before you can use the Pathfinder the way you want to. Check the status bar at the bottom of the Pathfinder's display area. If you are in a theme, the theme icon(s) will be displayed followed by a colon.

➤ **If you are using a T-TAM and some characters are not appearing on your computer screen:**

If you send chunks of text over 60 characters in length through the T-TAM, it will probably drop some of those characters. Turn Character Pacing on in the OUTPUT menu or send smaller chunks of text.

➤ **You get double-hits from keys or a 'bouncing' effect when you activate keys, or when you move across the keyboard to activate a key:**

The *Acceptance* and/or *Release* times are too short. Increase these times.

Open the SELECTION TECHNIQUE menu and activate your selection technique option. If the Pathfinder keys are bouncing as you move across the keyboard to your target key, **increase** the *Acceptance Time* and then try activating the keys. If you still experience bounce, increase the *Acceptance Time* a little more.

If you are getting double-activations, **increase** the *Release Time* and practice with the keys. If double-activations continue, increase *the Release Time* a little more.

➡ Your scanning speed is too slow


Go to the SELECTION TECHNIQUE menu, select your technique, then set the *Scanning Speed* option to a faster number.

To keep from scanning keys that you don't want to scan, read the pages on **Predictive Selection, Modify Predicted Keys** and **Linear Scan**

➡ About accessing masked messages

If you transfer vocabulary that contains masked messages from a DeltaTalker or Liberator into the Pathfinder, use the *OK* key on the Pathfinder to access the masked messages.

➡ You did something and an Error Box (below) appeared on your display.

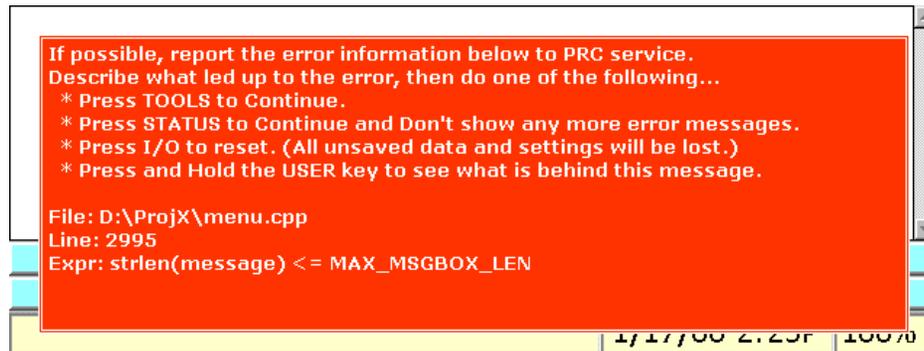

If you see this error box on your display, PRC would greatly appreciate it if you would take a few minutes to write down the error code at the bottom of the box (i.e., File; Line; Expr;) and to try to re-construct what you were doing just before this message appeared. Please phone the

PRC Service Department and give one of the service representatives the information, or mail the information to PRC. The information can be very helpful to our software engineers as they work to debug the software.

If you need to see the screen you were working on before the box appeared, press and hold the USER key on the front of the Pathfinder. The error box will disappear and you can see the display. When you release the USER key the error box will reappear.

We recommend that you do not choose to turn the error message off (the second option in the box). The information it contains can be very helpful if you have to call service or send the Pathfinder in to PRC.

About High-Efficiency Fluorescent Lights and the Use of Infrared:

There is a compatibility problem between many infrared controls and high-efficiency fluorescent lights.

If you are using the Pathfinder in a room that has high-efficiency fluorescent lighting, the infrared, including headpointing, computer access, remote control commands, etc., will not work. Turn the lights off and the IR will work fine.

A Note About Shipping Mode

When you put the Pathfinder in shipping mode (from the MAINTENANCE menu), wait 10 seconds after activating the option, for the shipping mode to take effect. Once the device is in shipping mode, you cannot activate keys.

Cleaning the Pathfinder

Regular Cleaning

Clean the case with a soft, damp cloth.

Clean the display screen with the same kind of product you use to clean glass. Use a soft cloth to wipe off the glass cleaner.

Disinfecting

To disinfect the equipment, mix $\frac{1}{4}$ cup vinegar with 1 cup water. Wipe all equipment with a cloth *dampened* in this mixture. Allow the equipment to air dry.

Index

Index

- 1-switch 67 85-100
- 2-switch 69, 85-100
- 32-location overlay 51, 137
- 64-location overlay 51, 137
- 1-switch scan mode 1 85
- 1-switch scan mode 2 86
- 2-switch scan mode 1 87
- 2-switch scan mode 2 88

- Abbreviations, speech list 265
- Acceptance time 91
- Activation beep 195
- Activation delay 93
- Activation feedback 195
- Activities
 - about 35, 151
 - create 151-153
 - copy 163
 - delete 155
 - delete dynamic 159
 - dynamic 35, 157-159
 - dynamic menu 157
 - linking 161
 - menu 151-153
 - modify 160
 - review dynamic 159
 - set dynamic 157
- Activity row 35
- Add tools 140, 173-182
 - to keys 140, 176, 178-179
 - to sequences 140, 182
- Add words to:
 - speech dictionary, *see Power Tools*
 - word prediction dictionary 118
- All Tools List 209-224
- Alternate pronunciations 263, *see also Power Tools*
- Appendices 209-283
- Area
 - about 47
 - blank area 49
 - copy an 53
 - erase an 54
 - overlay menu 51
 - rename an 52
 - select 49
 - switch areas 49
 - user area menu 49
- Arrow keys *see Power Tools*
- Assign different icon to a key 175
- Assign icon 177
- Assign key 173-180
- Assign key menu 173
- Assign tool
 - to a key 176-179
 - to a sequence 182
 - tools and text to a key 176-179, 184
 - tools and text to a sequence 183-185
- Assign-macros 184-185
- Audio feedback 195
- Auditory prompt menu 189
- Auditory prompts 189-193
 - about 189
 - category 191
 - delete 192
 - mid-sequence 193
 - on/off 189
 - storing 189-193
- Auto-line feed *see Power Tools*
- Auto repeat rate 94
- Auto spell 198
- Average pitch *see Power Tools*
- Averaging headpointer 77

- Backlight shut-off 25, 203
- Batteries
 - about 23

- display status 10, 23, 195
- disposal 285
- Battery charger
 - headpointing 82
 - Pathfinder 23, 195, 285
- Battery charger connector 21
- Battery door 21
- Battery indicator icon 12
- Battery status area 23, 195
- Baud rate *see Power Tool*
- Beep
 - Activation 195
 - Scan 195
- Blank notebook *see Power Tools*
- Blank overlays 51
- Block marker *see Power Tools*

- Calibrate display screen 203
- Calculator *see Power Tools*
- Caps/Lock 13
- Category prompts 189, 191
- Change a message 131, 134
- Change voices in a message *see Power Tools*
- Character, delete 180-181, 192, 213
- Character pacing *see Power Tools*
- Character/Word *see Power Tools*
- Charge (battery) 23, 81
- Charger (battery) 23, 81
- Charger light (LED) 12, 23, 81
- Choose activity 12, 35
- Cleaning the Pathfinder 290
- Clear display 180-181, 212, 219
- Clock
 - display 195, 204
 - set date 204
 - set time 204
 - set time/date format 204
- Close notebook *see Power Tools*
- Colors
 - Notebook status bar *see Power Tools*
- IR Headpointer pointer 81
 - set default 201
 - system 200
 - word prediction 122
- Column/row scan 57, 87-89
- Computer 13, 47, 55, 77, 128, *see also Power Tools*
 - 'backspace' command (word prediction) 121
 - overlay 55
 - icon, front of case 13
- Computer commands *see Power Tools*
- Computer emulation *see Power Tools*
- Connect to other devices *see Power Tools*
- Consonants 254
- Contrast, display 195
- Control (CTRL) commands *see Power Tools*
- Coordinates, key 37
- Copy activity 47, 163
- Copy calculator answer *see Power Tools*
- Copy text *see Power Tools*
- Copy user area 53
- Core vocabulary 105, 127, 129, 151
- Create activity 151-153
- CTRL *see Power Tools*
- CTRL serial marker *see Power Tools*
- Cut text *see Power Tools*

- Data logging, 167, 224-231
- DECTalk speech 259-268, *see also Power Tools*
- DECTalk voice *see Power Tools*
- Defaults
 - Reset 205
 - system colors 198, 200-201
- Delete character 180-181, 192, 213
- Delete last selection 213

Delete notebook *see Power Tools*
 Delete word 213
 Delete words from dictionary
 Speech dictionary *see Power Tools*
 Word prediction dictionary 118
 Dependent icon 271-274
 Dictionary
 speech *see Power Tools*
 word prediction 117-119
 Digitized speech
 erase 132
 store 131
 Direct selection
 Keyboard 61
 Headpointer 63-65
 Disinfecting Pathfinder 290
 Display battery status 195
 Display
 brightness/contrast 195
 calibration 204
 clear 180-181, 212, 219
 regular 11
 text size, 195
 word prediction 107-111
 Display clock 195
 Display contrast 195
 Display date 195
 Disposal, batteries 285
 Dual (2) switch scanning 69, 87-88
 Dynamic activity 157-159

 Earphone, external 17, 275
 Earphone jack 17, 275
 Earphone volume 43, 195
 Edit core 135
 Edit menu *see Power Tools*
 Edit mode *see Power Tools*
 Edit text *see Power Tools*
 Editing tools *see Power Tools*
 Emulation, computer *see Power Tools*
 English consonant phonemes 254
 English vowel phonemes 253
 Erase
 all vocabulary 203
 areas 54
 auditory prompts 192
 key assignment 181-182
 macros 185
 messages 132, 135
 notebooks *see Power Tools*
 recorded speech 132
 tools 181
 user area 54
 words from dictionary
 speech *see Power tools*
 word prediction 118
 workspace *see Power Tools*
 Error reporting 204
 Exit theme 145
 External earphone 17, 275
 External speaker 17, 275

 Feedback beep 195
 Feedback menu 195
 Flash predictor keys 115
 Formatted notebook *see Power Tools*
 Frequency/recency 113
 Front of Pathfinder 11-13

 Go to computer 55, 128, 214
 Go to Minspeak 128, 198, 212
 Go to spell 128, 198, 215
 Graphics, importing (*see MTI Manual*)

 Handle, Pathfinder's 18-19
 Hardware tests 204
 Headpointer 13, 63-65, 77-84
 Headpointing 13, 57, 63-65, 77-84
 averaging 77
 batteries 82
 horizontal/vertical adjustment 79
 mounting 63, *see also Mounting Instructions*

- pointer color 81
- predictive selection 95, 97, 198
- reset options to default 205
- sensor cord position 63-65
- Hide/show activity 161
- Hide/show key 187
- Hide/show notebook *see Power Tools*
- Highlight text *see Power Tools*

- IBM memory transfer *see MTI Manual*
- IBM computer commands 279
- I/O button 11, 25
- Icon
 - About 27-29, 33
 - assign different 177
 - delete 213
 - dependent 271-273
 - examples 26-33
 - import *see MTI Manual*
 - intensity 29
 - label 151-157
 - masking 288
 - Minspeak 27-33
 - Normal (ordinary) 274
 - one-shot 147
 - ordinary (normal) 274
 - pick 1 icon 139, 140-141
 - pick 2 Icon 139, 140-141
 - prediction 29, 95-97, 195
 - rename 176-177, 272-274
 - selecting 26-33
 - sequences 26-33
 - set 1 icon 139, 143
 - set 2 icon 137, 143
 - single meaning 147-149
 - tutor 199
- Infrared *see Power Tools*
- Insert tool 140, 173-182, 215
- IR *see Power Tools*
- IrDA *see Power Tools*

- IR window location 17, 21
- IRR *see Power Tools*

- Joystick 17, 57, 71, 75, 89-100

- Key assignments 173-182
- Key coordinates 37
- Key locations 37
- Key prompts 189-191
- Keyboard 15, 57, 61
- Keypad 15

- LAM instructions 165, 224-230
- Language acquisition monitor (LAM) 165, 224-230
- Lights (LEDs)
 - activation 195
 - battery 12
 - computer 13
 - intensity 196
 - shift/caps lock 13
 - spell/Minspeak 15
 - system status 12
 - tools 11
 - user area 13
 - user programmable key 12
- Linear scan cut-off 100
- Linking activities 161
- List notebooks *see Power Tools*
- Load graphics *see MTI Manual*
- Load memory *see MTI Manual*
- Locations, key 37
- Lock notebook *see Power Tools*
- Lock system 203
- Long "A" (create in sentence) 268

- Macintosh computer commands 281-282
- Macintosh memory transfer *see MTI Manual*
- Macros 183-185, 276-278

Maintenance menu 203
 Main voice *see Power Tools*
 Masking messages 288
 Memory card door 17
 Memory review 165
 Memory transfer interface *see MTI Manual*
 Menu choices 39-42
 Menu prompts 189-193
 Menu scanning (on-screen) 43, 75, 109
 Menu voice *see Power Tools*
 Menus
 activity 151-163
 assign key 173
 auditory prompt 189
 computer access (in User Area menu) 55
 edit 214, *see also Power tools*
 feedback 195-196
 IR setup *see Power Tools*
 maintenance 203-204
 notebook *see Power Tools*
 output *see Power Tools*
 review vocabulary 165, 232-234
 scanning (on-screen) 75, 109
 selection technique 59
 speech *see Power Tools*
 transfer memory *see MTI Manual*
 user area 47
 word prediction 111
 Merge/Load Module *See MTI Manual*
 Microphone 13
 Minserts *see Power Tools*
 Minspeak
 about 27-33
 icon, front of case 13
 and Word Prediction 103-105, 112
 Mispronounced words 263-264
 Modify activity 160
 Modify predicted keys 97
 Mounting headpointer *see Headpointer Mounting Instructions*
 Mounting plate 19
 MTI *See MTI Manual*
 Music 236-259
 Notebook menu *see Power Tools*
 Notebook, *see Power Tools*
 Number of rescans 75
 Numbers (DECtalk speech)
 dates 269-270
 math problems 270
 phone numbers 270
 One-shot icon 147
 On/Off
 character pacing *see Power Tools, see MTI Manual*
 Pathfinder 11, 25
 prompts 189
 speech 221
 toolbox 12
 word prediction 107
 on-screen scanning 75, 109
 Open/close notebook *see Power Tools*
 Ordinary icon 274
 Overlay
 128-location 51
 64-location 51
 32-location 51
 keyguards 15
 computer 55, *see also Power Tools*
 default 51
 protector 15
 toolbox 15
 Unity 14, 51
 Output menu *see Power Tools*
 Output method *see Power Tools*
 Output mode *see Power tools*
 Pathfinder 11-25, 290

Paste text *see Power Tools*
 Pause, create in speech 136, 268
 Phonemes, consonants 254
 Phonemes, vowel 253
 Pick 1 icon 139, 141
 Pick 2 icon 139, 141
 Pictures, importing *see MTI Manual*
 Pitch, average *see Power Tools*
 Pitch notes 255
 Pitch range *see Power Tools*
 Pitch value code 255-256
 Pitch values 257
 Pointer color 81
 Predictive selection 95
 Predictor keys (Word Prediction) 115
 Pre-stored songs 236, 237
 Print notebook *see Power Tools*
 Prompts, auditory 190-193
 Prompts, menu 189
 Prompt voice *see Power Tools*
 Pronunciation, alternate 263-264
 Pronunciation, speech dictionary *see Power Tools*

 Quarter row/column scan 51, 73

 Recency/frequency in Word Prediction 113
 Record
 message 129, 131
 erase 132
 store 131
 Regular icon 271
 Regular workspace *see Power Tools*
 Release time 95
 Rename icon 176-177, 272
 Rename notebook *see Power Tools*
 Rename user area 52
 Reset defaults in menus 203
 Resize text area 195
 Retrieve Minserts *see Power Tools*
 Return *see Power Tools*

 Review
 dynamic activities 159
 speech dictionary *see Power Tools*
 to display 165
 to display and speaker 165
 to output 165, 232
 word prediction dictionary 117-119
 vocabulary menu 165
 Row/column scan 57, 87-89
 Save notebook *see Power Tools*
 Save setup settings 283
 Scan
 1-switch 67, 85-86
 2-switch 69, 87-88
 32-64-128 location 51, 137
 acceptance time 91
 activation delay 93
 auto scan restart 75
 beep 195
 centering 99
 column/row 57, 87-89
 direction 89
 joystick 57, 71, 75, 99
 LED intensity 195
 linear scan cutoff 100
 menu scanning (on-screen) 75, 109
 mode
 1-1-switch 85
 1-2-switch 86
 2-1-switch 87
 2-2-switch 88
 on-screen 75, 109
 quarter/row column 57, 73
 rate 67, 85, 90
 release time 92
 row/column 57, 87-89
 selection type 75
 speed 90
 start position 99
 type 75
 menus 75, 109
 Screen as switch 75

Screen calibration 203
 Search and replace *see Power Tools*
 Select user area 49
 Select words in word prediction 109
 Selection technique menu 59
 Selection. techniques
 1-switch 57, 67, 85-86
 2-switch 69, 87-88
 headpointer 13, 57, 63-65,77-84
 joystick 17, 57, 71, 75, 99
 keyboard 57, 61
 quarter row/column 51, 73
 row/column 57, 87-89
 Sensor cord, headpointing 63-65
 Sequences, icon 27-31
 Serial marker *see Power Tools*
 Serial mouse *see Power Tools*
 Serial off *see Power Tools*
 Serial on *see Power Tools*
 Serial output *see Power Tools*
 Serial printer, connect to *see Power Tools*
 Serial port jack 17
 Set 1 Icon 139, 143
 Set 2 Icon 139, 143
 Set Activity 161, 219
 Set time/date on clock 203
 Settings menu 283
 Setup settings, save...283
 Setup word prediction 111
 Shift 13
 Shipping mode 203, 289
 Show status bar 195
 Shut-off, backlight 25, 203
 Shut-off with I/O 11
 Side panel on Pathfinder 21
 Singing 236-259
 pitch notes and values 255
 pitch value code 257
 store song
 pre-stored songs 236-237
 your own song 239-259
 time value of notes 259
 Single meaning icon 147
 Single switch scanning *see 1-switch scanning*
 Software version 12, 15
 Songs, pre-stored 236
 Songs, store 237-259
 Speak display 131, 218 *see also Power Tools*
 Speak notebook 218, *see also Power Tools*
 Speak when spelling 195, 286
 Speak word 194, 221, *see also Power Tools*
 Speaker, external 275
 Speaker grille 17
 Speaker volume down 43, 195
 Speaker volume up 43, 195
 Speaking tools *see Power Tools*
 Speech *see Power Tools*
 Speech dictionary *see Power Tools*
 Speech marker 221, *see also Power Tools*
 Speech menu *see Power Tools*
 Speech on/off 221
 Speech rate, notebook *see Power Tools*
 Speech tips (DECtalk) 261-269
 Spell mode 13
 Spell one word 98, 197
 Spell words
 regular vocabulary 128, 133
 speech dictionary *see Power Tools*
 word prediction dictionary 118-119
 Spell/Minspeak key 128 221
 Spoken prompts 189-193
 Status bar, notebook *see Power Tools*
 Status bar, show 195
 Store
 Activity 151-163
 auditory prompts 189-193

INDEX

computer commands *see Power Tools*
 core vocabulary 129, 133, 151
 digitized speech 129
 IR commands *see Power Tools*
 macros 183-185, 276
 minserts *see Power Tools*
 notebook *see Power Tools*
 songs 237, 239-259
 synthesized speech 133
 themes 133-145
 X-10 commands *see Power Tools*
 Storing tips 136
 Strap 19
 Stuffed Minserts *see Power Tools*
 Swap key 186
 Switch connectors 17
 Switch notebooks *see Power Tools*
 Switch user area 48-49
 Switches
 1-switch 59, 67, 85-86
 2-switch 69, 87-88
 joystick 17, 59, 71, 75, 99
 screen 75
 Synthesized speech
 edit 134
 erase 135
 store 133
 System colors 203
 System information
 key on front of case 12, 15
 maintenance menu 203
 status bar 11, 195
 System lock 203
 System status 12
 key on front of case 12, 15
 maintenance menu 203
 status bar 11, 195
 Table stand 19
 Text
 delete 135
 edit 135
 storing 129, 133
 Text and speech marker 220, *see also Power Tools*
 Text marker *see Power Tools*
 Text sizes 195
 Themes 139-145
 about 139
 exit 145
 pick 1 icon 137, 141
 pick 2 icon 137, 141
 set 1 icon 137, 143
 set 2 icon 137, 143
 store 137-145
 Time/Date, set 203
 Time values of notes 259
 Tips, storing 136
 Toolbox 11, 15, 39, 41, 175
 Tools key on front of case 11, 15
 Tools, list of 39-42, 209-223
 Transfer memory *See MTI Manual*
 Transfer memory menu *see MTI Manual*
 Manual
 Transfer notebook *see MTI Manual*
 Troubleshooting
 IR *see Power Tools*
 Pathfinder 285
 Tutor, icon 199
 Unity *see also Unity Tour Guide and Vocabulary Sort*
 User area 47-54
 copy 53
 erase 54
 icon on front of Pathfinder 12
 rename 52
 select 49
 switch 49
 User area menu 47
 User programmable key 12
 Using the calculator *see Power Tools*
 Using insert tool 140, 175-185, 215

Version number, software 12, 15
View battery status 23, 195

Vocabulary, Unity See *Unity
Pathfinder Manual*

Voice see *Power Tools*

Volume, earphone 43, 195
Volume, speaker 43, 195
Vowel phonemes 253

Wait any key See *Power Tools*

Word Prediction 103-122
and Minspeak 103-105
dictionary 117
add words 118
delete words 118
review words 119
display 111

flash predictor keys 115
frequency list 113
on/off 107
predict in Minspeak 112
predictor keys 115
recency and frequency 113
reset defaults 205
select words 109
serial delete (computer backspace)
121
set up 111
spell 109
turning on word prediction 107
word prediction menu 111

Workspace see *Power Tools*

X-10 see *Power Tools*

